"GOOD" ELECTRIC TOOTHBRUSH NZAOD 600 B1

for EU market:
HOYER Handel GmbH
Tasköprüstraße 3
22761 Hamburg
GERMANY

for GB market:
Lidl Great Britain Ltd
Lidl House · 14 Kingston Road
Surrey · KT8 9NU

Status of information · Informationernes stand · Version des informations
Stand van de informatie · Stand der Informationen:
12/2020 ID: NZAOD 600 B1_20_V1.2

GB **IE** **NI** Operating instructions

**"GOOD" ELECTRIC
TOOTHBRUSH**

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**ELEKTRISCHE ZAHNBÜRSTE
„BASIC“**

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**ELEKTRISK TANDBØRSTE
"BASIC"**

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**ELEKTRISCHE TANDENBORSTEL
'BASIC'**

IAN 352040_2007

IAN 352040_2007


English	2
Dansk	14
Français	26
Nederlands	38
Deutsch	50

Overview / Oversigt / Aperçu de l'appareil / Overzicht / Übersicht

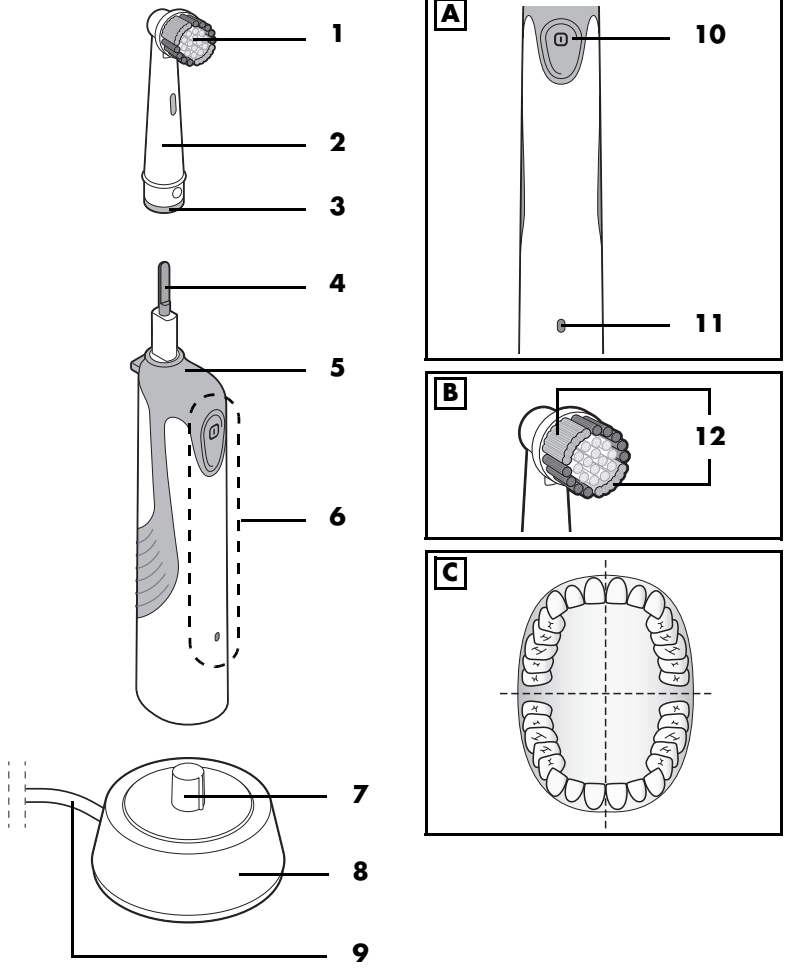Contents

1. Overview	2
2. Intended purpose	3
3. Safety information	3
4. Items supplied	5
5. Setting up the battery charger	5
6. Charging the battery	6
7. Toothbrush heads	6
8. Operation	6
8.1 Modes.....	7
8.2 Brushing teeth.....	7
9. Transport lock	7
10. Cleaning / Hygiene requirements	8
11. Storage	8
12. Disposal	8
13. Trouble-shooting	9
14. Ordering toothbrush heads	9
15. Technical specifications	10
16. Warranty of the HOYER Handel GmbH	11

1. Overview

- 1 Brush head
- 2 Toothbrush head (replaceable)
- 3 Marking ring
- 4 Shaft of the toothbrush
- 5 Toothbrush
- 6 Control elements
- 7 Peg of the battery charger
- 8 Battery charger
- 9 Mains power cable and mains plug

Figure A: control elements

- 10  On/off switch, switching between 2 modes, switching transport lock on and off
- 11 Charging indicator light

Figure B: brush head

- 12 Indicator bristles

Thank you for your trust!

Congratulations on the purchase of your new electric toothbrush.

For safe handling of the device and in order to get to know the entire range of features:

- **Thoroughly read these user instructions prior to initial use.**
- **Above all, observe the safety information!**
- **The device should only be used as described in this copy of the user instructions.**
- **Keep this copy of the user instructions for reference.**
- **If you pass the device on to someone else, please include this copy of the user instructions. The user instructions are a part of the product.**

We hope you enjoy using your new electric toothbrush!

All brand and manufacturer names are the property of their respective owners. They are specified solely in order to simplify the explanation of the points illustrated.

2. Intended purpose

The electric toothbrush is intended solely for tooth cleaning and oral hygiene in humans. The device is designed for private, domestic use. This device must not be used for commercial purposes.

3. Safety information

Warnings

If necessary, the following warnings are used in this copy of the user instructions:

DANGER! High risk: failure to observe this warning may result in injury to life and limb.

WARNING! Moderate risk: failure to observe this warning may result in injury or serious material damage.

CAUTION: low risk: failure to observe this warning may result in minor injury or material damage.

NOTE: circumstances and specifics that must be observed when handling the device.

Instructions for safe operation

- ⊙ This device can be used by children and people with restricted physical, sensory or intellectual abilities or people without adequate experience and/or understanding, if they are supervised or instructed on how to use this device and if they are aware of the resulting hazards. Cleaning and user maintenance must not be performed by children unless they are supervised.
- ⊙ Children must not play with the device.
-  ⊙ The battery charger is connected to the mains voltage, and must therefore not be used in the shower, in the bath, or over washbasins filled with water.
- ⊙ The mains power cable of this device cannot be replaced. If the cable should become damaged, the device itself must be disposed of.
- ⊙ This device contains a battery that cannot be replaced.

DANGER for children!

- ⊙ Packing materials are not children's toys. Children must not be allowed to play with the plastic bags. There is a risk of suffocation.

DANGER to and from pets and livestock!

- ⊙ Electrical devices can represent a danger to pets and livestock. In addition, animals can also cause damage to the device. For this reason you should keep animals away from electrical devices at all times.

DANGER! Risk of electric shock due to moisture!

- ⊙ The toothbrush, the battery charger, the mains power cable and the mains plug must not be immersed in water or any other liquids.

- ⊙ Protect the device from moisture, water drops and water splashes: Failure to observe this instruction may result in an electrical shock hazard.
- ⊙ Never touch the device with wet hands.
- ⊙ Make sure that the surface under the device is dry.
- ⊙ If the battery charger has fallen into water, pull out the mains plug immediately, and only then remove the battery charger from the water. In this case do not use the battery charger again, but have it checked by a specialist workshop.

DANGER! Risk of electric shock!

- ⊙ The installation of a residual current device with a nominal trip current not exceeding 30 mA in the building wiring offers additional protection. Contact a skilled electrician in this regard.
- ⊙ Do not start using the device if the device, the battery charger or the mains power cable shows visible damage, or if the device has previously been dropped.

- ⊙ Only connect the mains plug to a properly installed and easily accessible wall socket whose voltage corresponds to the specifications on the rating plate. The wall socket must still be easily accessible after the device is plugged in.
- ⊙ The hand unit of the toothbrush must only be charged using the battery charger supplied.
- ⊙ Ensure that the mains power cable cannot be caught or crushed, and that it cannot be damaged by sharp edges or hot points.
- ⊙ To fully disconnect the battery charger from the mains power supply, pull out the mains plug.
- ⊙ When removing the mains plug from the wall socket, always pull on the plug and never the cable.
- ⊙ Pull the mains plug out of the wall socket, ...
 - ... if there is a fault,
 - ... when you are not using the device,
 - ... before you clean the device and
 - ... during thunderstorms.
- ⊙ In order to avoid any hazards, do not make modifications to the device. Repairs must be carried out by a specialist workshop or the Service Centre.
- ⊙ When using an extension cable, this must correspond to the applicable safety regulations.

DANGER from batteries!

- ⊙ Protect the battery against mechanical damages. **Fire hazard!**
- ⊙ Do not expose the device to direct sunshine or heat. The ambient temperature should not be lower than 10 °C or higher than 40 °C.
- ⊙ The battery may only be charged with the original accessory (power unit).

CAUTION:

- ⊙ Users could possibly experience mild bleeding of the gums during the first days of electric toothbrush usage. Should the bleeding continue over several days, a dentist should be consulted.
- ⊙ Exert only light pressure when brushing to prevent damage to the gums and to increase the service life of the bristles.

4. Items supplied

- 1 toothbrush **5**
- 2 toothbrush heads **2**
- 1 battery charger **8**
- 1 copy of the user instructions

- Remove all packing materials.
- Check that the device is undamaged.

5. Setting up the battery charger

DANGER! Risk of electric shock due to moisture!

- ⊙ The battery charger **8** is connected to the mains voltage, and must therefore not be used in the shower, in the bath, or over washbasins filled with water.
- Never set up the battery charger **8** anywhere it could fall into water.

6. Charging the battery

NOTES:

- Fully charge the battery before using it for the first time.
- It takes approx. 16 hours to fully charge the battery.
- Place the toothbrush **5** on the peg of the battery charger **7** after use. This ensures that the toothbrush **5** always works at full power. The battery is not damaged by regular charging.

1. Plug the mains plug **9** of the battery charger **8** into a suitable wall socket. The wall socket must still be easily accessible after connection.
2. Place the switched-off toothbrush **5** onto the peg of the battery charger **7**. The charging indicator light **11** is lit red throughout the entire charging process. When the optimal charging state has been reached, the charging indicator light **11** turns green.

Battery signals: charging indicator light 11

Toothbrush **5** on the battery charger **8**:

- Charging indicator light **11** is lit green: the battery is fully charged.
- Charging indicator light **11** is lit red: the battery is being charged.

Toothbrush **5** switched on:

- The toothbrush **5** switches itself off and the charging indicator light **11** flashes red: the battery is weak. Charge the battery.

7. Toothbrush heads

The items supplied include 2 toothbrush heads **2**.

The various marking rings **3** allow you to assign the toothbrush heads **2** to different family members.

- Take a toothbrush head **2** and slide it onto the shaft of the toothbrush **4**.
- In order to remove, carefully pull the toothbrush head **2** off upwards.
- **Figure B:** the indicator bristles **12** on the brush head **1** slowly lose their colour during use. The decolouration shows when the toothbrush heads **2** should be replaced.

NOTE: we recommend that you replace the toothbrush heads **2** every 3 months or after infections (see "Ordering toothbrush heads" on page 9).

8. Operation

CAUTION:

- ⊙ Users could possibly experience mild bleeding of the gums during the first days of electric toothbrush usage. Should the bleeding continue over several days, a dentist should be consulted.
- ⊙ If you have had dental surgery in the last two months, consult a dentist before using the toothbrush.
- ⊙ Exert only light pressure when brushing to prevent damage to the gums and to increase the service life of the bristles.

8.1 Modes

The toothbrush **5** has 2 different modes:

- **Normal mode:**
When switching on with the on/off switch **10** the toothbrush starts with normal mode for cleaning teeth with high pulsation.
- **Sensitive mode:**
By pressing the on/off switch **10** again, the toothbrush switches to the Sensitive mode (gentle brushing) for cleaning teeth with sensitive gums.

NOTE: if you are not used to using an electric toothbrush and you have comparatively sensitive gums, you can use this mode for the first 1 to 2 weeks, and then switch to normal mode later.

8.2 Brushing teeth

1. Remove the toothbrush **5** from the battery charger **8**.
2. Moisten the bristles, squeeze some toothpaste on and switch on the toothbrush **5** only when the brush head **1** is already in your mouth.
3. Press the on/off switch **10** once to start the toothbrush in Normal mode.
4. Press the on/off switch **10** again to switch the toothbrush to Sensitive mode.
5. Slowly guide the brush head **1** from tooth to tooth without exerting too much pressure. Carefully clean all tooth surfaces, the inner surfaces of the teeth and the tooth edges.

Figure C: the recommended brushing time is at least 2 minutes, divided up into 30-second periods for each area of the mouth (bottom right, bottom left, top right, top left). The toothbrush **5** indicates these intervals with 2 brief interruptions. After 2 minutes there are several short interruptions to indicate

that the recommended brushing time has been reached.

6. Press the on/off switch **10** repeatedly until the toothbrush **5** is switched off. Remove the brush head **1** from the mouth only when the toothbrush **5** has been switched off.

NOTES:

- We recommend that you replace the toothbrush heads **2** every 3 months or after infections (see "Ordering toothbrush heads" on page 9).
 - **Figure B:** the indicator bristles **12** on the brush head **1** slowly lose their colour during use. The decolouration shows when the toothbrush heads **2** should be replaced.
-

9. Transport lock

To secure the toothbrush **5** against accidentally being switched on, activate the transport lock.

- Switching on the transport lock:
Press the on/off switch **10** until the charging indicator light **11** flashes red once.
- If you try to switch on the toothbrush **5** with the transport lock activated, the charging indicator light **11** will flash red once.
- Switching off the transport lock:
Press the on/off switch **10** repeatedly until the toothbrush **5** is switched on and the charging indicator light **11** flashes red once. Change the mode if required by pressing the on/off switch **10** again.

10. Cleaning / Hygiene requirements

DANGER! Risk of electric shock!

- ⊙ Before cleaning the battery charger **8**, pull the mains plug **9** out from the mains socket.
 - ⊙ Do not immerse the toothbrush **5** and battery charger **8** in water!
The battery charger **8** must be completely dry before you reconnect it to the wall socket.
-
- To ensure the hygienic functionality of the toothbrush heads **2** they must be rinsed thoroughly after every use and allowed to dry.
 - Remove the toothbrush head **2** once per week and clean it in hot water.
 - Rinse the toothbrush **5** as well and wipe it dry.
 - Wipe the battery charger **8** with a damp cloth if necessary, and rub it dry.

11. Storage

- Disconnect the device from the mains and store it in a dry location.

12. Disposal

This product is subject to the provisions of European Directive 2012/19/EC. The symbol showing a wheelee bin crossed through indicates that the product requires separate refuse collection in the European Union. This applies to the product and all accessories marked with this symbol. Products identified with this symbol may not be discarded with normal household waste, but must be taken to a collection point for recycling electric and electronic appliances.

This recycling symbol marks e.g. an object or parts of the material as valuable for recycling. Recycling helps to reduce the consumption of raw materials and protect the environment.

The battery integrated into this device must not be disposed of in household waste. The device must be disposed of properly with the battery installed.

Packaging

When disposing of the packaging, make sure you comply with the environmental regulations applicable in your country.

13. Trouble-shooting

If your device fails to function as required, please try this checklist first. Perhaps there is only a minor problem, and you can solve it yourself.

DANGER! Risk of electric shock!

- ⊙ Do not attempt to repair the device yourself under any circumstances.

Fault	Possible causes / Action
The device cannot be switched on.	<ul style="list-style-type: none">• Is the battery empty?• Is the transport lock activated? Press the on/off switch 10 repeatedly until the toothbrush 5 is switched on.
Charging indicator light 11 is not lit.	<ul style="list-style-type: none">• The battery charger 8 is not connected to the mains power.

14. Ordering toothbrush heads

You can order additional toothbrush heads **2** on the internet:

shop.hoyerhandel.com

1. Scan the QR code with your smartphone/tablet.
2. With the QR code you will access a website where you can submit your re-order.

15. Technical specifications

Model:	NZAOD 600 B1
Mains voltage:	100-240 V ~ 50/60 Hz
Protection class:	II 
Power rating:	1.4 W standby: 0.4 W
Battery:	1 x 3.7 V  Li ion 600 mAh 2.22 Wh
Charging time:	16 hours
Operating time:	approx. 90 min
Dimensions:	(W x H x D): approx. 65 x 250 x 65 mm
Weight:	approx. 244 g
Protection class of battery charger:	IPX7 (protected against the effects of temporary immersion)
Protection class of hand unit:	IPX4 (protected against water splashes)
LEDs:	RG 0
Ambient conditions:	10 °C to 40 °C humidity 40 to 85 % relative

Symbols used

	Double insulation
	Geprüfte Sicherheit (certified safety). Devices must comply with the generally acknowledged rules of technology and the German Product Safety Act (Produktsicherheitsgesetz - ProdSG). (Not for UK.)
	With the CE marking, HOYER Handel GmbH declares the conformity with EU guidelines.
	This symbol reminds you to dispose of the packaging in an environmentally friendly manner.
	Recyclable materials: cardboard (except corrugated cardboard)
	Direct current
	AC voltage

Subject to technical modifications.

NOTE: no action by the user is required to switch the product between 50 and 60 Hz. The product adjusts itself for both 50 and 60 Hz.

16. Warranty of the HOYER Handel GmbH

Dear Customer,
your device is provided with a 3 year warranty starting with the purchase date. In the event of product defects, you are entitled to statutory rights against the vendor. These statutory rights are not restricted by our warranty presented in the following.

Warranty conditions

The warranty period starts with the purchase date. Please keep the original purchase receipt in a safe place. This document is required to verify the purchase.

If within three years from the purchase date of this product a material or factory defect occurs, the product will be repaired or replaced by us – at our discretion – free of charge to you, or the purchase price will be refunded. This warranty implies that within the period of three years the defective device and the purchase receipt are presented, including a brief written description of the defect and the time it occurred.

If the defect is covered by our warranty, the repaired or a new product will be returned to you. No new warranty period starts with a repair or replacement of the product.

Warranty period and statutory claims for defects

The warranty period is not extended when the warranty has been claimed. This also applies to replaced and repaired parts. Any damages and defects already existing at the time of purchase must be reported immediately upon unpacking. Repairs arising after the expiration of the warranty period are subject to a charge.

Warranty coverage

The device was produced carefully according to strict quality guidelines and tested diligently prior to delivery.

The warranty applies to material or factory defects.

Excluded from the warranty are wear parts subject to normal wear and damages to fragile parts, e.g. switches, batteries, lamps or other parts manufactured from glass.

This warranty expires if the product is damaged, not used as intended or not serviced. For the proper operation of the product, all instructions listed in the operating instructions must be observed carefully. Any form of use and handling that is advised against in the operating instructions or warned against must always be avoided.

The product is only intended for private and not for commercial use. In the case of incorrect and improper treatment, use of force and interventions not performed by our authorised Service Centre, the warranty shall cease.

Handling in case of a warranty claim

In order to ensure prompt processing of your matter, please observe the following notes:

- Please keep the article number **IAN: 352040_2007** and the purchase receipt as a purchase verification for all inquiries.
- The article number can be found on the rating plate, an engraving, the title page of your instructions (in the bottom left), or as a label on the rear or underside of the device.
- If faulty operation or other defects occur, first contact the Service Centre listed in the following by **telephone** or **email**.
- Then, you are able to send a product reported as defective free of charge to the service address specified to you, including the purchase receipt and the information on the defect and when it occurred.

At www.lidl-service.com you can download this and many other manuals, product videos and installation software.

This QR code takes you directly to the Lidl Service website (www.lidl-service.com), where you can enter the article number (IAN) **352040_2007** to open your user instructions.

Supplier

Please note that the following address is **no service address**. First contact the aforementioned Service Centre.

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