## NACKEN-MASSAGEKISSEN MIT AKKU / CORDLESS NECK MASSAGER / COUSSIN MASSANT SANS FIL POUR NUQUE SMKA 2000 A4

DE AT CH

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Tel.: 0800562153

E-Mail: [owim@lidl.ch](mailto:owim@lidl.ch)

<b>List of pictograms used</b> .....	Page 16
<b>Introduction</b> .....	Page 16
Intended use.....	Page 16
Description of parts.....	Page 16
Technical data.....	Page 17
Scope of delivery.....	Page 17
<b>Safety advice</b> .....	Page 17
Safety instructions for batteries / rechargeable batteries.....	Page 19
<b>Before use</b> .....	Page 20
Charging the rechargeable battery.....	Page 20
<b>Operation</b> .....	Page 20
Switch off the product.....	Page 21
<b>Cleaning and care</b> .....	Page 21
<b>Storage</b> .....	Page 21
<b>Troubleshooting</b> .....	Page 22
<b>Disposal</b> .....	Page 22
<b>Warranty</b> .....	Page 23
Warranty claim procedure.....	Page 23
Service.....	Page 23

List of pictograms used			
	Only use in dry indoor rooms.		Safety class II
	Low/High rotation speed	<b>IP20</b>	Protection against penetration of solid objects with a diameter > 12.5 mm
	ON/OFF button		Clockwise/Anti-clockwise rotation
<b>CE</b>	CE mark indicates conformity with relevant EU directives applicable for this product.	<b>ta</b>	Highest rated ambient temperature
	Safety information Instructions for use		Do not use the charger if the pins of the plug are damaged.
	SMPS (switched-mode power supply)		Enclosed safety transformer

## Cordless neck massager SMKA 2000 A4

### ● Introduction

We congratulate you on the purchase of your new product. You have chosen a high quality product. The instructions for use are part of the product. They contain important information concerning safety, use and disposal. Before using the product, please familiarise yourself with all of the safety information and instructions for use. Only use the product as described and for the specified applications. If you pass the product on to anyone else, please ensure that you also pass on all the documentation with it.

### ● Intended use

The product is only designed for massaging the neck of people. The product is only for indoor household use. It is not intended for medical or commercial purposes and is not a substitute for medical treatment. Any other form of use is not permitted and may result in injury or damage of the product. The

manufacturer accepts no liability for any damage caused by improper use.

### ● Description of parts

- 1 Massage pillow
- 2 Massage heads
- 3 Adjustment strap with snap fasteners
- 4 Socket
- 5 Control panel
- 6 Clockwise/anti-clockwise rotation button with indicator light
- 7 ON/OFF button with indicator light
- 8 Low/High rotation speed button with indicator light
- 9 USB port
- 10 Plug-in mains adapter
- 11 USB type A plug
- 12 USB cable
- 13 USB type C plug

## ● Technical data

### Plug-in mains adapter

Manufacturer's name or trade mark, commercial registration number and address:	Xiamen Keli Electronics Co.,Ltd 91350211737865674J No.19, XingLin North 3rd Road, Jimei District, Xiamen, Fujian, China.
Model identifier:	SW-050100
Input voltage:	100-240V~
Input AC frequency:	50/60Hz
Output voltage:	5.0V===
Output current:	1.0A
Output power:	5.0W
Average active efficiency:	73.7%
No-load power consumption:	0.10W
Input current:	0.5A
Highest rated ambient temperature (ta):	45°C
Protection class:	II/ 

### Massage pillow

Input voltage:	5.0V===
Input current:	1.0A
Rechargeable battery:	3.7V/2000mAh Lithium-ion battery (non-replaceable)
Duration of use:	max. 15 min (automatic shut-off)

## ● Scope of delivery

Immediately after unpacking, please check the package contents for completeness and if the product and all parts are in good condition. Remove all packaging materials before use.

- 1 Cordless neck massager (with rechargeable battery)
- 1 Plug-in mains adapter

- 1 USB cable
- 1 Instruction manual

## Safety advice

Read all safety advices and instructions. Non-observance of the safety advices and instructions may cause electric shock, fire and/or bad injuries. No action is needed from users to shift the product between 50 Hz and 60 Hz. The product adapts itself for both 50 Hz and 60 Hz.

**KEEP ALL SAFETY NOTICES AND INSTRUCTIONS FOR FUTURE REFERENCE!**

-  **⚠ WARNING!**  
**DANGER TO LIFE AND RISK OF ACCIDENT FOR INFANTS AND CHILDREN!** Never leave children unsupervised with the packaging materials. The packaging material poses a suffocation hazard. Children frequently underestimate the dangers. Please always keep the product out of the reach of children.

**⚠ CAUTION!** This product is not a toy for children! Children are

not aware of the dangers associated with handling electrical products.

- This product can be used by children aged from 14 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the product in a safe way and understand the hazards involved. Children shall not play with the product. Cleaning and user maintenance shall not be made by children unless they are older than 14 and supervised.

**⚠ WARNING!** Strangulation hazard. The product contains adjustment strap with snap fasteners **3** that is placed around the neck. Make sure that the adjustment strap with snap fasteners **3** cannot be trapped when moving. Otherwise there is a risk of strangulation and injury.

- The output cord cannot be replaced. If the cord is damaged the product should be scrapped.

**⚠ WARNING!** Avoid the hazard of life by electric shock!

**⚠ WARNING!** For the purposes of recharging the battery, only use the detachable supply unit provided with this product.

- Do not insert needles or pointed objects into the massage pillow **1**.
- Never immerse the product into water or any liquid.
- Never touch the product, or the plug-in mains adapter **10** with wet hands or use it in very humid environments.
- Only use the provided plug-in mains adapter **10**.
- Make sure that the details on the rating label are in conformity with your mains voltage.
- In the event of malfunctions, remove the plug-in mains adapter **10** from the mains socket.
- Do not open the product or attempt to repair it yourself. Defective products should only be repaired by qualified personal.
- Always remove the plug-in mains adapter **10** from the mains socket after each use and before clean the product.
- Keep the USB cable **12** away from hot surfaces, sharp edges and mechanical forces.
- Allow the product to cool down before moving it to a different location.

## **WARNING!** **HAZARD OF FIRE!**

- Do not operate the product close to or below curtains, wall cupboards or other inflammable materials.
- Do not use the product on swollen, burnt, inflamed, diseased or injured regions of the skin and body areas. If in doubt, obtain medical advice before using the product.
- Never use the product if you have injuries on neck or neck pain.

## **Safety instructions for batteries/ rechargeable batteries**

- **DANGER TO LIFE!** Keep batteries/rechargeable batteries out of reach of children. If accidentally swallowed seek immediate medical attention.
-  **DANGER OF EXPLOSION!** Never recharge non-rechargeable batteries. Do not short-circuit batteries/rechargeable batteries and/or open them. Overheating, fire or bursting can be the result.
- Never throw batteries/rechargeable batteries into fire or water.

- Do not exert mechanical loads to batteries/rechargeable batteries.

## **Risk of leakage of batteries/ rechargeable batteries**

- Avoid extreme environmental conditions and temperatures, which could affect batteries/rechargeable batteries, e.g. radiators/direct sunlight.
- If batteries/rechargeable batteries have leaked, avoid contact with skin, eyes and mucous membranes with the chemicals! Flush immediately the affected areas with fresh water and seek medical attention!
-  **WEAR PROTECTIVE GLOVES!** Leaked or damaged batteries/rechargeable batteries can cause burns on contact with the skin. Wear suitable protective gloves at all times if such an event occurs.
- This product has a built-in rechargeable battery which cannot be replaced by the user. The removal or replacement of the rechargeable battery may only be carried out by the manufacturer or his customer service or by a similarly qualified person in order to avoid hazards. When disposing of the product, it should

be noted that this product contains a rechargeable battery.

## ● Before use

**Note:** Remove all packaging material from the product.

**Note:** The product must be powered using the built-in rechargeable battery.

## ● Charging the rechargeable battery

**Note:** Before first use, the rechargeable battery in the product should be charged for approximate 120 minutes. A charging time of approximate 120 minutes is also required subsequently.

**Note:** The product can be charged by using the plug-in mains adapter [10] and the USB cable [12] only.

**Note:** During charging, the ON/OFF button [7] flashes red. When the product is fully charged, the ON/OFF button [7] lights up continuously red.

**Note:** During charging, the product cannot be operated.

**Note:** When the rechargeable battery is almost empty, the massager heads [2] do not rotate and the ON/OFF button [7] lights red and flashes for around 15 seconds.

- To avoid injuries and damage to the product, always switch the product off before charging the battery. Press the ON/OFF button [7] once and then the green indicator light turns off (see Fig. B).
- Insert the USB type C plug [13] of the USB cable [12] into the socket [4] of the product (see Fig. B and D).
- Insert the USB type A plug [11] of the USB cable [12] into the USB port [9] of the plug-in mains adapter [10] (see Fig. C and D).
- Insert the plug-in mains adapter [10] (see Fig. C) into a properly installed socket which meets the product requirements (see "Technical data").

**The indicator light on On/Off button [7] shows the status of the rechargeable battery as follows:**

Battery status	Indicator light color
Charging	Red light flashes.
Fully charged	Red light keeps on.
High battery	Green light keeps on
Medium battery	Yellow light keeps on
Low battery	Red light flashes for around 15 seconds

## ● Operation

- Never use the product for more than 15 minutes. A longer massage may lead to tension caused by overstimulation of the muscles.
- The product switches off automatically after approximate 15 minutes. Cool it down for at least 15 minutes before switching on again to avoid overheating.
- Never use the product while sleeping or laying down. Extra mechanical load on the neck against the massage heads [2] may cause injuries or pain in the neck.
- Place the massage pillow [1] around your neck and fasten it using the adjustment strap with snap fasteners [3] (see Fig. A).
- Switch on the massage pillow [1] by pressing the ON/OFF button [7] on the control panel [5] for approximate 1 second (see Fig. B). The indicator light on ON/OFF button [7] turns on and lights green constantly. The massage heads [2] rotates in a clockwise direction.
- Press the clockwise/anti-clockwise rotation button [6] on the control panel [5] to reverse the direction of rotation of the massage heads [2] (see Fig. B). The massage heads [2] now turn anti-clockwise and the indicator light color turns from blue to green.
- Press the low/high rotation speed button [8] on the control panel [5] to increase the rotation speed of the massage heads [2] (see Fig. B). The speed of rotation of the massage heads [2]

can be adjusted in two levels. Blue indicator light color is low speed, green indicator light is higher speed.

## ● Switch off the product

- Press the ON/OFF button **7** on the control panel **5** to switch off the product (see Fig. B). All indicator lights on control panel **5** turns off and the massage heads **2** stop rotating. The product switches off automatically after approximate 15 minutes.
- Cool it down for at least 15 minutes before switching on again to avoid overheating.

**The indicator light on control panel **5** shows the operating status of the product as follows:**

Status	On/off button <b>7</b>
The product is switched on.	Green light turns on and lights constantly.
The product is switched off.	Green light turns off.

Status	Clockwise/anti-clockwise rotation button <b>6</b>
Clockwise rotating	Blue light
Anti-clockwise rotating	Green light

Status	Low/High rotation speed button <b>8</b>
Low speed	Blue light
High speed	Green light

the USB type C plug **13** from the socket **4** before cleaning.

- **⚠ WARNING!** Never immerse the product, or plug-in mains adapter **10** into water.
- Ensure the massage pillow **1** is completely dry before connecting the plug-in mains adapter **10**.
- Do not use abrasive, aggressive or chemical cleaning agents to clean the massage pillow **1**. It could damage the product.
- Never allow any liquid to enter the massage pillow **1**.
- Use a damp lint-free cloth to clean the product. Use a little mild detergent on the cloth if necessary.

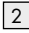
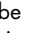
## ● Storage

- Let the massage pillow **1** completely cool down before storage.
- Do not place other objects on top of the massage pillow **1**.
- If the product is not used for a certain period of time, store the product in a clean and dry location.

## ● Cleaning and care

-  **⚠ WARNING! DANGER TO LIFE BY ELECTRIC SHOCK!**  
Disconnect the plug-in mains adapter **10** from the mains socket and unplug

## ● Troubleshooting

Error	Cause	Solution
The product does not work.	The battery may not be charged.	Charge the battery (see chapter “Charging the rechargeable battery”).
The product switches off during use.	The product switches off automatically after approximate 15 minutes to avoid overheating.	Let the product cool down for at least 15 minutes before switching it on again.
The massage heads  rotates slowly.	The mechanical load on the massage heads  is too high.	Do not press the neck against the massage heads  too hard.
The product cannot be switched on by pressing the ON/OFF button  .	The ON/OFF button  is not pressed long enough.	Press the ON/OFF button  for approximate 1 second.
The product cannot be charged.	The battery level is sufficient and thus no charging is required.	Charge the product only when the battery level is medium or low.

## ● Disposal

The packaging is made entirely of recyclable materials, which you may dispose of at local recycling facilities.

Observe the marking of the packaging materials for waste separation, which are marked with abbreviations (a) and numbers (b) with following meaning: 1–7: plastics/20–22: paper and fibre-board/80–98: composite materials.

Contact your local refuse disposal authority for more details of how to dispose of your worn-out product.

To help protect the environment, please dispose of the product properly when it has reached the end of its useful life and not in the household waste. Information on collection points and their opening hours can be obtained from your local authority.

The product incl. accessories and packaging materials are recyclable and are subject to extended producer responsibility. Dispose them separately, following the illustrated Info-tri (sorting information), for better waste treatment. The Triman logo is valid in France only.

Faulty or used batteries/rechargeable batteries must be recycled in accordance with Directive 2006/66/EC and its amendments. Please return the batteries/rechargeable batteries and/or the product to the available collection points.

### Environmental damage through incorrect disposal of the batteries/rechargeable batteries!

Batteries/rechargeable batteries may not be disposed of with the usual domestic waste. They may contain toxic heavy metals and are subject to hazardous waste treatment rules and regulations. The chemical symbols for heavy metals are as follows: Cd = cadmium, Hg = mercury, Pb = lead. That is why you should dispose of used batteries/rechargeable batteries at a local collection point.

## **This disposal is free of charge to you.**

Protect the environment and dispose of properly.

The product is recyclable, subject to extended manufacturer responsibility, and collected separately.

The built-in rechargeable battery cannot be removed for disposal. Return the product completely to a collection site for used electronics.

## **Warranty**

The product has been manufactured to strict quality guidelines and meticulously examined before delivery. In the event of material or manufacturing defects you have legal rights against the retailer of this product. Your legal rights are not limited in any way by our warranty detailed below.

The warranty for this product is 3 years from the date of purchase. The warranty period begins on the date of purchase. Keep the original sales receipt in a safe location as this document is required as proof of purchase.

Any damage or defects already present at the time of purchase must be reported without delay after unpacking the product.

Should the product show any fault in materials or manufacture within 3 years from the date of purchase, we will repair or replace it – at our choice – free of charge to you. The warranty period is not extended as a result of a claim being granted. This also applies to replaced and repaired parts.

This warranty becomes void if the product has been damaged, or used or maintained improperly.

The warranty covers material or manufacturing defects. This warranty does not cover product parts subject to normal wear and tear, thus considered consumables (e.g. batteries, rechargeable batteries, tubes, cartridges), nor damage to fragile parts, e.g. switches or glass parts.

## **Warranty claim procedure**

To ensure quick processing of your case, please observe the following instructions:

Please have the till receipt and the item number (IAN 435808\_2304) available as proof of purchase.

You will find the item number on the rating plate, an engraving on the front page of the instructions for use (bottom left), or as a sticker on the rear or bottom of the product.

If functional or other defects occur, please contact the service department listed either by telephone or by e-mail.

You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

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**OWIM GmbH & Co. KG**