SHIATSU MASSAGE CUSHION SHIATSU-MASSAGEKISSEN COUSSIN DE MASSAGE SHIATSU SSMK 18 A1

(GB) (IE)

SHIATSU MASSAGE CUSHION

Operating instructions

(FR) (BE)

COUSSIN DE MASSAGE SHIATSU

Mode d'emploi

(CZ)

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(SK)

MASÁŽNY VANKÚŠ SHIATSU

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(DK)

SHIATSU MASSAGEPUDE

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(HU)

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(DE) (AT) (CH)

SHIATSU-MASSAGEKISSEN

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(NL) (BE)

SHIATSU-MASSAGEKUSSEN

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(PL)

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(ES)

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
(IT)

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IAN 467756_2401

(DE) (FR) (BE) (NL)
(CZ) (PL) (SK)

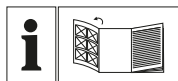GB IE

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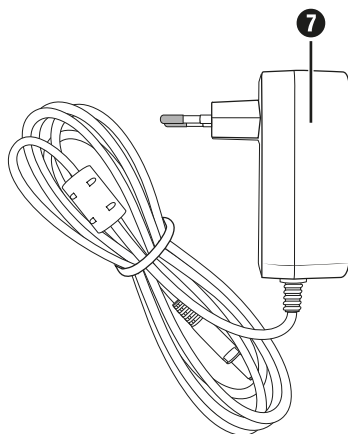
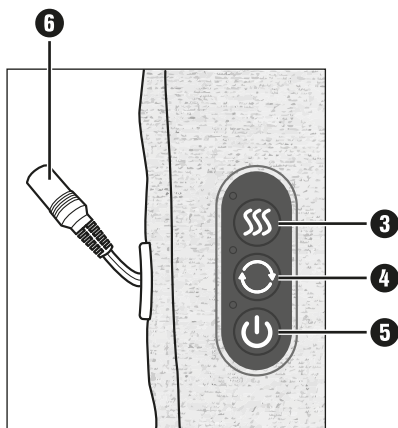
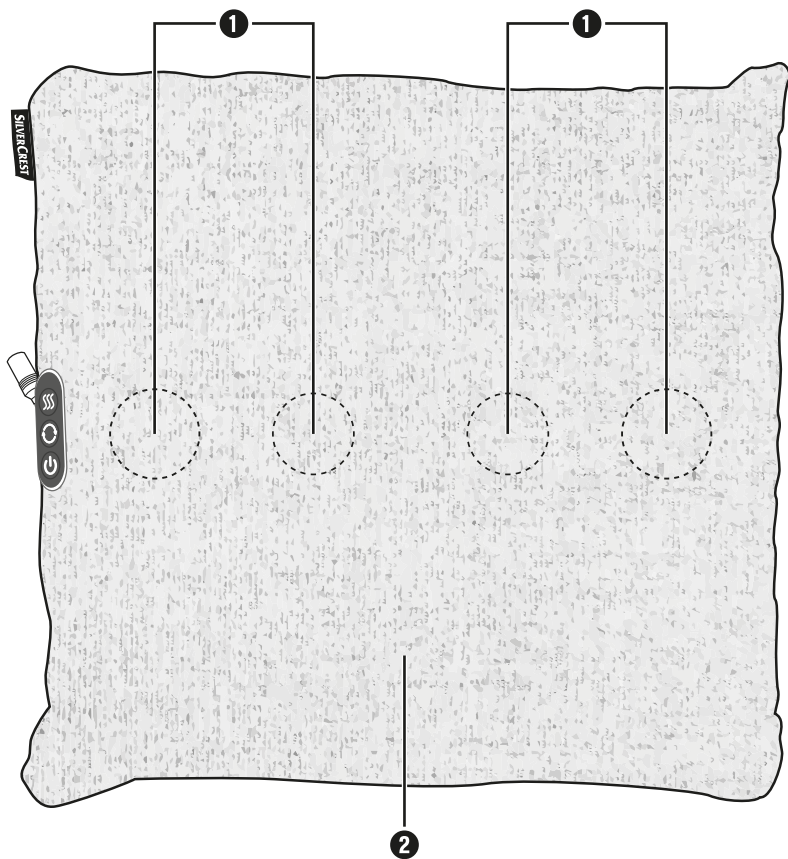
Contents

Introduction	2
Information about these instructions for use	2
Intended use	2
Warnings and symbols used	2
Safety information	3
Operating elements/parts	5
Before first use	5
Check package contents	5
Using the appliance	6
Operation and use	6
Switching the appliance on/off	6
LED operating light	6
Switching the heat function on/off	6
Changing the direction of the rotating massage heads	6
Cleaning	7
Cleaning the appliance	7
Washing the removable cover	7
Storage	8
Troubleshooting	8
Disposal	8
Disposal of the appliance	8
Disposal of the packaging	8
Appendix	9
Technical data	9
Kompernass Handels GmbH warranty	9
Service	10
Importer	10

Introduction

Information about these instructions for use

Congratulations on the purchase of your new appliance. You have selected a high-quality appliance. The instructions for use are part of this appliance. They contain important information about safety, usage and disposal. Before using the appliance, please familiarise yourself with all operating and safety instructions. Use this appliance only as described and for the areas of applications specified. Hand over all documents when passing the appliance on to third parties.

Intended use

This appliance is designed to relax and loosen muscles. The appliance is versatile and can be used on the shoulders, neck, back and legs. You can also activate the heat function for a more enhanced effect.

Any other use or use beyond that specified is considered improper use. The appliance is not intended for use in medical/therapeutic or commercial settings and cannot replace medical treatment. Claims of any type due to damage from improper use are excluded. The risk is borne solely by the user.

Warnings and symbols used

The following warnings and symbols are used in these instructions for use, on the packaging and on the appliance:

	DANGER! A warning with this symbol and the signal word "DANGER" indicates an imminently hazardous situation which will result in death or serious injury if not avoided.
	WARNING! A warning with this symbol and the signal word "WARNING" indicates a potentially hazardous situation which could result in death or serious injury if not avoided.

	CAUTION! A warning with this symbol and the signal word "CAUTION" indicates a potentially hazardous situation that could result in a minor or moderate injury if not avoided.
	ATTENTION! A warning with this symbol and the signal word "ATTENTION" indicates a potential situation which could result in property damage if not avoided.
	Note: A note provides additional information aimed at making it easier to use the appliance.
	Protection class II: Protection by double or reinforced insulation between live and touchable parts.
	Use the appliance indoors only.
	DC current/voltage
	AC current/voltage
	Only the removable cover is machine washable at 30 °C.
	Do not bleach
	Do not tumble dry
	Do not iron
	Do not dry clean

Safety information

⚠ DANGER!

- Do not allow children to play with the packaging material! Keep all packaging materials away from children. Risk of suffocation!
- Prior to use, check the appliance for visible external damage. Do not use an appliance that has been damaged or dropped.
- Only use the mains adapter in dry indoor areas.
- Connect the mains adapter only to correctly installed and earthed power outlets. Ensure that the power supply corresponds with the information on the type plate of the appliance/mains adapter.
- Only use the appliance with the mains adapter supplied and only use the mains adapter supplied with the appliance.
- Protect the appliance, mains adapter and permanently connected cord from moisture, dripping or splashing water. Never immerse the appliance, mains adapter or permanently connected cord in water or other liquids.
- Keep the appliance and mains adapter away from hot surfaces.
- Route the connection cord so that it cannot be pinched or damaged.
- The power outlet used must always be easily accessible so that the mains adapter can be quickly removed from the power outlet in a hazardous situation.
- Do not use the mains adapter if it has a damaged plug or connection cord.
- The mains adapter and permanently connected connection cord may not be repaired. In the event of a defect, the entire mains adapter must be replaced with another of identical design.
- Switch the appliance off immediately and disconnect the mains adapter from the appliance if you smell fumes or notice any smoke. Have the appliance checked by qualified specialists before using it again.
- Switch off the appliance after each use and before cleaning it. For additional protection, the appliance has an automatic switch-off function after 15 minutes.
- Hold the mains adapter to disconnect the appliance from the mains power. Do not pull on the permanently connected cord.

- Always unplug the mains adapter from the power outlet before cleaning.
- Do not carry out any repair work on the appliance. All repairs must be performed by our Customer Service or by a qualified specialist technician.
- Never open the zip of the white cushion cover. There are no parts requiring maintenance inside. It is only used for technical production reasons.
- People with pacemakers, artificial joints or electronic implants should consult their doctor before using the appliance.
- Do not use the appliance around your heart, face or genital area.
- Do not use the appliance if you suffer from one or more of the following illnesses or complaints: circulatory disorders, varicose veins, open wounds, bruises, skin tears, phlebitis, thrombosis, diabetes or slipped discs.

⚠ WARNING!

- This appliance may be used by children aged 8 years or older and by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, as long as they are being supervised or have been instructed on how to use the appliance safely and are aware of the potential risks.
- Do not allow children to use the appliance as a toy.
- Cleaning and user maintenance may not be carried out by children unless they are supervised.
- Never use the appliance on children, animals, during pregnancy, on your head, throat or other parts of the body that are exhibiting swelling, burns, inflammation, eczema, wounds or other sensitive areas.
- The appliance has a hot surface. Persons who are insensitive to heat must take care when using the appliance.
- Do not use the appliance under blankets or pillows. Never use the appliance in the vicinity of petrol or other inflammable materials.

⚠ CAUTION!

- Do not use the appliance while sleeping or after taking substances that may reduce your ability to react quickly (e.g. painkillers or alcohol).
- Never use the appliance directly on bones, e.g. on your spine, shins, the tops of your feet or similar places on your body.
- Use the appliance only on dry, clean parts of your body.
- Stop using the appliance immediately if you find the massage unpleasant or if you experience pain.
- Use the appliance only for its intended purpose.

⚠ ATTENTION!

- Only use the appliance indoors.
- Do not use any aggressive cleaning liquids or solvents to clean the appliance. They could damage the surface.

ⓘ Note:

- No action is required on the part of the user to set the appliance to 50 or 60 Hz. The appliance automatically sets itself to 50 or 60 Hz.

Operating elements/parts

(See fold-out page for illustrations)

- ❶ Rotating massage head 4×
- ❷ Removable cover
- ❸ SSS button (heat function)
- ❹ ⌚ button (internal and external rotation)
- ❺ ⏻ On/Off button
- ❻ Connection socket
- ❼ Mains adapter

Before first use**Check package contents**

⚠ DANGER! Risk of suffocation! Do not allow children to play with packaging materials.

- ◆ Remove all parts and the instructions for use from the box.
- ◆ Remove all packaging material.

The product includes the following components:

- Shiatsu massage cushion
- Mains adapter
- These instructions for use

ⓘ Note: Check the package for completeness and signs of visible damage. If the delivery is incomplete or damage has occurred as a result of defective packaging or during transport, contact the service hotline (see section **Service**).

Using the appliance

The appliance allows you to massage your shoulders, neck, back and legs to relax them and take stiffness away. You can also activate the heat function for a more enhanced effect. Depending on the extent of the tension in the part of the body to be treated, a maximum treatment time of about 15 minutes is recommended. The appliance has an automatic switch-off function, which switches the appliance off automatically after approx. 15 minutes.

To begin with, only press the part of the body you want to massage lightly against the appliance. Then check whether the massage position is comfortable for you. After this, gradually shift your body weight towards the appliance and increase the contact pressure.

A massage should always feel pleasant. A slight reddening of the treated areas and even slight bruising are a normal side effect, as the massage specifically promotes blood circulation. However, if you experience excessive skin reactions, discontinue use and consult a doctor. If you find the massage painful or unpleasant, change your position or the contact pressure or stop the treatment.

Operation and use

Switching the appliance on/off

- ◆ Insert the connector plug of the mains adapter ⑦ into the connection socket ⑥ on the appliance.
- ◆ Plug the mains adapter ⑦ into a power outlet.
- ◆ Press the  On/Off button ⑤ to start the appliance. The massage heads ① start rotating automatically. The LED operating light is switched on.
- ◆ Press the  On/Off button ⑤ to switch off the appliance. The LED operating light goes out.
- ① **Note:** The appliance switches off automatically after approx. 15 minutes.
- ◆ Unplug the mains adapter ⑦ from the power outlet when you are finished using the appliance.
- ◆ Disconnect the mains adapter plug ⑦ from the connection socket ⑥ on the appliance.

LED operating light

LED	Colour	State/function
	Red	Heat function switched on
	Off	Heat function switched off
	Green	External rotation
	Red	Internal rotation
	Green	Appliance switched on
	Off	Appliance switched off

Switching the heat function on/off

- ① **Note:** When you switch on the appliance, the heat function is switched off.
- ◆ Press the  button ③ to switch on the heat function. The light for the rotating massage heads ① lights up. The LED operating light of the  button ③ lights up red.
- ◆ Press the  button ③ to switch off the heat function. The light for the rotating massage heads ① goes out. The LED operating light on the  button ③ also goes out.

Changing the direction of the rotating massage heads

Two massage heads, each rotating in opposite directions in pairs, have internal and external rotation.

- ① **Note:** When you switch on the appliance, the massage heads ① rotate in an external rotation.
- ◆ Press the  button ④ to change the direction to internal rotation. The LED operating light of the  button ④ lights up red.
- ◆ Press the  button ④ to change the direction to external rotation. The LED operating light of the  button ④ lights up green.

Cleaning

Cleaning the appliance

⚠ DANGER! Risk of fatal electrical shock!

Switch off the appliance before cleaning. Always unplug the mains adapter **7** from the power outlet before cleaning. Disconnect the mains adapter plug **7** from the connection socket **6** on the appliance. Ensure that no moisture can get into the appliance during cleaning.

- ⓘ ATTENTION!** Damage to the appliance!
Do not use caustic, abrasive or solvent-based cleaning agents. They can damage the surfaces of the appliance.
- ◆ Wipe off the appliance using a damp cloth. Do not use any solvents or cleansers which can damage plastic. If necessary, put a little mild detergent on the cloth and wipe this off afterwards with clear water.
- ◆ Allow everything to dry thoroughly before storing the appliance or putting it back into operation.

Washing the removable cover

If it is very dirty, you can wash the removable cover **2** in a washing machine at 30 °C. Observe the symbols for care and cleaning on the sewn-on label on the removable cover **2**.

- ⓘ Note:** Never open the zip of the white cushion cover. It is only used for technical production reasons.
- ◆ Open the zip of the removable cover **2**.
- ◆ Remove the white cushion from the removable cover **2**.
- ◆ From the inside of the removable cover **2**, pull the hook-and-loop fastener off the control panel (see **1** Fig. 1).
- ◆ Remove the circuit board (see **2** Fig. 1).
- ◆ Pull the connection socket **6** through the recess in the removable cover **2** (see **3** Fig. 1).

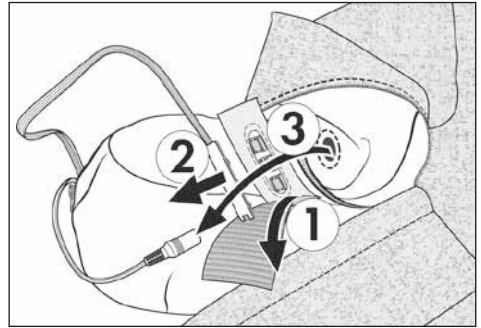
Fig. 1

- ◆ Place the removable cover **2** in the washing machine.
- ◆ Allow the removable cover **2** to dry thoroughly after the wash cycle.
- ◆ Place the removable cover **2** on the white cushion so that the connection socket **6** and the control panel are on the same side.
- ◆ Plug the circuit board onto the control panel. Observe the correct direction (see **1** Fig. 2).
- ◆ Attach the hook-and-loop fastener (see **2** Fig. 2).
- ◆ Guide the connection socket **6** through the recess on the removable cover **2** (see **3** Fig. 2).

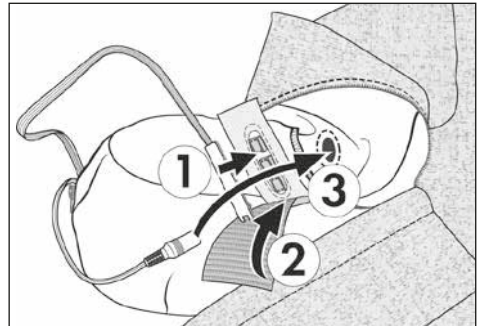
Fig. 2

- ◆ Put the removable cover **2** over the white cushion.
- ◆ Close the zip of the removable cover **2**.

Storage

- ◆ Unplug the mains adapter ⑦ from the power outlet when not in use.
- ◆ Disconnect the mains adapter plug ⑦ from the connection socket ⑥ on the appliance.
- ◆ When not in use, store the appliance in a clean, dry place out of direct sunlight.

Troubleshooting

Fault	Cause	Remedy
The massage heads ① rotate very slowly.	The massage heads ① are overloaded.	Reduce the load on the massage heads ①.
The massage heads ① do not move.	The appliance is not switched on or not connected to the mains.	Establish a connection to the power supply and/or switch on the appliance.

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
The product, its packaging and the operating instructions are recyclable. They are subject to an extended manufacturer responsibility and will be collected separately.

Disposal of the appliance

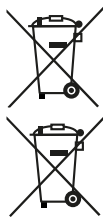The adjacent symbol of a crossed-out dustbin means that this appliance is subject to Directive 2012/19/EU. This directive states that this appliance may not be disposed of in the normal household waste at the end of its useful life, but must be taken to specially set-up collection locations, recycling depots or disposal companies.

The disposal is free of charge for the user. Protect the environment and dispose of this appliance properly.

If your old appliance has stored any personal data, you are responsible for deleting it yourself before returning it.

Your local community or municipal authorities can provide information on how to dispose of the worn-out product.

Disposal of the packaging

The packaging materials have been selected for their environmental friendliness and ease of disposal and are therefore recyclable. Dispose of packaging materials that are no longer needed in accordance with applicable local regulations.


Dispose of the packaging in an environmentally friendly manner. Note the labelling on the packaging and separate the packaging material components for disposal, if necessary.

The packaging material is labelled with abbreviations (a) and numbers (b) with the following meanings: 1-7: plastics, 20-22: paper and cardboard, 80-98: composites.

Appendix

Technical data

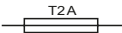Mains adapter	
Manufacturer	
Xiamen Invcy Electronics Co., Ltd. 1F/3F-5F, BLDG 5, No.943-5, Tonglong 2nd Rd, Torch Hi-tech Industrial Zone (Xiang'an), Xiamen, 361000, P.R. China Commercial register number: 91350200568403208H Importer: KOMPERNASS HANDELS GMBH, BURGSTRASSE 21, 44867 BOCHUM, GERMANY, Register court AG Bochum, Register number: HRB 4598	
Model	IVP1200-1500GA
Input voltage	100–240 V ~
Input AC frequency	50/60 Hz
Output voltage	12.0 V ===
Output current	1.5 A
Output power	18.0 W
Average efficiency in operation	86.6%
Low load efficiency (10%)	80.3%
Power consumption at no load	0.08 W
Power consumption	0.6 A
Protection class	II/□ (double insulation)
Polarity	 (positive inside, negative outside)
Efficiency class 6	
Mains adapter	

Appliance	
Input voltage	12.0 V ===
Input power	18.0 W
Polarity	 (positive inside, negative outside)
Removable power supply	

Kompernass Handels GmbH warranty

Dear Customer,

This appliance has a 3-year warranty valid from the date of purchase. If this product has any faults, you, the buyer, have certain statutory rights. Your statutory rights are not restricted in any way by the warranty described below.

Warranty conditions

The warranty period starts on the date of purchase. Please keep your receipt in a safe place. This will be required as proof of purchase.

If any material or manufacturing fault occurs within three years of the date of purchase of the product, we will either repair or replace the product for you or refund the purchase price (at our discretion). This warranty service requires that you present the defective appliance and the proof of purchase (receipt) within the three-year warranty period, along with a brief written description of the fault and of when it occurred.

If the defect is covered by the warranty, your product will either be repaired or replaced by us. The repair or replacement of a product does not signify the beginning of a new warranty period.

Warranty period and statutory claims for defects

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking. Repairs carried out after expiry of the warranty period shall be subject to a fee.

Scope of the warranty

This appliance has been manufactured in accordance with strict quality guidelines and inspected meticulously prior to delivery.

The warranty covers material faults or production faults. The warranty does not extend to product parts subject to normal wear and tear or to fragile parts which could be considered as consumable parts such as switches or parts made of glass.

The warranty does not apply if the product has been damaged, improperly used or improperly maintained. The directions in the operating instructions for the product regarding proper use of the product are to be strictly followed. Uses and actions that are discouraged in the operating instructions or which are warned against must be avoided.

This product is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications/repairs which have not been carried out by one of our authorised Service centres.

Warranty claim procedure

To ensure quick processing of your case, please observe the following instructions:

- Please have the till receipt and the item number (IAN) 467756_2401 available as proof of purchase.
- You will find the item number on the type plate on the product, an engraving on the product, on the front page of the operating instructions (below left) or on the sticker on the rear or bottom of the product.
- If functional or other defects occur, please contact the service department listed either by **telephone** or by **e-mail**.

- You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

You can download these instructions along with many other manuals, product videos and installation software at www.lidl-service.com.

This QR code will take you directly to the Lidl service page (www.lidl-service.com) where you can open your operating instructions by entering the item number (IAN) 467756_2401.

Service

GB Service Great Britain

Tel.: 0800 404 7657

E-Mail: kompernass@lidl.co.uk

IE Service Ireland

Tel.: 1800 101010

E-Mail: kompernass@lidl.ie

IAN 467756_2401

Importer

Please note that the following address is not the service address. Please use the service address provided in the operating instructions.

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