WATER JET FLOSSER / MUNDUSCHE HYDROPULSEUR SMD 24 A1

(GB) (IE)

WATER JET FLOSSER

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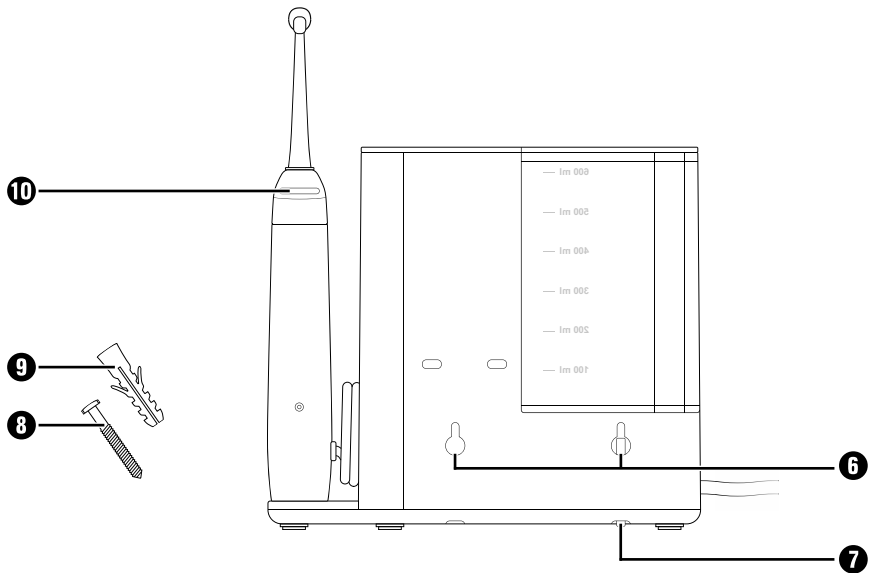
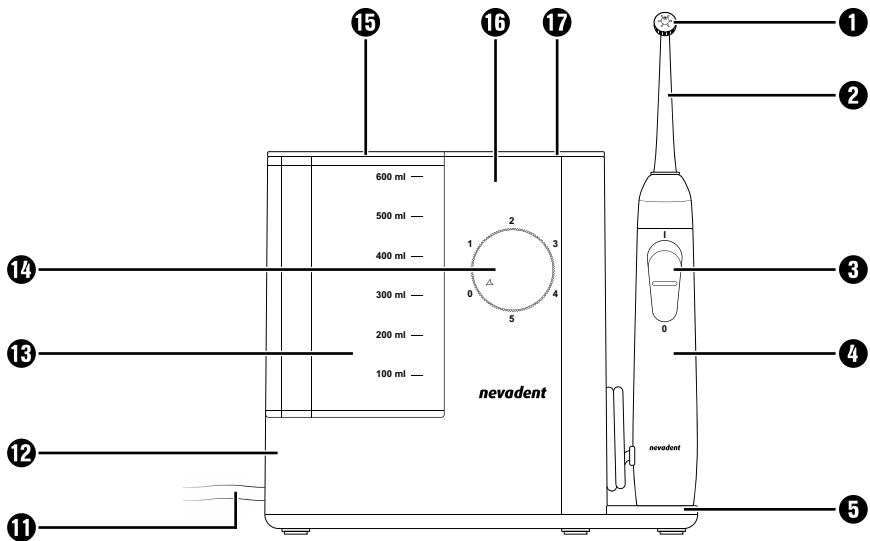
Contents

Introduction	2
Information about this short manual	2
Intended use	2
Safety information	3
Package contents and transport inspection	5
Description of the appliance	6
Before first use	7
Requirements for the set-up location	7
Wall mounting	7
Handling and operation	8
Attaching/removing nozzle attachments	8
Adjusting the nozzle control	8
Using the oral irrigator	9
Disposal	10
Disposal of the appliance	10
Disposal of the packaging	11
Appendix	11
Technical data	11
Service	11
Importer	12

Introduction

Information about this short manual

This document is a shortened print version of the complete operating instructions.

Scanning the QR code will take you directly to the Lidl Service page (www.lidl-service.com), where you can view and download the full version of the operating instructions by entering the article number (IAN) 385287_2107.

WARNING!

Observe the information in the complete operating instructions and the safety instructions to avoid personal injury and damage to property.

The short manual is an integral part of this product. Before using the product, please familiarise yourself with all operating and safety instructions. Keep the short manual in a safe place and hand over all documents when passing the product on to any future owner/user. I documents when passing the product on to any future owner/user.

Intended use

This appliance is intended exclusively for dental care and mouth hygiene in humans. The appliance is not intended for use in industrial or medical environments. It is not suitable for animal care. This appliance is intended solely for use in private households. It is not intended for any other purpose nor for use beyond the scope described.

Claims of any kind for damage resulting from misuse, incompetent repairs, unauthorised modification or the use of unauthorised replacement parts will not be accepted.

The operator bears sole liability.

Safety information

DANGER! ELECTRIC SHOCK!

- ▶ Connect the appliance only to a properly installed and easily accessible mains power socket supplying a mains power voltage of 100–240 V ~, 50/60 Hz.
- ▶ Ensure that the power cable is not damaged. Keep it out of hot areas and route it in such a way that it cannot be damaged.
- ▶ Before use, check the appliance for visible external damage. Never use the appliance if the cable, plug or housing are damaged.
- ▶ To avoid potential risks, if the appliance power cable is damaged, it must be replaced by the manufacturer, its customer service or by a qualified technician.
- ▶ Have all repairs carried out by a specialist workshop. Under no circumstances should you open the appliance yourself. Repairs that are not carried out by a specialist workshop could lead to physical injuries.
- ▶ Do not use the appliance in the vicinity of a bathtub/shower or a sink filled with water. The proximity of moisture presents a danger even when the appliance is switched off.
- ▶ Never touch the appliance if it has fallen into water. Always disconnect the power plug from the power socket first!
- ▶ Never immerse the connected appliance in water or other liquids.
- ▶ Never touch the power cable or the appliance with wet or damp hands.

⚠ DANGER! ELECTRIC SHOCK!

- ▶ Ensure that the power cable can never become wet or damp during operation. Route the power cable in such a way that it cannot become trapped or damaged in any other way.
- ▶ Do not use the power cable with an extension cable; connect the power cable directly to a power socket.

⚠ WARNING! RISK OF INJURY!

- ▶ Do not change any accessories while the appliance is switched on.
- ▶ Always store the appliance indoors. To prevent accidents, keep the appliance in a dry location when not in use.
- ▶ Oral irrigators may be used by children over the age of 8 and by persons with reduced physical, sensory or mental capabilities or lack of experience and/or knowledge if they are supervised. Do not allow children to use the appliance as a toy.
- ▶ Switch off the appliance before you put it down or change accessories.
- ▶ Lay the power cable so that nobody can step on it or trip over it.
- ▶ Always pull out the plug before cleaning the appliance. Do not pull on the power cable. Always grip the power plug to disconnect the appliance from the power supply.

ATTENTION! PROPERTY DAMAGE!

- ▶ Use only accessories recommended by the manufacturer.
- ▶ Do not make any unauthorised modifications or alterations to the appliance.
- ▶ Never put down the appliance next to radiators, ovens or other heated appliances or surfaces.
- ▶ The appliance may only be used with cold or lukewarm tap water up to a maximum of 40°C . If necessary, add a few drops of mouthwash. The water tank should never be completely filled with additives such as mouthwash or mouth rinse. These can damage the appliance.

Package contents and transport inspection

- ◆ Remove all parts of the appliance and the short manual from the carton.
- ◆ Remove all packaging materials and any films and labels.

DANGER

- ▶ Do not allow children to play with packaging materials. There is a risk of suffocation!

The package contents include the following components (see fold-out page for illustrations):

- Water jet flosser (main unit with hand unit and water tank)
- 4 nozzle attachments
- 2 wall plugs (∅ 5 mm)
- 2 screws (3 mm)
- Short manual (not shown)

NOTE

- ▶ Check the package for completeness and for signs of visible damage.
- ▶ If the delivery is incomplete or damage has occurred as a result of defective packaging or during transport, contact the service hotline (see section **Service**).

Description of the appliance

(See fold-out page for illustrations)

- ① Nozzle control
- ② Nozzle attachment
- ③ On/Off switch for the water jet I/O
- ④ Hand unit
- ⑤ Cradle for hand unit
- ⑥ Hangers
- ⑦ Locking pin
- ⑧ Screw
- ⑨ Wall plug
- ⑩ Release button for nozzle attachment
- ⑪ Power cable (with mains plug)
- ⑫ Main unit
- ⑬ Water tank
- ⑭ Rotary switch (On/Off switch and water pressure control)
- ⑮ Water tank lid
- ⑯ Storage compartment for nozzle attachments
- ⑰ Storage compartment lid

Before first use

Requirements for the set-up location

To ensure safe and trouble-free operation of the appliance, the set-up location must fulfil the following requirements:

- ▶ The appliance may only be operated free-standing if used in dry rooms. The set-up location must be a horizontal, level and stable surface.
- ▶ In damp rooms, especially in the vicinity of water, the appliance must be firmly fixed to the wall to protect it from falling into water.
- ▶ The mains power socket must be easily accessible so that the power cable can be quickly disconnected in an emergency.

Wall mounting

NOTE

- ▶ Before mounting to the wall, make sure that there are no electric cables, wires or pipes running behind the wall surface where you plan to hang the appliance.

1. Mark two drill holes on the wall 7.2 cm apart and drill the marked holes with a 5 mm bit. Ensure that there is a power socket in the immediate vicinity.
2. Push the two supplied wall plugs **9** into the drilled holes and screw the two supplied screws **8** into the wall plugs **9**. Take care that the screws **8** still protrude about 3 mm out of the holes.
3. On the back side of the appliance, the right hanger **6** has a locking pin **7** to prevent the appliance from coming free of the screws **8**. Press the locking pin **7** back somewhat until you hear a click. Then pull the locking pin **7** down slightly out of the hanger **6**.
4. Push the appliance onto the screws **8** so that the screws **8** enter into the hangers **6** on the back side of the appliance. Then press the appliance downward so that the screws **8** slide into the rails of the hangers **6** and the appliance rests firmly on the screws **8**.
5. Slide the locking pin **7** upwards until it audibly engages with the hanger **6**.
6. Plug the mains plug of the power cable **11** into a socket.

Handling and operation

Attaching/removing nozzle attachments

NOTE

- ▶ If used daily, replace the nozzle attachments **2** after no more than six months.
- ▶ Use only original nozzle attachments **2**.
- ◆ Make certain that the appliance is switched off before attaching or removing the nozzle attachment **2**.
- ◆ Fit a nozzle attachment **2** onto the hand unit **4**. Press the nozzle attachment **2** into the opening until the nozzle attachment **2** audibly engages in the hand unit **4**. Ensure that the guide tab on the lower end of the nozzle attachment **2** is aligned with the corresponding groove in the opening on the hand unit **4**.
- ◆ To remove the nozzle attachment **2**, press the release button **10** on the hand unit **4** and then pull the nozzle attachment **2** out of the hand unit **4**.

Adjusting the nozzle control

You can use the nozzle control **1** on the head of the nozzle attachment **2** to switch between two different nozzle functions:

- ◆ Turn the nozzle control **1** anti-clockwise with your fingers to obtain a single jet for targeted removal of food remnants from hard-to-reach areas.
- ◆ Turn the nozzle control **1** clockwise with your fingers to obtain multiple jets for massaging the gums.

Using the oral irrigator

⚠ CAUTION – PROPERTY DAMAGE!

- ▶ Fill the water tank **13** with pure tap water. If necessary, add a few drops of mouthwash. The water tank **13** should never be completely filled with additives such as mouthwash or mouth rinse.
- ▶ Use only cold or lukewarm water with a maximum temperature of 40°C. Do not fill the water tank **13** with hot water.

NOTE

- ▶ Use of the appliance is not a replacement for daily tooth cleaning with a toothbrush and toothpaste. Use the appliance only as a supplement to daily tooth cleaning with toothbrush and toothpaste.
 - ▶ Only use the appliance for a maximum of 2 minutes or only long enough to use up one full tank of water (600 ml). One full tank of water (600 ml) is sufficient to operate the appliance for about 150 seconds at level **1** or about 110 seconds at level **5**.
 - ▶ The unfamiliar stimulation may cause some gum bleeding at first, which is not unusual. If this continues for more than two weeks, you should consult your dentist.
 - ▶ Consult your dentist before using the appliance if you suffer from severe periodontitis, have injuries or abscesses in the mucous membrane of the oral cavity or underwent dental surgery within the last two months.
- ◆ Remove the water tank **13** from the main unit **12** and open the lid **15** of the water tank **13**. Fill the water tank **13** with cold or lukewarm water. Take care not to fill the water tank **13** past the maximum fill mark of **600 ml**. Then close the water tank **13** again with the lid **15** and place the water tank **13** in the main unit **12**.
 - ◆ Fit a nozzle attachment **2** onto the hand unit **4**. Press the nozzle attachment **2** into the opening until the nozzle attachment **2** audibly engages in the hand unit **4**. Ensure that the guide tab on the lower end of the nozzle attachment **2** is aligned with the corresponding groove in the opening on the hand unit **4**.
 - ◆ Switch on the main unit **12** by setting the rotary switch **14** to the desired water pressure level **1–5** (**1** = gentle, **5** = strong).
 - ◆ Open your mouth slightly and insert the nozzle attachment **2**. Bend over a washbasin to allow the water to run into the basin.

NOTE

- ▶ Begin with a lower level to allow your gums to become accustomed to the water jet. Increase the strength of the water jet in stages. If the pressure is uncomfortable, select a lower level again.
- ▶ To avoid splashes, do not switch on the hand unit ④ until it is in your mouth.
- ◆ Switch on the hand unit ④ by moving the On/Off switch for the water jet I/O ③ on the hand unit ④ to I.
- ◆ Direct the water jet at the gaps between your teeth and the edges of your teeth to clean them. Position the nozzle attachment ② so that the water jet can rinse the food remnants away from the gums. Move the nozzle attachment ② slowly along the edge of the teeth from gap to gap. Clean all the front and back sides of your teeth in this way.
- ◆ Start at the front side of your upper molars and move the nozzle attachment ② along the edges of the teeth from one side to the other. Repeat this process for the back of the upper teeth and for the back and front of the lower teeth.

Disposal

Disposal of the appliance

The adjacent symbol of a crossed-out wheelee bin means that this appliance is subject to Directive 2012/19/EU. This directive states that the product may not be disposed of in normal household waste at the end of its useful life, but must be taken to a specially set-up collection point, recycling depot or disposal company.

This disposal is free of charge for the user. Protect the environment and dispose of this appliance properly.

Your local community or municipal authorities can provide information on how to dispose of the worn-out product.


The product is recyclable, subject to extended producer responsibility and is collected separately.

Disposal of the packaging

The packaging materials have been selected for their environmental compatibility and ease of disposal and are therefore recyclable.

Dispose of packaging materials that are no longer needed in accordance with applicable local regulations.

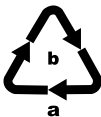Dispose of the packaging in an environmentally friendly manner.

Note the labelling on the packaging and separate the packaging material components for disposal, if necessary. The packaging material is labelled with abbreviations (a) and numbers (b) with the following meanings:

1-7: plastics, 20-22: paper and cardboard, 80-98: composites.

Appendix

Technical data

Power supply	100-240 V ~ (AC), 50/60 Hz
Power consumption	24 W
Protection class	II/□ (double insulation)
Protection type	IPX4 (protected against water splashes from any direction)

NOTE

- ▶ No user action is required to switch the product between 50 and 60 Hz. The product switches automatically to either 50 or 60 Hz.

Service

- GB Service Great Britain**
Tel.: 0800 404 7657
E-Mail: kompernass@lidl.co.uk

- IE Service Ireland**
Tel.: 1890 930 034
(0,08 EUR / Min., (peak))
(0,06 EUR / Min., (off peak))
E-Mail: kompernass@lidl.ie

IAN 385287_2107

Importer

Please note that the following address is not the service address. Please use the service address provided in the operating instructions.

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