WÄRMEKISSEN **HEAT PAD** **COUSSIN CHAUFFANT**

(DE) (AT) (CH)

WÄRMEKISSEN

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(FR) (BE)

COUSSIN CHAUFFANT

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(PL)

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(SK)

TERMOVANKŮŠ

Krátky návod

(DK)

VARMEPUDE

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(GB) (IE)

HEAT PAD

Quick Guide

(NL) (BE)

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(CZ)

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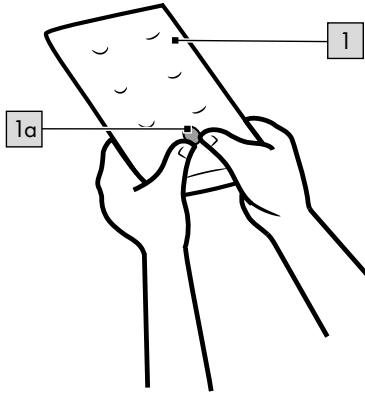
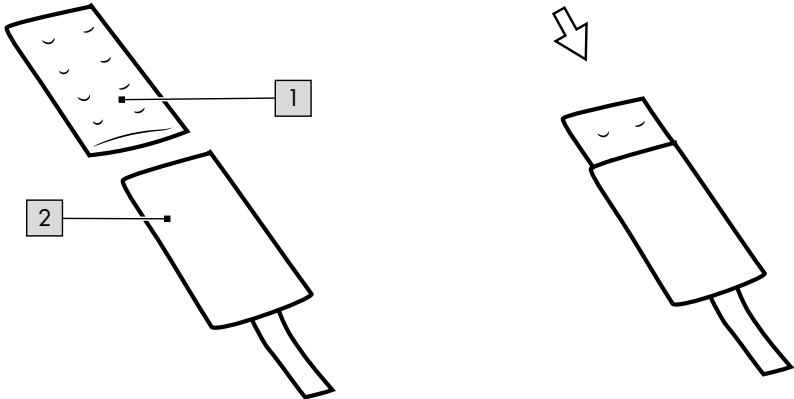
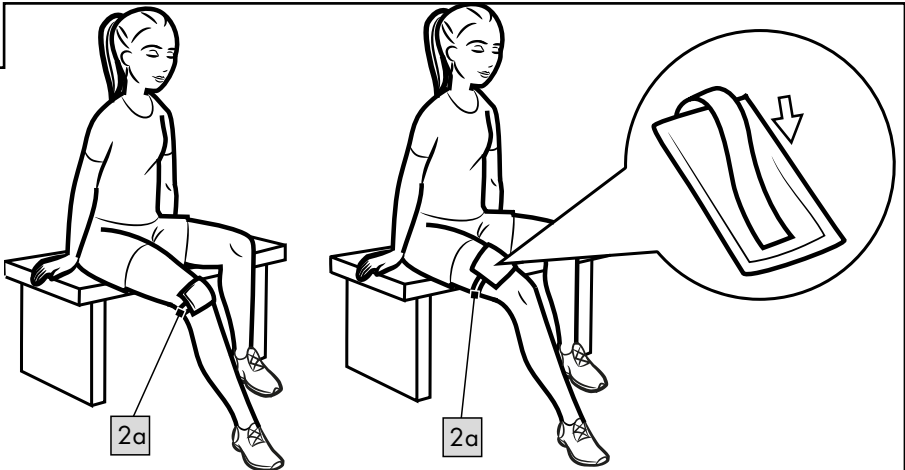
(ES)

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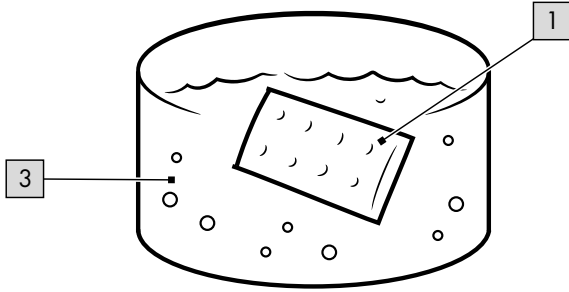
Guía rápida

IAN 376085_2101

(DE) (AT) (CH) (GB) (IE) (FR) (BE)
(NL) (PL) (CZ) (SK) (ES) (DK)

A**B****C**

D

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This document is an abridged print version of the full instructions for use. By scanning the QR code you gain direct access to the Lidl service website (www.lidl-service.com), and by entering the product number (IAN) 376085_2101 you can look at the full instructions for use and download them.

WARNING! Adhere to the full instructions for use and the safety information in order to avoid personal injuries and damage to property. The quick guide is an inherent part of this product. Before using the product familiarise yourself with all usage and safety information. Keep the quick guide handy and pass it on along with all documents when handing over the product to third parties.

Scope of delivery

2 x heat pad (1)
2 x cover (2)
1 x instructions for use

Technical data

Dimensions: approx. 29 x 14.5cm (l x w)/
approx. 36.5 x 14.5cm (l x w)
Heating temperature: max. 55 °C

 Date of manufacture (month/year):
08/2021

Intended use

The product has been designed as a heating pad for private use.
Not suitable for charging in the microwave.

Safety instructions

Danger to life!

- Never leave children unsupervised with the packaging material. Risk of suffocation.
- **WARNING.** Not a toy! Keep out of the reach of children! Use only under supervision.
- Check the product for damage and wear before each use.
The product may only be used in good order and condition!

Risk of injury!

- If handled carefully and correctly, no liquid should leak out. Should this happen, contact with the skin may cause skin irritations. Rinse such areas thoroughly with water.
- If the liquid makes contact with the eyes or other mucous membranes, rinse them thoroughly and consult a doctor.

Danger of scalding!

- The product heats up to a temperature of 55 °C. Do not wear the product directly against the skin. Instead, use the cover supplied.

Avoiding material damage!

- Keep the product away from pointed and sharp objects. If in exceptional circumstances liquid leaks out, avoid contact with your garments because the liquid can be removed only with difficulty with standard commercial detergents.
- Do not charge the product in the microwave because the product could explode.

Use

Two heating pads (1) and two covers (2) are included in the delivery. To illustrate, here is an explanation of how to use a heating pad and a cover. The heating pads and covers differ only in size.

1. Remove the packaging and spread out the product (1) and the cover (2).
2. Bend the metal platelet (1a) that is located in the product (1) (fig. A). The liquid begins to harden and the product becomes warm.

Note: If the liquid does not harden, bend the metal platelet again. If you cannot recognise the metal platelet because of the product's colour, hold the product up against the light.

3. When the product has become warm, place it in the cover (2) (fig. B) and close the cover at the opening using the hook-and-loop fastener.
4. Position the product on your knee or your thigh, as shown in fig. C.

5. Lead the strap (2a) (fig. C) through under your knee or thigh and attach the ends on the cover (2).

Note: Make sure that the strap is not fastened too tightly. The product gives off heat for up to 60 min.

Charging

You can charge the product around 1,000 times.

1. Remove the product from the cover (2).
2. Put the product in boiling water (3) for ten minutes until the inside of the product becomes clear and liquid again (fig. D).
3. Let the water cool down.
4. Remove the product and dry it off. The product is now ready for use again.

Storage, cleaning

When not in use, always store the product clean and dry at room temperature. Wipe clean with a damp cloth only.

IMPORTANT! Never clean the product with harsh cleaning agents.

The cover (2) can be machine-washed at 40 °C. Before washing, take the product (1) out of the cover and close the hook-and-loop fasteners!

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