## KIDS' ELECTRIC TOOTHBRUSH NKZ 3 A1

(GB) (IE) (NI)

### KIDS' ELECTRIC TOOTHBRUSH

Operating instructions

(NL) (BE)

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
(GB) (IE) (NI)  
(DK) (NL)

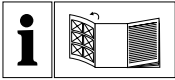GB IE NI

Before reading, unfold the page containing the illustrations and familiarise yourself with all functions of the device.

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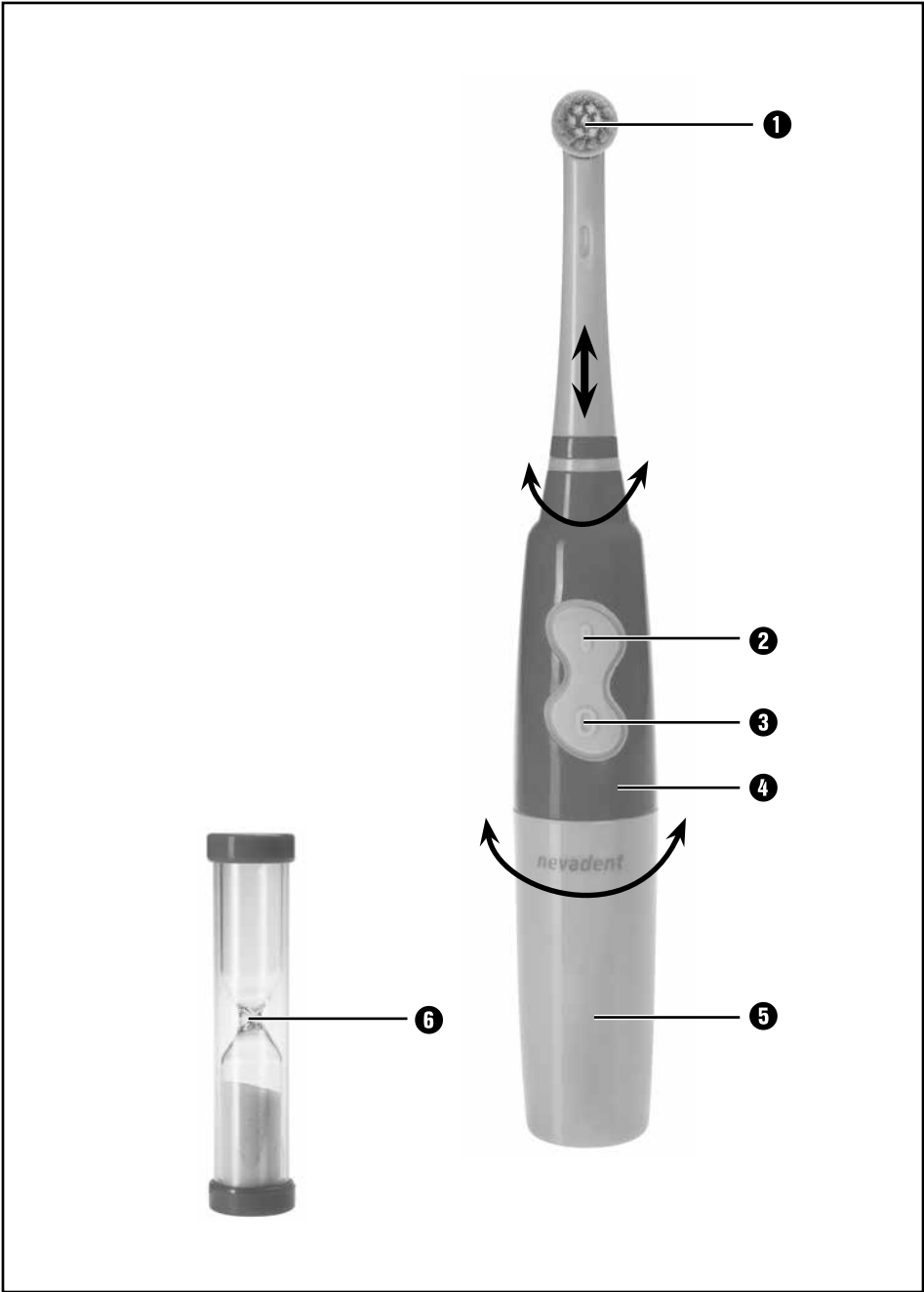
# Contents

<b>Introduction</b> .....	<b>2</b>
<b>Intended use</b> .....	<b>2</b>
<b>Package contents</b> .....	<b>2</b>
<b>Description of the appliance</b> .....	<b>2</b>
<b>Technical data</b> .....	<b>2</b>
<b>Safety instructions</b> .....	<b>3</b>
<b>Preparation</b> .....	<b>4</b>
<b>Fitting/removing brush heads</b> .....	<b>4</b>
<b>Brushing teeth</b> .....	<b>4</b>
<b>Cleaning and care</b> .....	<b>5</b>
<b>Disposal</b> .....	<b>6</b>
<b>Kompernass Handels GmbH warranty</b> .....	<b>6</b>
Service .....	7
Importer .....	7
<b>Ordering replacement parts</b> .....	<b>8</b>

## Introduction

Congratulations on the purchase of your new appliance.

You have selected a high-quality product. The operating instructions are part of this product. They contain important information on safety, usage and disposal. Before using the product, familiarise yourself with all operating and safety instructions. Use the product only as described and for the specified areas of application. Please also pass on these operating instructions to any future owner.

## Intended use

This appliance is intended exclusively for personal dental care at home. This appliance is intended solely for use in private households. The appliance is not intended for use in industrial or medical environments. It is not suitable for animal care. Do not use it for commercial purposes.

## Package contents

Hand element  
4 brush heads  
2 batteries  
1 sticker sheet with 4 stickers  
Sandglass  
Storage box  
Operating instructions

### NOTE

Check contents for completeness and damage immediately after unpacking. Contact Customer Services if necessary.

Remove all packaging materials and plastic film from the appliance.

## Description of the appliance

- 1 Brush head
- 2 On switch
- 3 Off switch
- 4 Hand element
- 5 Battery compartment cap
- 6 Sandglass

## Technical data

Voltage supply	2 x 1.5 V $\text{---}$ (direct current)
Battery type	1.5 V, micro/AAA/LR03
Protection class	IPX4 (protected against water splashes)

## Safety instructions

### **⚠ WARNING! RISK OF INJURY!**

- ▶ To avoid damage, do not allow liquids to penetrate the appliance and do not submerge it in water.
- ▶ Ensure that no foreign objects penetrate the appliance.
- ▶ In the event of malfunctions or obvious damage, please contact our Customer Service department.
- ▶ Toothbrushes may be used by children and by persons with reduced physical, sensory or mental capabilities or lack of experience and/or knowledge if they are supervised or have been told how to use the appliance safely and are aware of the potential risks.
- ▶ Cleaning and user maintenance should not be performed by children unless they are under supervision.
- ▶ Children must not use the appliance as a plaything.
- ▶ Do not allow the appliance to come into contact with hot surfaces. The appliance could be irreparably damaged!
- ▶ Do not open the appliance (except the battery compartment) yourself. Risk of injury if appliance is not properly assembled!

### **HANDLING BATTERIES**

- ▶ Do not throw batteries into a fire. Do not recharge batteries. There is a risk of explosion and injury!
- ▶ Never open batteries and never solder or weld batteries. There is a risk of explosion and injury!
- ▶ Check the condition of the batteries at regular intervals. Leaking batteries can cause damage to the appliance.

## HANDLING BATTERIES

- ▶ If you do not intend to use the appliance for an extended period, remove the batteries.
- ▶ If the batteries have leaked, use protective gloves. Clean the battery compartment and the battery contacts with a dry cloth.
- ▶ Keep batteries out of the reach of children. Children might put batteries into their mouth and swallow them. If a battery is swallowed, medical assistance must be sought immediately.

## Preparation

Unscrew the battery compartment cap ⑤.

- 1) Insert the batteries as per the polarity diagram in the battery compartment.
- 2) Screw the battery compartment cap ⑤ back on.
- 3) Select a sticker and stick it on the front of the appliance.
- 4) Before initial use, clean the new brush heads ① with fresh water.

## Fitting/removing brush heads

- Push a brush head ① onto the hand element ④ with a slight twisting motion.
- To remove the brush head ①, turn it a little anticlockwise while pulling it upwards and off.

### NOTE

To avoid damage to the appliance, always use original brush heads.

(See section "Order replacement parts" for replacement brush heads)

## Brushing teeth

This appliance is not intended for children under 3 years of age.

### ⚠ WARNING!

Due to the unaccustomed effect of the toothbrush, there may be some initial gum bleeding. This is normal. If this continues for more than two weeks, you should consult your dentist.


- 1) Place the brush head ① onto the drive shaft on the handset ④.
- 2) Moisten the brush under flowing water.
- 3) Put a little toothpaste onto the wet brush.

### NOTE

To avoid spray, do not switch on the appliance until it is in your mouth

- 4) Turn the sandglass ⑥ over and press the on switch ②.

- 5) Divide your mouth into four parts (quadrants): upper right, upper left, lower right and lower left.

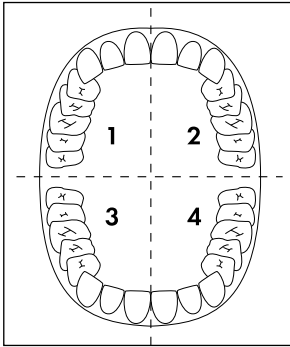- 6) Start in one quadrant:  
Start with the brush head **1** on your gums, and move the brush upwards/downwards from the gums to the crown of the tooth. Proceed in this manner on both the inside and the outside of the tooth. Start again at the gums for each individual tooth. Do not forget to brush the back of the rear molars. To clean the chewing surfaces, guide the brush over them with a little pressure.
- 7) After about 30 seconds, move on to the next quadrant. Clean each jaw quadrant using the same technique.
- 8) After about two minutes, the sand in the sand-glass **6** will have run through, indicating that you have completed the two minutes of cleaning recommended by dental professionals.
- 9) Press the off switch **3**. The toothbrush stops.
- 10) Then take the toothbrush out of your mouth.

- 11) After use, rinse the brush head **1** thoroughly under running water.
- 12) Remove the brush head **1** from the handset **4**. To clean the drive shaft and the handset **4**, rinse them briefly with water. Hold the handset **4** so that the drive shaft is at a downward angle.

## Cleaning and care

### **WARNING! PROPERTY DAMAGE!**

-  Never submerge the hand element **4** in water or other liquids!
- Do not use detergents or solvents. They can damage the plastic surfaces.

- 1) Remove the brush head **1** from the handset **4**. After use, clean the brush heads **1** thoroughly under running water. We recommend cleaning the brush heads **1** in a standard commercial denture cleaning fluid every 2-4 weeks. A new brush head **1** should be fitted after about 3 months, and the brush head must always be replaced when the bristles start to bend outwards.
- 2) To clean the drive shaft and the handset **4**, rinse them briefly with water. Hold the handset **4** so that the drive shaft is at a downwards angle.

## Disposal

**Never dispose of the appliance in the normal domestic waste. This product is subject to the provisions of European Directive 2012/19/EU.**

Dispose of the appliance via an approved waste disposal company or your municipal waste disposal facility. Please comply with all applicable regulations. Please contact your waste disposal centre if you are in any doubt.

Your local community or municipal authorities can provide information on how to dispose of the worn-out product.


### Disposal of (rechargeable) batteries

Used batteries may not be disposed of in the household waste. Consumers are legally obliged to dispose of (rechargeable) batteries at a collection point in their community/city district or at a retail store. The purpose of this requirement is to ensure that batteries are disposed of in an environmentally friendly manner. Only dispose of batteries when they are fully discharged.

The packaging is made from environmentally friendly material and can be disposed of at your local recycling plant.

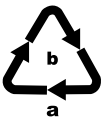Dispose of the packaging in an environmentally friendly manner. Note the labelling on the packaging and separate the packaging material components for disposal if necessary. The packaging material is labelled with abbreviations (a) and numbers (b) with the following meanings: 1-7: Plastics, 20-22: Paper and cardboard, 80-98: Composites.

## Kompernass Handels GmbH warranty

Dear Customer,

This appliance has a 3-year warranty valid from the date of purchase. If this product has any faults, you, the buyer, have certain statutory rights. Your statutory rights are not restricted in any way by the warranty described below.

### Warranty conditions

The validity period of the warranty starts from the date of purchase. Please keep your original receipt in a safe place. This document will be required as proof of purchase.

If any material or production fault occurs within three years of the date of purchase of the product, we will either repair or replace the product for you at our discretion. This warranty service is dependent on you presenting the defective appliance and the proof of purchase (receipt) and a short written description of the fault and its time of occurrence.

If the defect is covered by the warranty, your product will either be repaired or replaced by us. The repair or replacement of a product does not signify the beginning of a new warranty period.

### Warranty period and statutory claims for defects

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking. Repairs carried out after expiry of the warranty period shall be subject to a fee.

## Scope of the warranty

This appliance has been manufactured in accordance with strict quality guidelines and inspected meticulously prior to delivery.

The warranty covers material faults or production faults. The warranty does not extend to product parts subject to normal wear and tear or fragile parts such as switches, batteries, baking moulds or parts made of glass.

The warranty does not apply if the product has been damaged, improperly used or improperly maintained. The directions in the operating instructions for the product regarding proper use of the product are to be strictly followed. Uses and actions that are discouraged in the operating instructions or which are warned against must be avoided.

This product is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications/repairs which have not been carried out by one of our authorised Service centres.

## Warranty claim procedure

To ensure quick processing of your case, please observe the following instructions:

- Please have the till receipt and the item number (e.g. IAN 12345) available as proof of purchase.
- You will find the item number on the type plate, an engraving on the front page of the instructions (bottom left), or as a sticker on the rear or bottom of the appliance.
- If functional or other defects occur, please contact the service department listed either by **telephone** or by **e-mail**.
- You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

You can download these instructions along with many other manuals, product videos and software on [www.lidl-service.com](http://www.lidl-service.com).

## Service

### GB Service Great Britain

Tel.: 0800 404 7657

E-Mail: [kompennass@lidl.co.uk](mailto:kompennass@lidl.co.uk)

### IE Service Ireland

Tel.: 1890 930 034

(0,08 EUR/Min., (peak))

(0,06 EUR/Min., (off peak))

E-Mail: [kompennass@lidl.ie](mailto:kompennass@lidl.ie)

IAN 312116

## Importer

Please note that the following address is not the service address. Please use the service address provided in the operating instructions.

KOMPENASS HANDELS GMBH

BURGSTRASSE 21

DE-44867 BOCHUM

GERMANY

[www.kompennass.com](http://www.kompennass.com)

## Ordering replacement parts

You can order the following replacement parts for product NKZ 3 A1:

► 8-pack of toothbrush heads (blue)

► 8-pack of toothbrush heads (pink)

Order the replacement parts via the service hotline (see section "Service") or simply visit our website at [www.kompernass.com](http://www.kompernass.com).

### NOTE

- Have the IAN number of the appliance ready (can be found on the cover of these operating instructions) when you place your order.

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