PORE CLEANSER / PORENREINIGER / NETTOYEUR DE PORES SPR 600 B1

(GB) (IE)

PORE CLEANSER

Operating instructions

(FR) (BE)

NETTOYEUR DE PORES

Mode d'emploi

(CZ)

ČISTIČ PÓRŮ

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(SK)

ČISTIČ PÓROV

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(DK)

PORERENSER

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(HU)

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(DE) (AT) (CH)

PORENREINIGER

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(NL) (BE)

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(PL)

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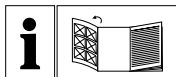
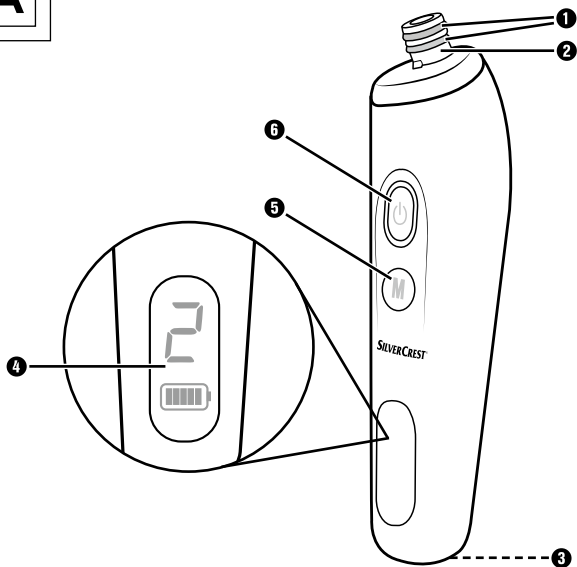
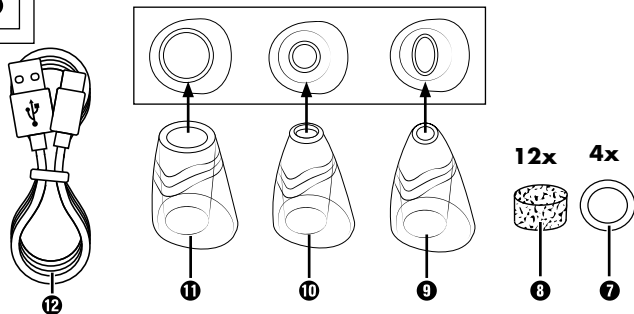
A**B**

Contents

Introduction.....	2
Proper use	2
Warnings and symbols used	2
Notes on trademarks	3
Safety information	4
Package contents	9
Appliance description.....	10
Charging the integrated battery.....	10
Before first use	11
Preparing the appliance.....	11
Checking compatibility.....	12
Handling and operation.....	12
Selecting an attachment.....	12
Fitting an attachment with filter.....	13
Intensity levels.....	14
Before use	14
Treatment of the face	14
Areas of application.....	16
After use.....	17
Cleaning and care.....	17
Storage	18
Ordering replacement parts	19
Disposal.....	19
Disposal of the appliance.....	20
Disposal of the packaging	20
Technical specifications.....	21
Kompernass Handels GmbH warranty	21
Service	23
Importer	23

Introduction

Congratulations on the purchase of your new appliance.

You have selected a high-quality product. The operating instructions are part of this product. They contain important information about safety, usage and disposal. Before using the product, please familiarise yourself with all operating and safety instructions. Use the product only as described and for the specified range of applications. Please also pass these operating instructions on to any future owner.

Proper use

This appliance uses vacuum technology for deep cleaning of facial pores and treating blackheads. It also removes dead skin cells and stimulates circulation. The appliance is not suitable for treating sensitive areas around the eyes or lips.

The appliance is intended solely for personal, private use. Do not use the appliance on animals. It is not intended for use in medical/therapeutic or commercial fields.

It is not intended for any other purpose nor for use beyond the scope described. If you are uncertain whether the appliance is appropriate for you, please consult a physician.

Warnings and symbols used

The following warnings and symbols are used in these operating instructions, on the packaging and on the appliance (where applicable):

DANGER! A warning with this symbol and the signal word "DANGER" indicates an imminently hazardous situation that will result in death or serious injury if not avoided.

WARNING! A warning with this symbol and the signal word "WARNING" indicates a potentially hazardous situation that could result in death or serious injury if not avoided.

	CAUTION! A warning with this symbol and the signal word "CAUTION" indicates a potentially hazardous situation that could result in a minor or moderate injury if not avoided.
	ATTENTION! A warning with this symbol and the signal word "ATTENTION" indicates a potential situation that could result in property damage if not avoided.
	A note identifies additional information that facilitates the use of the appliance.
	Read the instructions.
	DC current/voltage
	Removable power supply

Notes on trademarks

USB® is a registered trademark of USB Implementers Forum, Inc.

The SilverCrest trademark and the retail name are the property of their respective owner.

All other names and products may be trademarks or registered trademarks of their respective owner.

Safety information

DANGER!

- Connect the mains adapter only to correctly installed and earthed mains sockets. Ensure that the power supply corresponds with the information on the type plate of the appliance/mains adapter.
- To avoid potential risks, damaged mains cables should be replaced by the manufacturer, their customer service department or someone similarly qualified.
- Make sure that the charging cable does not become wet or damp. Route the cable in such a way that it cannot be caught or damaged.
- Keep the appliance and the charging cable away from hot surfaces.
- Do not perform any repairs to the appliance. All repairs must be performed by our customer service department or by a qualified specialist technician.
- Never open the appliance housing. There are no parts inside the appliance that require user maintenance. In addition, this will invalidate the warranty.
- Switch off the appliance after each use and before cleaning it.
- Never immerse the appliance in water or other liquids!

- ▶ Never reach for the appliance if it has fallen into water.
- ▶ Do not use the appliance while taking a bath or shower. Do not store or keep the appliance in places from where it can fall into the bath tub, sink or other containers filled with water or other liquids.
- ▶ Do not use the appliance if the power cable is damaged, if the mains adapter or accessories are defective, if the appliance is not working properly, or if it has sustained a fall.

⚠ WARNING!

- ▶ This appliance may be used by children aged 8 over and by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, provided that they are under supervision or have been told how to use of the appliance safely and are aware of the potential risks. Children must not play with the appliance. Cleaning and user maintenance tasks must not be carried out by children unless they are supervised.
- ▶ Children must be supervised to make sure that they do not play with the appliance or its parts.
- ▶ Never leave children unattended with the packaging material. Keep children away from the packaging and the appliance. There is a risk of suffocation!

- ▶ Do not use the appliance if you suffer from one or more of the following illnesses or complaints: heart disease, multiple sclerosis, diabetes, autoimmune illnesses, skin cancer, circulatory disorders, or skin conditions such as herpes, vascular disorders, facial erythrosis, rosacea, skin growths, cracked skin, irritated skin.
- ▶ Never insert the suction head or a connected attachment into your eyes!

⚠ WARNING!

- ▶ This product is equipped with or supplied with a lithium-ionen battery.
- ▶ This appliance contains a battery that cannot be replaced.
- ▶ The appliance has an integrated lithium-ionen battery. Incorrect handling can lead to fire, explosions, leakage of dangerous substances or other hazards!
- ▶ Do not throw the appliance into an open fire.
- ▶ Do not open or short-circuit the battery. This could cause the battery to overheat and explode.
There is a risk of fire!

- ▶ Never store the device when it is discharged. Prolonged storage in a discharged state can cause permanent damage to the battery. If the device is stored for a longer period of time, the charge level of the battery must be checked regularly. The optimum charge level is between 50% and 80%.

⚠ CAUTION!

- ▶ Never use the appliance on animals, during pregnancy, after cosmetic surgeries on the face, or on parts of the body that are exhibiting swelling, burns, sunburn, inflammation, eczema, wounds or cracked skin.
- ▶ Do not use the appliance if you are taking a blood thinner or medications containing steroids.
- ▶ Stop treatment with the appliance immediately if you experience any discomfort or pain. Excessive strain on the skin can lead to allergic reactions, bruises, swelling or broken capillaries. In case of doubt, consult your physician to determine if the treatment is suitable for you.
- ▶ Do not use the appliance on sensitive, movable parts of the eyes or on the lips.
- ▶ Do not use the appliance on nipples or in the genital area.

- ▶ Never use the appliance without the attachment.
- ▶ Uninterrupted application of the appliance may not last any longer than 5 minutes. The first treatments may not exceed 3 minutes. Do not treat the same area for too long or with too much force. Do not treat the same area several times in succession.
- ▶ Do not insert your fingers or any other objects into the opening for the attachment or into the suction head. Keep long hair away from the appliance.

! ATTENTION!

- ▶ Use the appliance only with the supplied accessories.
- ▶ When using a USB power adapter, the mains socket used must always be easily accessible so that the USB power adapter can be quickly removed from the socket in a hazardous situation. You should also observe the operating instructions for the USB power adapter.
- ▶ Only charge the appliance with a USB power adapter with protection class II that is approved for use with household appliances and has an output voltage of 5 V $\equiv$ and an output current of max. 2 A.
- ▶ Use only the supplied charging cable to charge the appliance.

- ▶ Keep the appliance away from heat sources (radiators, ovens, etc.) and do not expose it to extreme temperature fluctuations.
- ▶ There are no serviceable parts or parts requiring cleaning inside the appliance. Penetrating moisture can cause damage to the appliance.
- ▶ Do not use any aggressive cleaning liquids or solvents to clean the appliance. They could damage the surface.

Package contents

The appliance is supplied with the following components as standard:

- Pore cleanser
- USB charging cable (USB-A to USB-C)
- 3x Attachments
- 12x Filters
- 6x Locking rings
- Operating instructions

⚠ DANGER! Do not allow children to play with packaging materials. There is a risk of suffocation.

- 1) Remove all parts of the appliance and the operating instructions from the box.
- 2) Remove all packaging material.

i Note:

Check the package for completeness and signs of visible damage. If the delivery is incomplete or damage has occurred as a result of defective packaging or during transport, contact the customer service hotline (see **Service** section).

Appliance description

Figure A:

- ❶ 2x Locking rings (on the appliance)
- ❷ Suction head
- ❸ USB-C port
- ❹ LED display
- ❺ Intensity button 
- ❻ On/Off button 

Figure B:

- ❼ 4x Spare locking rings
- ❽ 12x Filters
- ❾ Oval attachment
- ❿ Small attachment
- ⓫ Large attachment
- ⓬ USB charging cable (USB-A to USB-C)

Charging the integrated battery

 **DANGER!** Risk of fatal electric shock! Only charge the appliance indoors in dry locations away from water. For example, do not charge the appliance next to or above a sink filled with water.

 Note:

Before using the appliance, the integrated battery must be completely charged.

The USB charging cable ❶ and the USB-C port ❸ are used only to charge the internal battery. They are not suitable for data transmission.

Use only the supplied USB charging cable ❶ to charge the integrated battery. Disconnect the USB charging cable ❶ when charging is complete.

- 1) Make sure that the appliance is switched off. Hold down the On/Off button  **6** for about 2 seconds to switch off the appliance.
- 2) Connect the USB-A plug of the USB charging cable **12** to a suitable mains adapter.
- 3) Connect the USB-C plug of the USB charging cable **12** to the USB-C port **3** of the appliance.
- 4) Plug the power adapter into a mains socket.
In the LED display **4**, the battery symbol  lights up, and the charging progress is indicated by flashing of the last bar.
When all bars light up continuously, the battery is charged.

Battery symbol	Charge level
Frame and first bar flashing.	Low battery. Recharge immediately.
Bars are partially lit continuously, last bar flashes.	Battery is charging. Bars indicate charging progress.
All bars light up continuously.	Battery is fully charged.

- 5) After charging, unplug the power adapter from the mains socket and the USB charging cable **12** from the appliance.

Before first use

Preparing the appliance

- 1) Charge the battery as described in the section **Charging the integrated battery**.
- 2) We recommend cleaning the appliance and the attachments **9** - **11** as described in the section **Cleaning and care** before using them for the first time. This will remove any residues left from the production process.

The appliance is now ready for use.

Checking compatibility

- We recommend first testing the appliance in a less sensitive area. Check how your skin reacts, for example, by testing it on the inside of your lower arm. First read the section **Handling and operation** so that you know how to use the appliance safely and avoid injuries.
- Clean the skin and make sure that no cream or moisturizer has been applied to the area.
- Only continue with the application on your face if no skin reactions appear on the test area after about 30 minutes.
- In case of doubt, consult your physician to determine if the treatment is suitable for you.

Handling and operation

Selecting an attachment

The appliance comes with 3 attachments ⑨ - ⑪ for treating the skin.

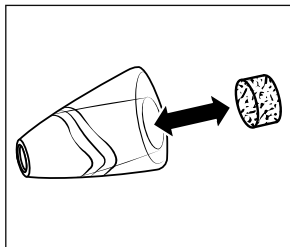
You can:

- Clean the pores of your face and treat blackheads.
- Remove dead skin cells.
- Stimulate blood circulation in the skin.

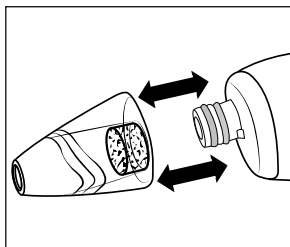
Attachment	Opening shape	Areas of application
Small attachment ⑩		This attachment treats a small area, making it suitable for hard-to-reach locations, such as on the nose. The attachment is also suitable for initial applications or for sensitive skin.

Attachment	Opening shape	Areas of application
Large attachment 11		This attachment treats a large area. The attachment can be applied to all areas of the face. Pay attention to the unique aspects of the eyes and mouth! This attachment is also suitable for use on blackheads.
Oval attachment 9		This attachment treats a somewhat large area and is particularly suitable for the sides of the nose and around the mouth. We recommend the intensity level "Soft".

Fitting an attachment with filter

- 1) If already fitted: remove the fitted attachment **9 - 11** from the suction head **2**.
- 2) To change the filter **8**, pull the current filter **8** out of the suction channel of the attachment **9 - 11**. Use tweezers, if necessary.
- 3) Insert a filter **8** into the suction channel of the selected attachment **9 - 11**.

- 4) Fit the attachment **9 - 11** onto the suction head **2** so that the small rectangular tab on the lower part of the suction head **2** engages precisely into the matching recess in the suction channel of the attachment **9 - 11**.
- 5) The filter **8** must be resting directly on the suction head **2** for the appliance to function properly. As soon as you start the appliance, the filter **8** is sucked in.

Intensity levels

You can choose between 3 intensity levels. The set level is shown on the LED display ④:

- "1": Light suction
- "2": Moderate suction
- "3": Strong suction

Before use

- Before beginning a treatment, clean your skin to remove any powder, cream or make-up. Ensure that your skin is completely dry.
- Avoid using a sauna or tanning bed on the day before and the day after a treatment. Protect your skin from the sun.
- For the first applications or for sensitive skin, use the small attachment ⑩, which treats only a small area at a time. Also select the intensity level "1" with the lowest suction power. Sensitive skin can react with swelling. Perform a compatibility test, if necessary: Before treating your face, check how your skin reacts at another location, such as the inside of your lower arm. Wait about 30 minutes for a possible skin reaction.
- Tip: Opening the pores before treatment can improve the result. Facial steamers or facial saunas are typically available for this purpose. Alternatively, a conventional facial steam bath can be done with hot water in a bowl and a towel over the head.

Treatment of the face

⚠ CAUTION! There is a risk of bruises, swelling and skin irritation.

Always begin at the lowest intensity level and gradually increase the intensity, if needed.

Take special care and use a low intensity level when treating sensitive areas or sensitive skin.

Do not use the appliance on the sensitive skin immediately around the eyes or on the lips.

Everyone's body reacts differently. If you find an application unpleasant, reduce the intensity level, switch immediately to a different location, reduce the pressure you exert on your skin with the appliance or stop the treatment completely.

- 1) Fit the selected attachment ⑨ - ⑪ with filter ⑧.
- 2) Start the appliance by pressing the On/Off button ⑥ for about 2 seconds.
Intensity level "2" lights up in the LED display ④.
- 3) Briefly press the Intensity button ⑤ to set the intensity level. The LED display ④ shows you the set intensity level.

⚠ CAUTION! There is a risk of bruises, swelling and skin irritation.

Uninterrupted application of the appliance may not last any longer than 5 minutes. The first treatments may not exceed 3 minutes. Do not treat the same area for too long or with too much force. Do not treat the same area several times in succession.

Do not use the appliance more often than twice per week.

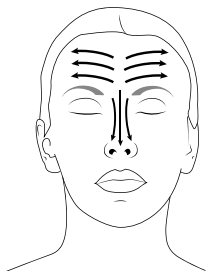
- 4) Place the opening of the attachment ⑨ - ⑪ against the skin of your face. Move the appliance as follows:
Press only hard enough for the attachment ⑨ - ⑪ to make full contact with the skin so that the appliance cannot draw in outside air. You will feel a light suction force on the area being treated. Never remain at a fixed location. Always keep the appliance in motion. Pull the appliance over each area of skin only a single time!
- 5) Treat the areas of the face as described in the section **Areas of application**.
- 6) Press the On/Off button ⑥ for about 2 seconds to switch off the appliance.

Areas of application

i Note:

Take care that your skin is firm so that the attachment **9 - 11** can make a seal against the skin and not draw in outside air.

If necessary, use one hand to pull the treated area and tighten it while moving the appliance over the area with your other hand.

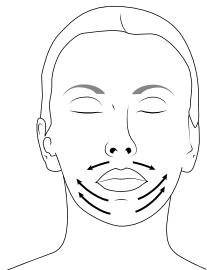
Forehead: Place the appliance in the middle of your forehead and pull it outward to just before the hairline.

Nose: Place the appliance on the upper part of your nose and pull it down the sides and middle of the nose.

Cheeks: Place the appliance next to your nose and pull it over your cheeks.

Temples: Place the appliance next to your eye. Maintain a distance of about 1 – 2 cm from the eye. Then pull the appliance towards the hairline.

Mouth and chin: Place the appliance above the mouth and pull it to the side towards your cheek. Do not touch the lips directly. Place the appliance on your chin and pull it up towards the cheeks. Stop before you reach your cheek.

After use

- Clean your face with water. Do not use any irritating care products. Apply soothing skincare products or a moisturizing cream.
- You may observe some reddening of the skin after the treatment. However, this is temporary and should go away within one or two days, at most. If you experience any issues, please consult your physician.
- Avoid using a sauna or tanning bed on the day before and the day after a treatment. Protect your skin from the sun.

Cleaning and care

⚠ DANGER! Risk of fatal injury from electrocution! Before every cleaning, disconnect the USB charging cable ⑫ from the appliance and pull the mains adapter out of the mains socket.

Never immerse the appliance, the USB power cable ⑫ or the mains adapter in water or other liquids!

⚠ ATTENTION! Do not use abrasive or caustic cleaning agents to clean the appliance. They can corrode the surface and damage the appliance beyond repair.

Never open the housing of the appliance. There are no serviceable parts or parts requiring cleaning inside the appliance. Penetrating moisture can cause damage to the appliance.

- Make sure that the appliance is switched off and the USB charging cable 12 is not connected to the appliance.
- Remove the fitted attachment 9 - 11 and take out the used filter 8. You can use the filters 8 multiple times. After use, wash the filter in lukewarm water and allow it to dry before the next use. If the filter is still dirty or defective, use a new filter 8.
- Clean the attachment 9 - 11 with water and washing-up liquid. If necessary, use a moistened cotton swab to clean the inside of the attachments 9 - 11. We recommend additionally applying a suitable disinfectant. Allow the attachments 9 - 11 to dry before the next use.

i Note:

If the locking rings 1 break during cleaning, replace them with the provided spare locking rings 7. If necessary, you can order more filters 8, locking rings 1 and attachments 9 - 11 (see section **Ordering replacement parts**).

- In case of soiling, pull the locking rings 1 off the suction head 2. Use tweezers, if necessary. Clean the locking rings 1 in lukewarm water with some washing-up liquid.
- Wipe the area around the suction head 2 with a slightly damp cloth. Then wipe it again with a dry cloth. Clean the other surfaces of the appliance only with a dry cloth.
- Place the locking rings 1 back on the suction head 2 and slide them into the provided grooves.

Storage

- Store the appliance and its accessories in a dry and dust-free location.
- To ensure a long battery life, charge the battery at regular intervals.
- Never store the device when it is discharged. Prolonged storage in a discharged state can cause permanent damage to the battery. If the device is stored for a longer period of time, the charge level of the battery must be checked regularly. The optimum charge level is between 50% and 80%.

Ordering replacement parts

You can order replacement parts for this product on the Internet at www.kompernass.com.

Scan the QR code with your smartphone/tablet. You can use this QR code to go directly to our website to view and order the available spare parts.

Note:

If you have problems with your online order, you can contact our service centre by phone or e-mail.

Please always quote the article number (IAN) 472246_2407 in your order.

Please note that online ordering of replacement parts is not possible for all countries.

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FR

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FR

The product, its packaging and the operating instructions are recyclable. They are subject to an extended manufacturer responsibility and will be collected separately.

Disposal of the appliance

The adjacent symbol of a crossed-out dustbin means that this appliance is subject to Directive 2012/19/EU. This directive states that this appliance may not be disposed of in the normal household waste at the end of its useful life, but must be taken to specially set-up collection locations, recycling depots or disposal companies.

The disposal is free of charge for the user. Protect the environment and dispose of this appliance properly.

If your old appliance has stored any personal data, you are responsible for deleting it yourself before returning it.

If it is possible to do so without destroying the old appliance, remove the old batteries or rechargeable batteries before returning the appliance for disposal and take them to a separate collection point. In the case of permanently installed rechargeable batteries, you must indicate during disposal that the appliance contains a battery.

The integrated rechargeable battery in this appliance cannot be removed for disposal.

Your local community or municipal authorities can provide information on how to dispose of the worn-out product.

Disposal of the packaging

The packaging materials have been selected for their environmental friendliness and ease of disposal and are therefore recyclable. Dispose of packaging materials that are no longer needed in accordance with applicable local regulations.

Dispose of the packaging in an environmentally friendly manner. Note the labelling on the packaging and separate the packaging material components for disposal, if necessary. The packaging material is labelled with abbreviations (a) and numbers (b) with the following meanings: 1-7: plastics, 20-22: paper and cardboard, 80-98: composites.

Applies only to Spain:

The packaging contains paper and/ or cardboard components.

The packaging contains plastic and/ or metal components.

Technical specifications

Input voltage	5.0 V ===
Input current	max. 2 A
Integrated battery (Li-Ionen)	600 mAh / 3.7 V === / 2.22 Wh
Protection class	III /  (protection through low voltage)

Kompernass Handels GmbH warranty

Dear Customer,

This appliance has a 3-year warranty valid from the date of purchase. If this product has any faults, you, the buyer, have certain statutory rights. Your statutory rights are not restricted in any way by the warranty described below.

Warranty conditions

The warranty period starts on the date of purchase. Please keep your receipt in a safe place. This will be required as proof of purchase.

If any material or manufacturing fault occurs within three years of the date of purchase of the product, we will either repair or replace the product for you or refund the purchase price (at our discretion). This warranty service requires that you present the defective appliance and the proof of purchase (receipt) within the three-year warranty period, along with a brief written description of the fault and of when it occurred.

If the defect is covered by the warranty, your product will either be repaired or replaced by us. The repair or replacement of a product does not signify the beginning of a new warranty period.

Warranty period and statutory claims for defects

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking. Repairs carried out after expiry of the warranty period shall be subject to a fee.

Scope of the warranty

This appliance has been manufactured in accordance with strict quality guidelines and inspected meticulously prior to delivery.

The warranty covers material faults or production faults. The warranty does not extend to product parts subject to normal wear and tear or to fragile parts which could be considered as consumable parts such as switches or parts made of glass.

The warranty does not apply if the product has been damaged, improperly used or improperly maintained. The directions in the operating instructions for the product regarding proper use of the product are to be strictly followed. Uses and actions that are discouraged in the operating instructions or which are warned against must be avoided.

This product is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications/repairs which have not been carried out by one of our authorised Service centres.

Warranty claim procedure

To ensure quick processing of your case, please observe the following instructions:

- Please have the till receipt and the item number (IAN) 472246_2407 available as proof of purchase.
- You will find the item number on the type plate on the product, an engraving on the product, on the front page of the operating instructions (below left) or on the sticker on the rear or bottom of the product.

- If functional or other defects occur, please contact the service department listed either by **telephone** or by **e-mail**.
- You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

You can download these instructions along with many other manuals, product videos and installation software at www.lidl-service.com.

This QR code will take you directly to the Lidl service page (www.lidl-service.com) where you can open your operating instructions by entering the item number (IAN) 472246_2407.

Service

GB Service Great Britain

Tel.: 0800 404 7657

E-Mail: kompernass@lidl.co.uk

IE Service Ireland

Tel.: 1800 101010

E-Mail: kompernass@lidl.ie

IAN 472246_2407

Importer

Please note that the following address is not the service address. Please use the service address provided in the operating instructions.

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