DIAGNOSEWAAGE SPWD 180 I1

GB IE NI

DIAGNOSTIC SCALES

Short manual

FR BE

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Guide succinct

NL BE

DIAGNOSE WEEGSCHAAL

Beknopte handleiding

DE AT

DIAGNOSEWAAGE

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

TARGA GMBH
Coesterweg 45
59494 Soest
GERMANY

Last Information Update - Version des informations
Datum nieuwste versie - Stand der Informationen:
04 / 2022 - Ident.-No.: SPWD 180 I1 042022-1

IAN 498992_2204

IAN 498992_2204

English..... 1
 Français.....12
 Nederlands.....24
 Deutsch.....36

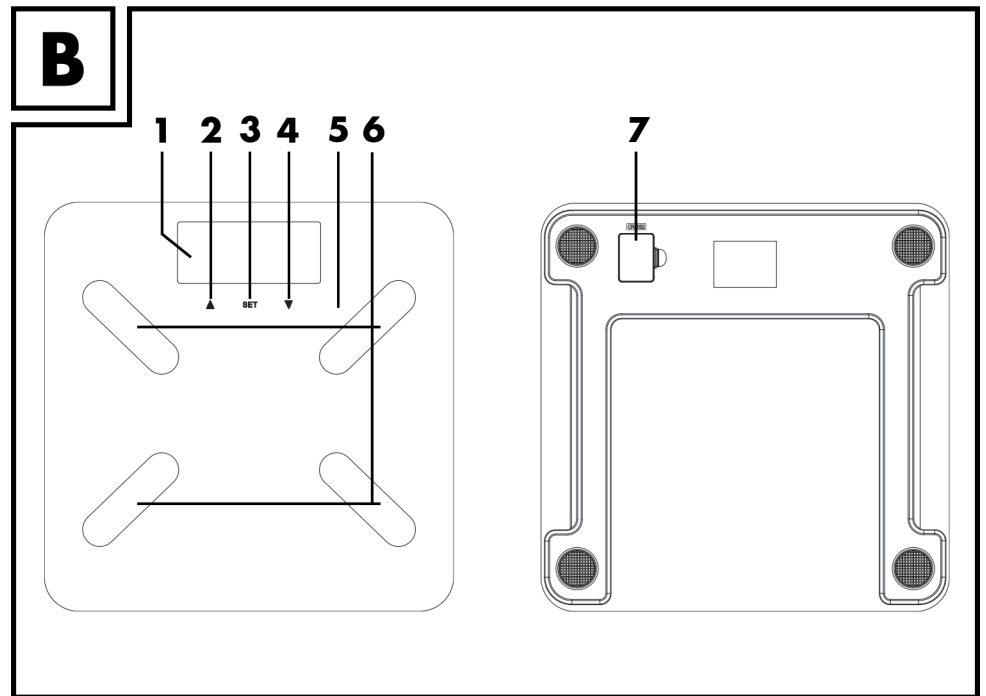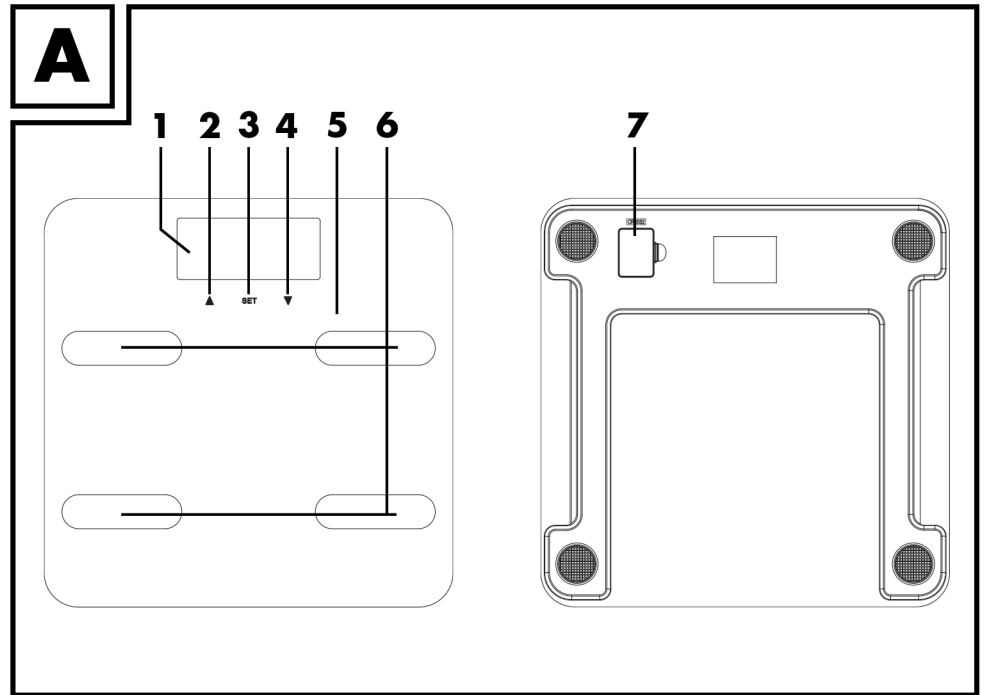This document is an abbreviated print version of the complete operating instructions. Scanning the QR code will take you directly to the Lidl service page (www.lidl-service.com) and by entering the first 6 digits of the item number (IAN) 498992_2204 you can view and download the full operating instructions.

WARNING! Observe the complete operating instructions and the safety instructions to avoid personal injury and damage to property. The Quick Start Guide is an integral part of this product. Familiarise yourself with all operating and safety instructions before using the product. Keep the Quick Start Guide in a safe place and hand over all documents when passing the product on to third parties.

Intended use

These diagnostic scales are used to determine a person's weight. Furthermore, it is possible to measure body fat, water, muscle and bone percentages and to determine your recommended daily calorie intake. These diagnostic scales may not be used outdoors or in tropical climates. These diagnostic scales have not been designed for corporate or commercial use. Use the diagnostic scales in domestic environments for private applications only. Any other use does not correspond to the intended use. These diagnostic scales fulfill all relevant norms and standards concerning CE conformity. Any modifications to the diagnostic scales other than recommended changes by the manufacturer may result in these norms no longer being met. The manufacturer accepts no liability for any damage or malfunctions resulting from this. Only use accessories recommended by the manufacturer.

Observe the regulations and laws in the country of use.


Package contents

- SilverCrest SPWD 180 I1 diagnostic scales
- 2x CR2032 lithium batteries
- These operating instructions

Overview

- 1 Display
- 2 ▲ button
- 3 SET button
- 4 ▼ button
- 5 Tread
- 6 Electrodes
- 7 Battery compartment

Display

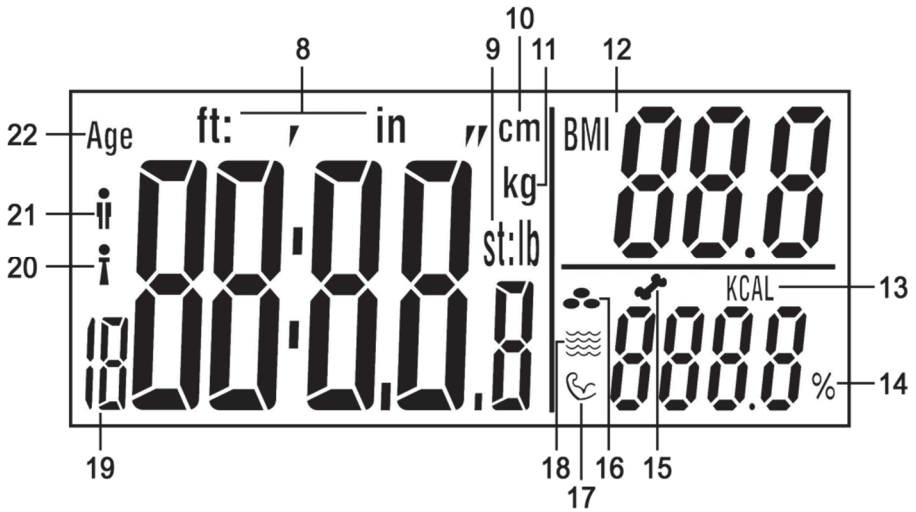- 8 Body height in feet:inches
- 9 **st:lb** Unit of weight in stones:pounds
lb Unit of weight in pounds
- 10 Body height in cm
- 11 Unit of weight in kilograms
- 12 BMI (Body Mass Index)
- 13 Recommended calorie intake per day
- 14 Percentage (of body fat, water, muscle, bones)
- 15 Bone percentage
- 16 Body fat percentage
- 17 Muscle percentage
- 18 Water percentage
- 19 User profile memory
- 20 Female
- 21 Male
- 22 Age

Technical specifications

Operating voltage:	6 V 
Battery:	2 Lithium battery, 3 V, CR2032
Units of weight:	kg / st / lb (adjustable)
Maximum weight:	180 kg / 396 lb / 28 st
Graduation:	0.1 kg / 0.2 lb
Measuring tolerance:	+/- (1% +0.1 kg/ 0.2 lb) for 30 kg to 180 kg
Auto Power Off function:	yes
User profile memories:	12
Gender:	female / male
Age range for diagnostic functions:	8 - 80 years
Levels of activity:	6 (none, little, normal, moderate, intense, athletic)
Supported body heights:	100 - 220 cm
Dimensions (W x H x D):	approx. 30 x 2,2 x 30 cm
Weight:	approx. 1.6 kg incl. batteries
Operating temperature:	+10°C to +40°C

The technical data and design may be changed without prior notice.

Safety instructions

Before you use the diagnostic scales for the first time, please read the following notes and heed all warnings, even if you are familiar with handling electronic devices. Keep these operating instructions in a safe place for future reference. If you sell the diagnostic scales or give them away, pass on the operating instructions at the same time.

WARNING! This signal word indicates a danger with a moderate risk, which can lead to death or major injuries if not avoided.

This symbol indicates that one must not step on the diagnostic scales with wet or damp feet. The same applies if the tread is wet or damp. You could slip and injure yourself.

This symbol indicates that the diagnostic scales cannot bear weights over 180 kg or 396 lb. Otherwise, the diagnostic scales could be damaged and there is a risk of injury.

This symbol denotes further information on the topic.

⚠ This diagnostic scales may be used by children of 8 years or over or by persons with physical, sensory or mental impairments or those with no knowledge or experience providing they are supervised or they have been given instructions on the proper use of the diagnostic scales and they understand the associated risks. Do not allow children to play with the diagnostic scales. The product must not be cleaned or maintained by unsupervised children. Keep the packaging materials out of the reach of children. There is a risk of suffocation.

⚠ Do not use the diagnostic scales under any circumstances if you wear a pacemaker. Never use the diagnostic scales in the vicinity of anyone who wears a pacemaker or other medical devices.

⚠ The diagnostic scales are not suitable for pregnant women.

- ⚠ Do not step on the diagnostic scales with wet or damp feet or when the tread is damp. You could slip and injure yourself.
- ⚠ Do not overload the diagnostic scales with more than 180 kg or 396 lb. Otherwise, the diagnostic scales could be damaged and there is a risk of injury.
- ⚠ Do not jump on the tread of the diagnostic scales. Otherwise, they could be damaged and there is a risk of injury.
- ⚠ If smoke, unusual noises or smells are noticed, remove the batteries from the diagnostic scales immediately. If this occurs, the diagnostic scales should not be used before it they have been inspected by authorised service personnel. Never inhale smoke from a possible device fire. If you do inadvertently inhale smoke, seek medical attention. Smoke inhalation can be damaging to your health.
- ⚠ The diagnostic scales must not be damaged in order to avoid any further risks.
- ⚠ Make sure that no fire hazards (e.g. burning candles) are placed on or near the diagnostic scales. Fire hazard!
- ⚠ Do not expose the diagnostic scales to any direct heat sources (e.g. heaters), direct sunlight or artificial light. The product must not be exposed to spray or dripping water or abrasive liquids. Do not use the diagnostic scales near water. The diagnostic scales must especially never be immersed (do not place any containers containing liquids such as drinks, vases, etc. on the diagnostic scales). Take care that the diagnostic

scales are not subject to excessive impacts or vibrations. Do not introduce any foreign bodies into the device. as these could damage the diagnostic scales.

 Do not short-circuit the batteries. Remove the batteries from the device when you are not going to use it for a long period of time. If the batteries have leaked, remove them to prevent damage to the diagnostic scales. You can wipe off any leaked battery fluid by using a dry, absorbing cloth. Always wear protective gloves! Avoid any skin contact. If the battery fluid comes into contact with the skin or eyes, rinse immediately with plenty of water and seek medical aid.

 Non-rechargeable batteries must not be recharged.

 Different types of batteries or new and old batteries must not be used together.


 Do not short-circuit any battery contacts.

 Batteries must be inserted respecting the correct polarity.

 Empty batteries must be removed from the device and must be disposed of properly.

Before you start

Activating the batteries

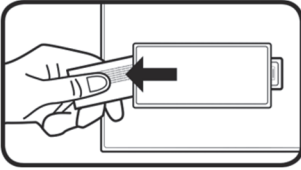Remove the strip from the battery compartment (7) on the bottom of the diagnostic scales. The diagnostic scales are now ready to use.


Replacing the batteries

Replace the batteries when they are empty. To replace the lithium batteries CR2032, remove the cover of the battery compartment (7) at the rear of the diagnostic scales by inserting a finger into the recess of the battery compartment (7). Then lift the battery compartment (7) cover upwards and remove it.

To remove the batteries, push them down a bit (see the direction of the arrow in the illustration). You can then remove the batteries from the holder.

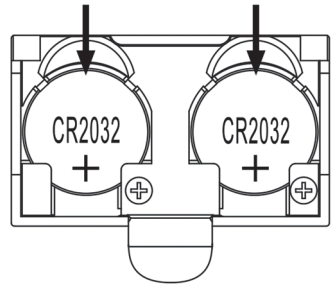Then insert two new CR2032 lithium batteries. Always pay attention to the polarity. This is indicated in the battery compartment (7), and the + symbol on the CR2032 lithium batteries must be visible, as shown in the adjoining illustration.

Next insert the battery compartment (7) cover back to its original position. The diagnostic scales are now ready to use.

Environmental regulations and disposal information

	<p>Devices marked with this symbol are subject to the European Directive 2012/19/EU. All electrical and electronic devices must be disposed of separately from household waste at official disposal centres. Avoid damage to the environment and risks to your personal health by disposing of the device properly. For further information about proper disposal, contact your local government, disposal bodies or the shop where you bought the device.</p>
	<p>The symbol of the crossed-out wheeled bin on conventional and rechargeable batteries indicates that they cannot be disposed of in normal domestic waste but must be disposed of separately.</p> <p>Where batteries contain toxic materials, the chemical symbol of the toxic material is shown underneath the symbol, with meanings as follows:</p> <ul style="list-style-type: none"> - Pb: Battery contains lead - Cd: Battery contains cadmium - Hg: Battery contains mercury <p>You are obliged by law to return used batteries. Old batteries can contain toxic materials which can be damaging to health or the environment if not properly stored or disposed of. Batteries also contain important raw materials, such as iron, zinc, manganese and nickel, which can be reused.</p> <p>After use, you can return the batteries to us or to a local collection point (e.g. in retail outlets or a local collection point) free of charge. Please note that batteries must be disposed of fully discharged at appropriate collection points for old batteries. If disposing of batteries which are not fully discharged, take precautions to prevent short circuits.</p>
	<p>Dispose of all packaging in an environmentally friendly manner. Cardboard packaging can be taken to paper recycling centres or public collection points for recycling. Any film or plastic contained in the packaging should be taken to your public collection points for disposal.</p>
 ES/PT	
<p><u>Only relevant for France:</u></p>  <p>“Sorting made easy”</p> <p>The product is recyclable, subject to advanced manufacturer responsibility and is sorted and collected separately.</p>	

	<p>Please note the markings on the packaging material when disposing of it, it is labelled with abbreviations (a) and numbers (b), the meanings of which are as follows:</p> <p>1-7: plastic / 20-22: paper and cardboard / 80-98: composite materials.</p>
	

Symbol	Material	Contained in the following packaging elements for this product
	Polyethylene terephthalate	Battery isolation strip
	High density polyethylene	–
	Polyvinyl chloride	–
	Low density polyethylene	–
	Polypropylene	–
	Polystyrene	–
	Other plastics	–
	Corrugated cardboard	Sales packaging and inner box
	Other cardboard	–
	Paper	Tissue paper to protect the surface
	Paper/cardboard/plastic	–
	Paper/cardboard/plastic/ aluminium	–

Conformity notes

The product complies with the requirements of the applicable European and national directives. Evidence of conformity has been provided. The manufacturer has the relevant declarations and documentation.

The product complies with the requirements of the applicable national directives of the Republic of Serbia.

The product complies with the requirements of the applicable national directives for Great Britain.

The complete EU Declaration of Conformity is available for download from this link:

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