
ELEKTRISCHE KINDER-ZAHNBÜRSTE NZAODK 600 A1
KIDS' ELECTRIC TOOTHBRUSH NZAODK 600 A1
BROSSE À DENTS ÉLECTRIQUE ENFANT NZAODK 600 A1

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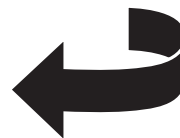DEUTSCHLAND

Contents

1. Overview	17
2. Intended purpose	18
3. Safety information	18
4. Items supplied	20
5. Setting up the battery charger	20
6. Charging the battery	21
7. Toothbrush heads	21
7.1 Selecting a toothbrush head	21
7.2 Replacing the toothbrush head	21
8. Operation	22
8.1 Modes	22
8.2 Brushing teeth	22
8.3 This is how your child brushes properly	23
9. Transport when travelling	23
9.1 Transport lock	23
9.2 Protective cap for the brush head	23
10. Cleaning / Hygiene requirements	24
11. Storage	24
12. Disposal	24
13. Trouble-shooting	25
14. Ordering toothbrush heads	25
15. Technical specifications	25
16. Warranty of the HOYER Handel GmbH	26

1. Overview

- 1 Protective cap for the toothbrush head
- 2 Large brush head
- 3 Small brush head
- 4 Opening
- 5 Toothbrush heads (replaceable)
- 6 Marking rings
- 7 Shaft of the toothbrush
- 8 Toothbrush
- 9 Control elements
- 10 Peg of the battery charger
- 11 Battery charger
- 12 Mains power cable and mains plug
- 13 Ridge (inside the protective cap)

Figure A: control elements

- 14  On/off switch + transport lock for on and off switch
- 15 Charging indicator light
- 16 Left LED: normal mode
- 17 Middle LED: sensitive mode
- 18 Right LED: massage mode

NOTE: the device is available in two different colours. This has no effect on the function or operation of the device.

Thank you for your trust!

Congratulations on the purchase of your new children's toothbrush.

For safe handling of the device and in order to get to know the entire range of features:

- Thoroughly read these user instructions prior to initial use.
- Above all, observe the safety information!
- The device should only be used as described in this copy of the user instructions.
- Keep this copy of the user instructions for reference.
- If you pass the device on to someone else, please include this copy of the user instructions. The user instructions are a part of the product.

We wish you a lot of joy with your new children's toothbrush!

All brand and manufacturer names are the property of their respective owners. They are specified solely in order to simplify the explanation of the points illustrated.

2. Intended purpose

The children's toothbrush is intended solely for tooth cleaning and oral hygiene for humans.

The device is designed for private, domestic use. This device must not be used for commercial purposes.

3. Safety information

Warnings

If necessary, the following warnings are used in this copy of the user instructions:

DANGER! High risk: failure to observe this warning may result in injury to life and limb.

WARNING! Moderate risk: failure to observe this warning may result in injury or serious material damage.

CAUTION: low risk: failure to observe this warning may result in minor injury or material damage.

NOTE: circumstances and specifics that must be observed when handling the device.

Instructions for safe operation

- ⦿ This device can be used by children and people with restricted physical, sensory or intellectual abilities or people without adequate experience and/or understanding, if they are supervised or instructed on how to use this device and if they are aware of the resulting hazards. Cleaning and user maintenance must not be performed by children unless they are supervised.
- ⦿ Children must not play with the device.
-  ⦿ The battery charger is connected to the mains voltage, and must therefore not be used in the shower, in the bath, or over washbasins filled with water.
- ⦿ The mains power cable of this device cannot be replaced. If the cable should become damaged, the device itself must be disposed of.

DANGER as small parts may be swallowed

- ⦿ The marking rings on the toothbrush heads are small parts that can be swallowed. Keep away from children under 36 months of age.

DANGER for children

- ⦿ Packing materials are not children's toys. Children should not be allowed to play with the plastic bags: there is a risk of suffocation.

DANGER to and from pets and livestock

- ⦿ Electrical devices can represent a danger to pets and livestock. In addition, animals can also cause damage to the device. For this reason you should keep animals away from electrical devices at all times.

DANGER! Risk of electric shock due to moisture

- ⦿ The toothbrush, the battery charger, the mains power cable and the mains plug must not be immersed in water or any other liquids.

- ⦿ Protect the device from moisture, water drops and water splashes: failure to observe this instruction may result in an electrical shock hazard.
- ⦿ Never touch the device with wet hands.
- ⦿ Make sure that the surface under the device is dry.
- ⦿ If the battery charger has fallen into water, pull out the mains plug immediately, and only then remove the battery charger from the water. In this case do not use the battery charger again, but have it checked by a specialist workshop.

DANGER! Risk of electric shock

- ⦿ The installation of a residual current device with a nominal trip current not exceeding 30 mA in the building wiring offers additional protection. Contact a skilled electrician in this regard.
- ⦿ Do not start using the device if the device, the battery charger or the mains power cable shows visible damage, or if the device has previously been dropped.
- ⦿ Only connect the mains plug to a properly installed and easily accessible wall

socket whose voltage corresponds to the specifications on the rating plate. The wall socket must still be easily accessible after the device is plugged in.

- ⊙ The hand unit of the toothbrush must only be charged using the battery charger supplied.
- ⊙ Ensure that the mains power cable cannot be caught or crushed, and that it cannot be damaged by sharp edges or hot points.
- ⊙ To fully disconnect the battery charger from the mains power supply, pull out the mains plug.
- ⊙ When removing the mains plug from the wall socket, always pull on the plug and never the cable.
- ⊙ Pull the mains plug out of the wall socket, ...
 - ... if there is a fault,
 - ... when you are not using the device,
 - ... before you clean the device and
 - ... during thunderstorms.
- ⊙ In order to avoid any hazards, do not make modifications to the device. Repairs must be carried out by a specialist workshop or the Service Centre.
- ⊙ When using an extension cable, this must correspond to the applicable safety regulations.

DANGER from batteries

- ⊙ Protect the battery against mechanical damages. Fire hazard!
- ⊙ Do not expose the device to direct sunshine or heat. The ambient temperature should not be lower than 10 °C or higher than 40 °C.
- ⊙ The battery may only be charged with the original accessory (power unit).

CAUTION

- ⊙ Users could possibly experience mild bleeding of the gums during the first days of electric toothbrush usage. Should the bleeding continue over several days, a dentist should be consulted.
- ⊙ Exert only light pressure when brushing to prevent damage to the gums and to increase the service life of the bristles.

4. Items supplied

- 1 toothbrush 8
- 2 + 2 toothbrush heads 5
(2x large brush head 2,
2x small brush head 3)
- 2 protective caps 1
- 1 battery charger 11
- 1 copy of the user instructions

- Remove all packing materials.
- Check that the device is undamaged.

5. Setting up the battery charger

DANGER! Risk of electric shock due to moisture!

- ⊙ The battery charger 11 is connected to the mains voltage, and must therefore not be used in the shower, in the bath, or over washbasins filled with water.

- Never set up the battery charger 11 anywhere it could fall into water.

6. Charging the battery

NOTES:

- Fully charge the battery before using it for the first time.
 - It takes approx. 16 hours to fully charge the battery.
 - Place the toothbrush 8 on the peg 10 of the battery charger 11 after use. This ensures that the toothbrush 8 always works at full power. The battery is not damaged by regular charging.
-

1. Plug the mains plug 12 of the battery charger 11 into a suitable wall socket. The wall socket must still be easily accessible after connection.
2. Place the switched-off toothbrush 8 onto the peg of the battery charger 10. The charging indicator light 15 is lit red throughout the entire charging process. When the optimal charging state has been reached, the charging indicator light 15 is lit green.

Battery signals: charging indicator light 15

Toothbrush 8 on the battery charger 11:

- Charging indicator light 15 is lit green: the battery is fully charged.
- Charging indicator light 15 is lit red: the battery is being charged.

Toothbrush 8 switched on:

- The toothbrush 8 switches itself off and the charging indicator light 15 flashes red: the battery is weak. Charge the battery.

7. Toothbrush heads

7.1 Selecting a toothbrush head

The items supplied include 4 toothbrush heads 5:

- 2 with small brush heads 3 (recommended for children aged 4 and above), and
 - 2 with large brush heads 2 (recommended for children aged 7 and above).
- Check which size is suitable for the mouth of your child.

7.2 Replacing the toothbrush head

- Take a toothbrush head 5 and slide it onto the shaft of the toothbrush 7.
- In order to remove, carefully pull the toothbrush head 5 off upwards.

NOTES:

- The various marking rings 6 allow you to assign the toothbrush heads 5 to different family members.
 - We recommend that you replace the toothbrush heads 5 every 3 months or after infections (see "Ordering toothbrush heads" on page 25).
-

8. Operation

CAUTION:

- ⊙ Users could possibly experience mild bleeding of the gums during the first days of electric toothbrush usage. Should the bleeding continue over several days, a dentist should be consulted.
 - ⊙ If your child has had dental surgery in the last two months, consult a dentist before using the toothbrush.
 - ⊙ Exert only light pressure when brushing to prevent damage to the gums and to increase the service life of the bristles.
-

8.1 Modes

The toothbrush 8 has 3 different modes. Depending on the desired mode, press the on/off switch  14 1, 2 or 3 times while switched off:

- 1x: normal mode
 - The left LED 16 is lit.
 - For cleaning teeth with high pulsation
 - 2x: sensitive mode
 - The middle LED 17 is lit.
 - Gentle brushing for cleaning teeth with sensitive gums
-

NOTE: if your child is not used to using an electric toothbrush and has comparatively sensitive gums, your child can use this mode for the first 1 to 2 weeks and then switch to normal mode later.

- 3x: massage mode
 - The right LED 18 is lit.
 - Pulsing brushing to massage your gums
-

8.2 Brushing teeth

1. Remove the toothbrush 8 from the battery charger 11.
 2. Moisten the bristles, apply some toothpaste and switch on the toothbrush 8 only when the brush head 2/3 is already in your mouth.
 3. Depending on the desired mode, press the on/off switch  14 1, 2 or 3 times:
 - left LED 16: normal mode
 - middle LED 17: sensitive mode
 - right LED 18: massage mode
 4. Slowly guide the brush head 2/3 from tooth to tooth without exerting too much pressure. Carefully clean all tooth surfaces, the inner surfaces of the teeth and the tooth edges.
 5. The recommended brushing time is at least 2 minutes, divided into intervals of 30 seconds each (see “This is how your child brushes properly” on page 23). The toothbrush indicates these intervals by beeping and short interruptions of the vibration. After 2 minutes, a tune will play to indicate that the recommended brushing time has been reached.
 6. Press the on/off switch  14 repeatedly until the toothbrush 8 is switched off. Remove the brush head 2/3 from your mouth only when the toothbrush 8 has been switched off.
-

NOTE: we recommend that you replace the toothbrush heads 5 every 3 months or after infections (see “Ordering toothbrush heads” on page 25).

8.3 This is how your child brushes properly

NOTE: the toothbrush is also suitable for cleaning orthodontic braces. In this case, however, the bristles wear out more quickly.

- The recommended brushing time is 2 minutes. When the time set has elapsed, a tune plays.
- Figure B: the brushing time is divided into 4 identical intervals. (For a brushing time of 2 minutes, this is 30 seconds each per mouth section: bottom right, bottom left, top right, top left.)
- The toothbrush indicates these intervals by beeping and short interruptions of the vibration. After 2 minutes, a tune will play to indicate that the recommended brushing time has been reached.
- It is recommended that children under the age of 8 should be supervised for the first 1 to 2 weeks while brushing their teeth with the electric toothbrush.

How to explain it to your child

- Figure B: explain to your child the 4 areas in the mouth (bottom right, bottom left, top right, top left).
- Tell your child that each time the toothbrush beeps and the vibration is briefly interrupted, they should brush another area of the mouth.
- Tell your child that they should brush until the tune plays.

9. Transport when travelling

9.1 Transport lock

To secure the toothbrush against accidentally being switched on, activate the transport lock.

- Switching on the transport lock:
Press the on/off switch  14 until the left LED 16 flashes once.
- If you try to switch on the toothbrush with the transport lock activated, the left LED 16 will flash.
- Switching off the transport lock:
Press the on/off switch  14 until the toothbrush is switched on. Change the mode if required by pressing the on/off switch  14 again.

9.2 Protective cap for the brush head

- Figure C: for transport, use gentle pressure to place the protective cap 1 over the toothbrush head 5. The small ridge 13 will click into the opening 4.
- The brush head 2 / 3 should be as dry as possible. If the brush head 2 / 3 has to be transported when it is still damp, remove the protective cap 1 as quickly as possible after you have transported the toothbrush and allow the brush head 2 / 3 to dry.

10. Cleaning / Hygiene requirements

DANGER! Risk of electric shock!

- ⦿ Before cleaning the battery charger 11, pull the mains plug 12 out from the mains socket.
- ⦿ Do not immerse the toothbrush 8 and battery charger 11 in water!
The battery charger 11 must be completely dry before you reconnect it to the wall socket.

- To ensure the hygienic functionality of the toothbrush heads 5 they must be rinsed thoroughly after every use and allowed to dry.
- Remove the toothbrush head 5 once per week and clean it in hot water.
- Rinse the toothbrush 8 as well and wipe it dry.
- Wipe the battery charger 11 with a damp cloth if necessary and rub it dry.

11. Storage

- Disconnect the device from the mains and store it in a dry location.

12. Disposal

This product is subject to the provisions of European Directive 2012/19/EC. The symbol showing a wheelie bin crossed through indicates that the product requires separate refuse collection in the European Union. This applies to the product and all accessories marked with this symbol. Products identified with this symbol may not be discarded with normal household waste, but must be taken to a collection point for recycling electric and electronic appliances.

This recycling symbol marks e.g. an object or parts of the material as valuable for recycling. Recycling helps to reduce the consumption of raw materials and protect the environment.

The battery integrated into this device must not be disposed of in household waste. The device must be disposed of properly with the battery installed.

Packaging

When disposing of the packaging, make sure you comply with the environmental regulations applicable in your country.

13. Trouble-shooting

If your device fails to function as required, please try this checklist first. Perhaps there is only a minor problem, and you can solve it yourself.

DANGER! Risk of electric shock!

- ⊙ Do not attempt to repair the device yourself under any circumstances.

Fault	Possible causes / Action
The device cannot be switched on.	• Is the battery empty? • Is the transport lock activated? Press the on/off switch  14 until the toothbrush 8 is switched on.
Charging indicator light 15 is not lit.	• The battery charger 11 is not connected to the mains power.

14. Ordering toothbrush heads

You can order additional toothbrush heads 5 on the internet:

shop.hoyerhandel.com

1. Scan the QR code with your smart-phone/tablet.
2. With the QR code you will access a website where you can submit your re-order.

15. Technical specifications

Model:	NZAODK 600 A1
Mains voltage:	100 – 240 V ~ 50/60 Hz
Protection class:	II 
Power rating:	1.4 W Standby: 0.4 W
Battery:	1x 3.7 V Li ion 600 mAh 2.22 Wh
Charging time:	16 hours
Operating time:	approx. 90 min
Dimensions:	(W x H x D): approx. 50 x 250 x 65 mm
Weight:	approx. 270 g

Protection class of battery charger:	IPX7 (protected against the effects of temporary immersion)
Protection class of hand unit:	IPX4 (protected against water splashes)
LEDs:	RG 0
Ambient conditions:	10 °C to 40 °C Humidity 40 to 85% relative

Symbols used

	Certified Safety. Devices must comply with the generally acknowledged rules of technology and the German Product Safety Act (Produktsicherheitsgesetz - ProdSG).
	With the CE marking, HOYER Handel GmbH declares the conformity with EU guidelines.
	This symbol reminds you to dispose of the packaging in an environmentally friendly manner.
	Recyclable materials: cardboard (except corrugated cardboard)

Subject to technical modifications.

16. Warranty of the HOYER Handel GmbH

Dear Customer,
your device is provided with a 3 year warranty starting with the purchase date. In the event of product defects, you are entitled to statutory rights against the vendor. These statutory rights are not restricted by our warranty presented in the following.

Warranty conditions

The warranty period starts with the purchase date. Please keep the original purchase receipt in a safe place. This document is required to verify the purchase.

If within three years from the purchase date of this product a material or factory defect occurs, the product will be repaired or replaced by us – at our discretion – free of charge to you, or the purchase price will be refunded. This warranty implies that within the period of three years the defective device and the purchase receipt are presented, including a brief written description of the defect and the time it occurred.

If the defect is covered by our warranty, the repaired or a new product will be returned to you. No new warranty period starts with a repair or replacement of the product.

Warranty period and statutory claims for defects

The warranty period is not extended when the warranty has been claimed. This also applies to replaced and repaired parts. Any damages and defects already existing at the time of purchase must be reported immediately upon unpacking. Repairs arising after the expiration of the warranty period are subject to a charge.

Warranty coverage

The device was produced carefully according to strict quality guidelines and tested diligently prior to delivery.

The warranty applies to material or factory defects.

Excluded from the warranty are wear parts subject to normal wear and damages to fragile parts, e.g. switches, batteries, lamps or other parts manufactured from glass.

This warranty expires if the product is damaged, not used as intended or not serviced. For the proper operation of the product, all instructions listed in the operating instructions must be observed carefully. Any form of use and handling that is advised against in the operating instructions or warned against must always be avoided.

The product is only intended for private and not for commercial use. In the case of incorrect and improper treatment, use of force and interventions not performed by our authorised Service Centre, the warranty shall cease.

Handling in case of a warranty claim

In order to ensure prompt processing of your matter, please observe the following notes:

- Please keep the article number IAN: 320207 and the purchase receipt as a purchase verification for all inquiries.
- The article number can be found on the rating plate, an engraving, the title page of your instructions (in the bottom left), or as a label on the rear or underside of the device.
- If faulty operation or other defects occur, first contact the Service Centre listed in the following by telephone or email.
- Then, you are able to send a product reported as defective free of charge to the service address specified to you, including the purchase receipt and the information on the defect and when it occurred.

At www.lidl-service.com you can download this and many other manuals, product videos and installation software.

This QR code takes you directly to the Lidl Service website (www.lidl-service.com), where you can enter the article number (IAN) 320207 to open your user instructions.

Service Centre

(GB) Service Great Britain
Tel.: 0800 404 7657
E-Mail: hoyer@lidl.co.uk

(IE) Service Ireland
Tel.: 1890 930 034
(0,08 EUR / Min., (peak))
(0,06 EUR / Min., (off peak))
E-Mail: hoyer@lidl.ie

(MT) Service Malta
Tel.: 80062230
E-Mail: hoyer@lidl.com.mt

(CY) Service Cyprus
Tel.: 8009 4409
E-Mail: hoyer@lidl.com.cy

IAN: 320207

Supplier

Please note that the following address is no service address. First contact the aforementioned Service Centre.

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