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[www.healthforyou.lidl/privacy.html](http://www.healthforyou.lidl/privacy.html)

<b>List of pictograms used</b> .....	Page 18
<b>Introduction</b> .....	Page 18
Intended use .....	Page 19
Parts description .....	Page 19
Technical data .....	Page 19
<b>Safety notes</b> .....	Page 20
<b>Measurement principle</b> .....	Page 22
<b>Initial use</b> .....	Page 22
Inserting / replacing batteries .....	Page 22
Setting up the diagnostic scales .....	Page 23
Set-up with the "HealthForYou" app .....	Page 23
Set-up without the app (alternative) .....	Page 24
<b>Operation</b> .....	Page 24
Measuring .....	Page 25
Weight only .....	Page 25
Switching off the diagnostic scales .....	Page 25
Configuring users .....	Page 25
Transferring assigned measurements .....	Page 25
Retrieving saved measurements from the scales .....	Page 26
Deleting scales data .....	Page 26
<b>Cleaning and care</b> .....	Page 26
<b>Disposal</b> .....	Page 26
<b>Warranty / Service</b> .....	Page 27
<b>Guidelines / Conditions</b> .....	Page 27
Declaration of conformity .....	Page 27
Privacy policy .....	Page 27

List of pictograms used			
	Please read the instructions for use!		Risk of tipping: Always step onto the middle of the platform.
	Direct current/voltage		Place the scales on a level surface, not on carpet.
	Observe the warnings and safety notes!		Danger of slipping: Never step on the scales with wet feet.
	Safety notes		Dispose of the scales and packaging in an environmentally friendly manner!
	Instructions		Do not dispose of batteries in the normal domestic waste!
	Keep children away from the scales!		Product meets the applicable product-specific European directives
	Caution! Explosion hazard!		Not suitable for persons with medical implants (e.g. pacemaker)
	Wear safety gloves!		Packaging from responsible sources.
	Data transfer via <i>Bluetooth</i> <sup>®</sup> low energy technology		Do not load the scales beyond 180 kg / 396 lb / 28 st.
	3-year warranty		
	Manufacturer		The UKCA symbol confirms that the product conforms to the applicable norms in the United Kingdom.

## Body Analyser Scales SBF 77

### ● Introduction

This quick start guide forms an integral part of the instruction manual, which enables you to start up your product immediately. You may download the full instruction manual from this website <http://www.lidl-service.com>.

Before using the product, read the instruction manual and pay particular attention to the safety notes contained in it. Keep this quick start guide in a safe place. When passing this product on to third parties, be sure to include all documentation.

## ● Intended use

These scales are intended to determine the weight, body fat, body water and muscle percentage, and the bone mass in persons, and to compute the basal and active metabolic rate. The scales should not be used by pregnant people, dialysis patients, people with edemas, as well as people who have pacemakers or other medical implants. The scales are only intended for personal use, not for commercial purposes. Any use of the scales other than previously mentioned or any product modification is prohibited and may result in injuries and/or product damage. The manufacturer accepts no liability for any damage caused by improper use. The scales are not intended for commercial use.

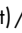
**Note:** The values and recommendations provided are not medically relevant. If your body fat or body water percentage is over or under the normal range on a long term basis, please consult a doctor.

## ● Parts description

- 1 LC-display
- 2 ▼ button
- 3 SET button
- 4 ▲ button
- 5 Mode button
- 6 Battery compartment
- 6a Battery compartment closure
- 7 Weighing platform
- 8 Electrode

### LC-Display

- 9 Weight
- 10 Unit of measurement (st, lb or kg)
- 11 BMR basal metabolic rate
- 12 AMR active metabolic rate
- 13 Metabolic rate in kcal
- 14 Age
- 15 Values in % or weight in lb or kg
- 16 Body data values
- 17 User initials, unknown value “- - -”
- 18 Gender
- 19 Date / time / birthday
- 20 Synchronisation symbol when transferring measurements

- 21 Bluetooth® symbol for connecting the diagnostic scales with a smartphone
- 22 Activity level
- 23 Body data mode BMI (Body Mass Indicator) / BF (body fat) /  body water percentage /  (muscle percentage) /  (bone mass)

## ● Technical data

Platform dimensions:	approx. 30 x 30 x 2.5 cm (W x D x H)
Battery:	3 x 1.5V battery (type AAA)
Maximum capacity:	180 kg / 396 lb / 28 st
Increments (weight, bone mass):	0.1 kg / 0.2 lb / 0.2 lb
Increments (body fat, body water, muscle percentage):	0.1 %
Increments (calories):	1 kcal
Material:	Safety glass, plastic, metal
Data transfer via Bluetooth® low energy technology	

This product uses Bluetooth® low energy technology, frequency band 2.402–2.480 GHz, transmission power max. 5 dBm. For a list of compatible smartphones, information about the “HealthForYou” app and software, and further details about the devices please visit: [www.healthforyou.lidl](http://www.healthforyou.lidl)

System requirements of the “HealthForYou” web version	MS Internet Explorer MS Edge Mozilla Firefox Google Chrome Apple Safari (respective current versions)
System requirements of the “HealthForYou” app	iOS version 12.0 and above, Android™ Version 8.0 and above, Bluetooth® 4.0 and above.

Technical specifications subject to change without notice due to updates.

# IMPORTANT INSTRUCTIONS – KEEP FOR FUTURE REFERENCE!

## Safety notes

**⚠ WARNING!** Failure to observe the following instructions may result in personal injury or property damage. The following safety notes and hazard warnings are not only intended to protect your own health and the health of others, but also to prevent damage to the scales. Please observe these safety notes and be sure to include them when passing the scales on to others.

■  **⚠ WARNING!**  
**DANGER TO LIFE AND RISK OF ACCIDENTS FOR INFANTS AND CHILDREN!**

Never leave children unsupervised with the packaging materials and the device. Otherwise there is a risk to life from suffocation. Always keep children away from the device. This product is not a toy.

■ This product may be used by children aged 8 years and up, as well as by persons with physical, sensory or mental impairments or lacking experience and/or knowledge, so long as they are

supervised or instructed in the safe use of the product and understand the associated risks. Do not allow children to play with the device. Cleaning and user maintenance must not be performed by children without supervision.

■ **⚠ WARNING!** This product is not suitable for persons with medical implants (e.g. pacemaker). The measuring voltage for the bioelectrical impedance analysis of the scales may interfere with the operation of such implants.

■ Do not use if you are pregnant.

■ **CAUTION! RISK OF INJURY!**

Do not load the scales beyond 180 kg/396 lb/28 st. Otherwise the scales might be damaged.

■ **CAUTION! RISK OF IN-**

**JURY!** Do not step on to the scales with wet feet and do not use the scales if the platform is moist. Otherwise you may slip.

■ **CAUTION! RISK OF INJURY!** Do not step

onto only one outside edge of the diagnostic scales: Tipping hazard!

- Protect your scales from impact, moisture, dust, chemicals and sudden changes in temperature, and keep them away from sources of heat (ovens, heating radiators).

- Do not expose the product to
  - strong mechanical stresses,
  - direct sunlight,
  - strong electromagnetic fields (e.g. mobile).

The product may otherwise be damaged.

- Check the product for damage before use. Do not use the product if damaged.
- Remove the batteries if the product malfunctions or is defective. Otherwise there is an increased risk of battery leaks.

**⚠ WARNING!** Never open the product housing. Doing so may result in personal injury or product damage. Always have any repairs carried out at the service centre or by an electrical equipment repair specialist.

- **CAUTION!** Do not store the scales in passageways. Otherwise you or others may trip over them.
- Do not allow the connection terminals to short circuit.

## Dangers associated with handling batteries

### CAUTION! DANGER OF EXPLOSION!

- Never recharge, short-circuit, mechanically deform, crush or open single-use batteries. This could cause batteries to overheat, catch fire or explode. Never use incorrect or differing battery types. Never throw batteries into fire or water and never place them on hot surfaces. The batteries may explode. Avoid extreme conditions and temperatures over 50 °C which may affect batteries, e.g. on radiators, fires, ovens. Otherwise there is a danger of explosions or a danger of flammable gas or liquid leaking.
- Be sure to insert the batteries correctly (polarity)!
- Batteries must not be disposed of in the normal domestic waste!
- Every consumer is legally obliged to dispose of batteries properly!
- Keep batteries out of the reach of children.
- Avoid contact with the skin, eyes or mucous membranes if fluid has leaked from the batteries. In the event of contact with battery acid, thoroughly flush the affected

area with clean water and seek immediate medical attention.

- Remove the batteries from the product if it will not be used for a long period of time.
- Remove drained batteries from the product immediately. There is an increased risk of leaks!
- Failure to observe the above instructions may result in the batteries discharging beyond their end voltage. This poses a risk of leaks. If the batteries have leaked inside your product, immediately remove the batteries to prevent product damage!
-  Leaked or damaged batteries may cause burns on contact with the skin. To avoid this, always wear suitable protective gloves in these circumstances!
- New and used batteries must not be used together.

## ● Measurement principle

These diagnostic scales use the principle of bioelectrical impedance analysis (BIA). This enables body parts to be measured within a few seconds using a non-perceivable, completely harmless and safe current.

By measuring the electrical resistance (impedance) and factoring in constants or individual values (age, height, gender, activity level), the body fat percentage and other values can be determined approximately. Muscle tissue and water have good electrical conductivity, hence little resistance. Bone and fatty

tissue on the other hand have a low conductivity, as the fat cells and bones hardly conduct the current due to their high resistance.

Please note, the values determined by the diagnostic scales are only approximate values of the body's medical, actual analytical values. Only a medical specialist can determine the exact body fat or water content percentages, the muscle mass or the bone structure through medical procedures (e.g. computer tomography).

Please note, due to the technology used measuring tolerances may occur, as these are not calibrated diagnostic scales for professional, medical purposes.

This product is not suitable for the following persons, as the values may be incorrect or implausible when determining the body fat percentage and other values:

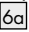
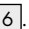
- Children under approx. 10 years  
**Note:** For children from 8 to 10 years, the device is only suitable for measuring weight.
- Competitive athletes and bodybuilders,
- Persons with a fever, receiving dialysis, oedema symptoms or osteoporosis,
- Persons taking cardiovascular medications (relating to the heart and blood vessels),
- Persons taking vasodilatory or vasoconstrictive medications,
- Persons with significant anatomical variances of the legs related to total body height (considerably shorter or longer leg).

## ● Initial use

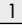
- Verify the parts listed in the scope of delivery are complete.
- Remove all packaging material.

## ● Inserting / replacing batteries

The batteries must be inserted prior to initial use.

- Open the battery compartment cover by pushing the battery compartment closure  in and up with your finger (see Fig. A).
- Insert the three included batteries (with the correct polarity) in the battery compartment .

Please note the marking inside the battery compartment .

**Note:** The LC-display  will show "USE APP" for approx. 10 minutes. If the product does not respond, remove and reinsert the batteries.

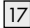
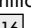
**Note:** When replacing the battery later, all the saved data will be retained.

## ● Setting up the diagnostic scales

-  Place the product on a level, solid surface. The product is not suitable for use on carpeted floors. A solid floor covering is required to ensure accurate values.

## ● Set-up with the "HealthForYou" app

Download the free app "HealthForYou" from the Apple App Store and Google Play.

- Activate *Bluetooth*® in the smartphone settings.
- To be able to connect to *Bluetooth*®, be sure to keep your smartphone/tablet within the *Bluetooth*® range of the product.
- Install "HealthForYou" from the Apple App Store / Google Play.
- Start the app and follow the instructions. Select the scales in the settings menu of the app. A randomly generated six-digit PIN code will be shown on the display on the scales (see fig. C). This is made up of the initials display  and measurements display . At the same time, an entry field (pairing request) will appear on your smartphone, where you should enter this six-digit PIN code.

- You will need to enter the following information in the "HealthForYou" app:

User data	Configuration Options
Initials	up to 3 letters or numbers (for up to 8 users)
Body height	100 to 220 cm (3' 3.5" to 7' 2.5")
Birthday	Year, month, day (age: 10-100 years)
Gender	Male (♂), female (♀)
Activity level (✖)	1 to 5 The medium- and long-term information is important when determining the activity level.

Activity level	Physical activity
1	None.
2	Low: Little and light exercise (e.g. walks, light gardening, gymnastic exercises).
3	Moderate: Exercising at least 2 to 4 times per week for 30 minutes at a time.
4	High: Exercising at least 4 to 6 times per week for 30 minutes at a time.
5	Very high: Intensive exercise, training or hard physical labour for at least 1 hour every day.

Set up the users as prompted.

- The first measurement must be matched with your personal user information for automatic user recognition. Follow the instructions in the app. To do so, step onto the diagnostic scales barefoot and be sure to remain still, standing on the electrodes with both legs, weight distributed evenly.

## ● Set-up without the app (alternative)

### Note:

We recommend setting up and implementing all settings using the "HealthForYou" app.

Only carry out this step if you want to set up entirely without the app.

If you have set up with the app, go straight to the chapter "Operation".

- Press the SET button [3] while "USE APP" is flashing to set the year/month/day and the hour/minute. If "USE APP" is not flashing, briefly load the scales.
- To do this, press the buttons ▼ [2] and ▲ [4]. Confirm your settings by pressing the SET button [3]. The LC-display [1] will show "0.0".
- Press the SET button [3] again. The LC-display [1] will show "P-01" or the number of the user memory to last use the scales.
- Use the buttons ▼ [2] and ▲ [4] to select your user memory (P-01 to P-08) and press the SET [3] button to confirm.
- Enter your settings, such as gender, initials, height, date of birth (year/month/day) and your activity level (see "Set-up using the "HealthForYou" app" (table)).
- Confirm these by pressing the SET button [3] each time.

Once you have finished making your entries, the 4 digit user PIN will appear. You might need this to log into the app later. The LC-display [1] will then show "0.0".

**Note:** For automatic user recognition, the first measurement must be assigned to your personal user information. Ideally, stand barefoot on the diagnostic scales and be sure to stand still on the electrodes with your weight distributed equally between both feet.

**Note:** You can deactivate Bluetooth® on the scales. To do so, when the product is switched on, press the mode button [5] for approx. 6 seconds until the Bluetooth® symbol [21] and "ON" appear in the LC-display [1]. Within 3 seconds, use the ▲ button [4] to switch to the deactivated status.

"OFF" will appear in the LC-display [1]. After approx. 3 seconds, the scales will switch off automatically with this change.

### Note for logging into the app later:

Should you have set up the scales without the app, you can add the existing user memory of your scales to the "HealthForYou" app. For this you will need to enter your 4 digit user PIN. This process is described in the app.

## ● Operation

- To obtain comparable results, where possible use the scales at the same time of day (preferably in the morning) after using the toilet, on an empty stomach and wearing clothing of a similar weight.
- Allow several hours following unusual physical strain.
- Wait approx. 15 minutes after getting up to allow water to distribute throughout the body.
- Important for measurements: the body fat percentage can only be determined if you are barefoot.

**Note:** Completely dry or very calloused feet may provide poor results, as the conductivity will be too low.

**Note:** Do not allow the feet, legs, calves or thighs to touch. This will interfere with measurement.

**Note:** Wearing socks will provide false results.

- Stand upright and still during measurement, with the weight distributed equally between both legs.
- Please note, only the long-term trend is important. Sudden weight changes within a few days are typically caused by loss of fluids; however, hydration is important for your well being.

**Note:** The units of measure for the scales can be changed in the app, under device settings (kg/lb/st). Alternatively, by briefly pressing the Mode button [5] on the back of the scales as many times as needed, you can select the desired unit.

When selecting the weight unit lb, the time mode automatically changes to 12-hour format (a.m./p.m.).

## ● Measuring

- Step on the diagnostic scales barefoot and be sure to stand still on the electrodes [8] with the weight distributed equally between both legs.

The diagnostic scales will start to measure immediately. The results will appear shortly after.

If a user has been assigned (the user initials [17] will appear, e.g. LEE), BMI, BF, water, muscle, bone, BMR and AMR will be displayed.

The following information will be displayed:

- Weight in kg with BMI (Fig. D)
- Body fat in % (BF) (Fig. E)
- Body water in % (Fig. F)
- Muscle mass in % (Fig. G)
- Bone mass in kg (Fig. H)
- Basal metabolic rate in kcal (BMR) (Fig. I)
- Active metabolic rate in kcal (AMR) (Fig. J)

**Note:** If the measurement value can be associated with several possible user profiles, the user display switches to LC display [1]. Select your user by pressing the SET button [3].

If a user is not recognised, only the weight will be displayed along with “-.-”. You can then use the ▲ button [4] to select your user. This is activated after approx. 3 seconds. If you do not do so, this measurement will not be saved.

**Note:** Alternatively, you can also select your user before taking the measurement.

Briefly step on the diagnostic scales to activate them. Select your user using the ▲ button [4]. The selected user is activated after approx. 3 seconds.

## ● Weight only

- Stand on the diagnostic scales wearing shoes. Stand still on the diagnostic scales with the weight distributed evenly between both legs.

The diagnostic scales will start to measure immediately. If a user has been assigned (the user initials [17] will be displayed), the weight and BMI will be displayed and saved.

**Note:** For two or more possibly recognised users, the user changes in the display. Select your user by pressing the SET button [3].

If a user is not recognised, only the weight will be displayed along with “-.-”. You can then select

your user using the buttons ▼ [2] and ▲ [4] and confirm using the SET button [3]. If you do not do so, this measurement will not be saved.

## ● Switching off the diagnostic scales

The diagnostic scales feature an auto shut-off and will automatically switch off after displaying the results.

## ● Configuring users

The weighing scales are able to save measurements from up to 8 different users. Upon starting a new measurement, the diagnostic scales will automatically assign the measurement to the created user if the last saved measurement was within +/- 3kg of the current value.

If the measurement value can be associated with several possible user profiles, the user display switches to LC display [1]. Select your user by pressing the SET button [3].

**Note:** Alternatively, you can also select your user before taking the measurement.

Briefly step on the diagnostic scales to activate them. Select your user using the ▲ button [4]. The selected user is activated after approx. 3 seconds.

## ● Transferring assigned measurements

If the “HealthForYou” app is open and actively connected via Bluetooth® to the diagnostic scales (indicated by the Bluetooth® symbol [21] in the LC-display [1]), newly assigned measurements will be immediately transferred to the app. The synchronisation symbol [20] appears in the LC-display [1] when new measurements are being transferred to the app. It does not appear if no new measurements were saved or recorded for this user. In this case the values will not be saved to the diagnostic scales. If the app is not open, newly assigned values will be saved to the diagnostic scales. The diagnostic scales can save up to 30 values per user. The saved

values are automatically transferred to the app after opening the app within *Bluetooth*<sup>®</sup> range. Data can automatically be transferred with the diagnostics scales off, provided you are within *Bluetooth*<sup>®</sup> range.

## ● Retrieving saved measurements from the scales

- When the product is switched on, select the user by pressing the ▲ button [4] as many times as necessary. The selected user is activated after approx. 3 seconds.
- Press the ▼ button [2] as many times as necessary to retrieve past values.
- Measurements that have already been transferred can no longer be viewed on the LC-display [1].

## ● Deleting scales data

- To completely delete all values and outdated user data from the diagnostic scales, switch on the diagnostic scales and hold the Mode button [5] for approx. 11 seconds. The LC-display [1] will read "CLr" for a few seconds.
- To delete just one user, switch the diagnostic scales on and press the ▲ button [4] to select the user. This will be activated after approx. 3 seconds and the user's initials along with "0.0" will appear in the LC-display [1]. Then press the Mode button [5] for approx. 6 seconds. The LC-display [1] will read "dEL" for a few seconds.

## ● Cleaning and care

- Do not clean with chemical cleaners or scouring agents.
- Clean the product with a damp cloth. For more stubborn dirt you may also apply a weak soap solution to the cloth.
- Never immerse the product in water or other liquids. The devices may otherwise be damaged.

- Ensure water does not enter the product. If this does occur, wait for the devices to dry completely before using it again.

## ● Disposal

The packaging is made entirely of recyclable materials, which you may dispose of at local recycling facilities.

Observe the marking of the packaging materials for waste separation, which are marked with abbreviations (a) and numbers (b) with following meaning: 1-7: plastics / 20-22: paper and fibre-board / 80-98: composite materials.

The product and packaging materials are recyclable, dispose of it separately for better waste treatment. The Triman logo is valid in France only.

Contact your local refuse disposal authority for more details of how to dispose of your worn-out product.

To help protect the environment, please dispose of the product properly when it has reached the end of its useful life and not in the household waste. Information on collection points and their opening hours can be obtained from your local authority.

Per Directive 2006/66/EC and its amendments, defective or used batteries must be recycled. Return batteries and/or the product to the provided recycling facilities.

**Improper disposal of batteries can harm the environment!**

Never dispose of batteries in your household bin. They may contain toxic heavy metals and are subject to hazardous waste regulations. The chemical symbols of the heavy metals are: Cd = cadmium,

Hg = mercury, Pb = lead. Therefore dispose of used batteries through your local collection facilities.

## ● Warranty / Service

The scales are covered by a 3 year warranty from the date of purchase for defects in material and manufacture. The warranty does not cover:

- Damage due to improper use.
- Wear parts.
- Damage caused by the customer.
- If the device has been opened by an unauthorised repair shop.

This warranty does not affect the customer's legal rights. The customer must provide proof of purchase in order for any warranty claims within the warranty period to be honoured. The warranty claim must be filed with Hans Dinslage GmbH, Uttenweiler, Germany, within 3 years from the date of purchase.

**For complaints, please contact our Service Department at:**

**Service Hotline (free of charge):**

GB

☎ **0800 93 10 319 (Mon–Fri: 8 am–6 pm)**

✉ **service-uk@sanitas-online.de**

IE

☎ **1800 93 92 35 (Mon–Fri: 8 am–6 pm)**

✉ **service-ie@sanitas-online.de**

If we ask you to return the defective diagnostic scales, return the product to:

NU Service GmbH  
Lessingstraße 10b  
89231 Neu-Ulm  
GERMANY

**IAN 366720\_2101**

## ● Guidelines / Conditions

### ● Declaration of conformity CE

We hereby confirm that this product complies with the European RED Directive 2014/53/EU.

The CE declaration of conformity for this product is available at:

[www.sanitas-online.de/web/de/landingpages/cedeclarationofconformity.php](http://www.sanitas-online.de/web/de/landingpages/cedeclarationofconformity.php)

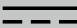
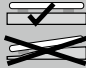
### ● Privacy policy

The "HealthForYou" privacy policy can be viewed via the following link:

[www.healthforyou.lidl/privacy.html](http://www.healthforyou.lidl/privacy.html)

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The *Bluetooth*<sup>®</sup> word mark and logos are registered trademarks owned by *Bluetooth*<sup>®</sup> SIG, Inc. and any use of such marks by Hans Dinslage GmbH is under license. Other trademarks and trade names are those of their respective owners.

Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries.

Google Play and the Google Play logo are trademarks of Google LLC.

Android is a trademark of Google LLC.

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Riedlinger Str. 28  
88524 Uttenweiler  
GERMANY

Stand der Informationen · Last Information Update · Version des  
informations · Stand van de informatie · Stan informacii · Stav informací  
Stav informácií · Estado de las informaciones · Tilstand af information:  
06/2021 · Ident.-No.: 749.969\_SBF77\_2021-03-08\_01\_IM-QG\_Lid\_LB8\_OS

IAN 366720\_2101