MUNDDUSCHE NMD 3.7 B4
WATER JET FLOSSER NMD 3.7 B4
HYDROPULSEUR NMD 3.7 B4

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MUNDDUSCHE

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IRYGATOR DO ZĘBÓW

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ÚSTNA SPRCHA

GB IE Operating instructions

WATER JET FLOSSER

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MONDDOUCHÉ

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ÚSTNÍ SPRCHA

HOYER Handel GmbH

Tasköprüstraße 3
DE-22761 Hamburg
GERMANY

Stand der Informationen · Status of information · Version des informations · Stand van de informatie
Stan informáci · Stav informací · Stav informácií:

07/2019 ID: NMD 3.7 B4_19_V1.2

IAN 326029_1904

IAN 326029_1904







Deutsch	2
English	16
Français	30
Nederlands	48
Polski	62
Česky	76
Slovenčina	90

Übersicht / Overview / Aperçu de l'appareil / Overzicht / Przegląd / Přehled / Prehľad

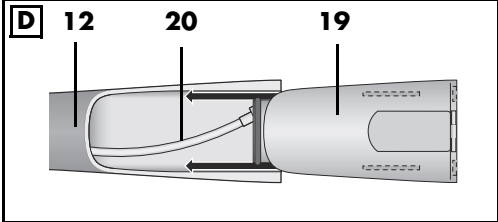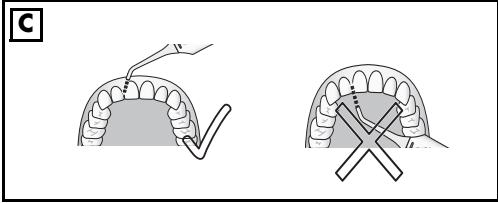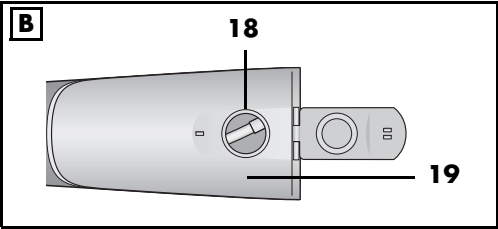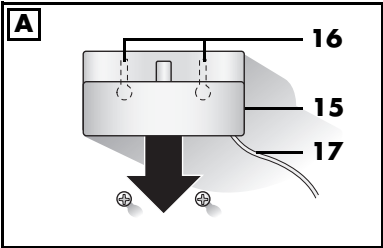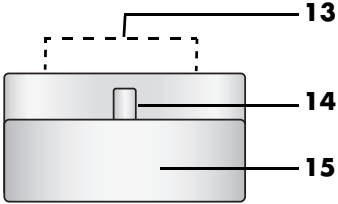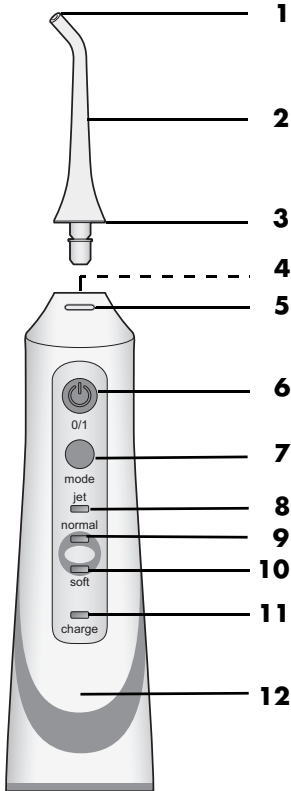Inhalt

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DEUTSCHLAND

Contents

1. Overview	17
2. Intended purpose	17
3. Safety information	18
4. Items supplied	20
5. Battery charger	20
6. Charging the battery	20
7. Nozzle attachments	21
8. Operation	21
8.1 Filling the water tank	21
8.2 Setting the water jet intensity	22
8.3 Before initial use	22
8.4 Using the dental water jet	22
8.5 After every use	23
9. Cleaning / Hygiene requirements	23
10. Storage	24
11. Disposal	24
12. Troubleshooting	25
13. Ordering accessories	25
14. Technical specifications	26
15. Warranty of the HOYER Handel GmbH	27

1. Overview

1	Nozzle
2	Nozzle attachment (replaceable)
3	Marking ring
4	Holder (for nozzle attachment)
5	Release button (for nozzle attachment)
6	⏻ On/off switch
7	mode Select water jet intensity
8	jet LED: setting: strong water jet
9	normal LED: setting: normal water jet
10	soft LED: setting: soft water jet
11	charge Charging indicator light
12	Dental water jet
13	Holders (for storing the nozzle attachments)
14	Peg of the battery charger
15	Battery charger
16	Holes for wall mounting (on the back)
17	Mains power cable and mains plug
18	Filling opening
19	Water tank
20	Suction hose (in the water tank)

(2 raw|plugs and 2 screws without image)

Thank you for your trust!

Congratulations on the purchase of your new dental water jet.

For a safe handling of the device and in order to get to know the entire scope of features:

- **Thoroughly read these user instructions prior to initial use.**
- **Above all, observe the safety information!**
- **The device should only be used as described in these user instructions.**
- **Keep these user instructions for reference.**

- **If you pass the device on to someone else, please include this copy of the user instructions. The user instructions are a part of the product.**

We hope you enjoy using your new dental water jet!

2. Intended purpose

The dental water jet is intended to be used exclusively for cleaning teeth and oral hygiene.

The device is designed for private, domestic use. This device must not be used for commercial purposes.

3. Safety information

Warnings

If necessary, the following warnings will be used in these user instructions:

DANGER! High risk: failure to observe this warning may result in injury to life and limb.

WARNING! Moderate risk: failure to observe this warning may result in injury or serious material damage.

CAUTION: low risk: failure to observe this warning may result in minor injury or material damage.

NOTE: circumstances and specifics that must be observed when handling the device.

Instructions for safe operation

- ⊙ This device can be used by children from the age of 8 and people with restricted physical, sensory or intellectual abilities or people without adequate experience and/or understanding if they are supervised.
- ⊙ Cleaning and user maintenance must not be performed by children unless they are supervised.
- ⊙  The battery charger is connected to the mains voltage, and must therefore not be used in the shower, in the bath, or over wash basins filled with water.
- ⊙ Charge the hand unit of the dental water jet only with the original battery charger.
- ⊙ The dental water jet may only be used with cold or warm water. Use only clear tap water without additives.
- ⊙ Children must not play with the device.
- ⊙ The mains power cable of this device cannot be replaced. If the cable should become damaged, the device itself must be disposed of.

DANGER as small parts may be swallowed

- ⊙ The marking rings on the nozzle attachments are small parts which can be swallowed. Keep away from children under 36 months of age.

DANGER for children

- ⊙ Children must not play with packing material. Do not allow children to play with plastic bags. There is a risk of suffocation.

DANGER to and from pets and livestock

- ⊙ Electrical devices can represent a hazard to pets and livestock. In addition, animals can also cause damage to the device. For this reason you should keep animals away from electrical devices at all times.

DANGER! Risk of electric shock due to moisture

- ⊙ The dental water jet, the battery charger, the mains power cable and the mains plug must not be immersed in water or any other liquids.
- ⊙ If the battery charger has fallen into water, pull out the mains plug immediately, and only then remove the battery charger from the water. In this case do not use the battery charger again, but have it checked by a specialist workshop.

DANGER! Risk of electric shock

- ⊙ The installation of a residual current device with a nominal trip current not exceeding 30 mA in the building wiring offers additional protection. Contact a skilled electrician in this regard.
- ⊙ Do not start using the device if the device, the battery charger or the mains power cable shows visible damage, or if the device has previously been dropped.
- ⊙ The mains power cable of this device cannot be replaced. If the cable should become damaged, the device itself must be disposed of.
- ⊙ Only connect the mains plug to a properly installed and easily accessible wall socket whose voltage corresponds to the specifications on the rating plate. The wall socket must continue to be easily accessible after the device is plugged in.
- ⊙ Ensure that the mains power cable cannot be caught or crushed, and that it

cannot be damaged by sharp edges or hot points.

- ⊙ To fully disconnect the battery charger from the mains power supply, pull out the mains plug.
- ⊙ When removing the mains plug from the wall socket, always pull the plug and never the cable.
- ⊙ Disconnect the mains plug from the wall socket ...
 - ... if there is a fault,
 - ... when you are not using the device,
 - ... before you clean the device and
 - ... during thunderstorms.
- ⊙ To avoid any risk, do not make modifications to the device. Repairs must be carried out by a specialist workshop or the Service Centre.
- ⊙ When using an extension cable, this must correspond to the applicable safety regulations.

DANGER from battery

- ⊙ Do not expose the device to direct sunshine or heat. The ambient temperature should not be lower than 0 °C or higher than +40 °C.
- ⊙ Protect the battery against mechanical damages. **Fire hazard!**
- ⊙ Should the electrolyte solution leak out of the battery, avoid contact with eyes, mucous membranes and skin. Rinse immediately with plenty of clean water and consult a doctor. The electrolyte solution can cause irritation.
- ⊙ The device is provided with a lithium-ion battery.
 - The battery cannot be removed!
 - This device must not be opened!
 - The device with the installed battery must be disposed of properly.

WARNING! Risk of material damage

- ⊙ If you wish to attach the battery charger to the wall, check what fixing material

is suitable. If in doubt, have the assembly performed by an expert.

- ⊙ Before drilling, check that there are no water pipes or power cables in the wall at the relevant point.
- ⊙ The battery charger is fitted with non-slip silicon feet. As furniture is coated with a variety of different varnishes and acrylics, and treated with different cleaning agents, it cannot be completely ruled out that components are contained in some of these substances that could damage and soften the plastic feet. If necessary, place a non-slip mat under the battery charger.
- ⊙ Never place the device on hot surfaces (e.g. hot plates) or near heat sources or open fire.
- ⊙ Do not use any astringent or abrasive cleaning agents.

CAUTION

- ⊙ In the first few days of use of a dental water jet, some users may experience slight gum bleeding. If the bleeding lasts longer than 14 days, a dentist should be consulted.
- ⊙ If you suffer from severe periodontitis, have injuries or ulcers in the oral mucosa or have had dental surgery during the last two months, please consult a dentist before using the dental water jet.

4. Items supplied

- 1 dental water jet **12**
- 2 nozzle attachments **2**
- 1 battery charger **15**
- 2 rawplugs, diameter 5 mm
screws 3 x 25 mm
- 1 copy of the user instructions

- Remove all packing material.
- Check whether all accessories are present and undamaged.

5. Battery charger

The battery charger **15** can either be placed on a surface or attached to a wall.

- Insert the mains power cable **17** through one of the 3 recesses in the housing of the battery charger **15**.

Wall mounting

NOTES:

- Please note the rawplug size when selecting the drill.
 - If your wall material requires it, replace the fixing material supplied with a suitable one.
-
- **Figure A:** mount the battery charger **15** to a wall according to the holes **16** on the back as shown in the illustration.

6. Charging the battery

NOTES:

- Fully charge the battery before using it for the first time. This takes 14 - 16 hours.
 - Place the dental water jet **12** on the peg **14** of the battery charger **15** after use. This ensures that the dental water jet **12** always works at full power. The battery is not damaged by regular charging.
-
1. Plug the mains plug **17** of the battery charger **15** into a suitable wall socket. The wall socket must continue to be easily accessible after the device is plugged in.
 2. Place the switched-off dental water jet **12** on the peg **14** of the battery charger **15**. The charging indicator light **charge 11** is lit in red throughout the entire charging process. When the optimal charging

state has been reached, the charging indicator light **charge 11** is lit green.

Battery signals: charging indicator light charge 11

- Dental water jet **12** on battery charger **15**: charging indicator light **charge 11** is lit green: the battery is fully charged.
- Dental water jet **12** on battery charger **15**: charging indicator light **charge 11** is lit red: the battery is being charged.
- Dental water jet **12** switches itself off: charging indicator light **charge 11** flashes red. Charge the battery.
- The battery becomes weak during use: charging indicator light **charge 11** flashes red. Charge the battery.

7. Nozzle attachments

The items supplied include 2 nozzle attachments **2**.

The various marking rings **3** allow you to assign the nozzle attachments **2** to different family members.

The nozzle attachments **2** can be inserted into the holders **13** of the battery charger **15** for storage.

Attaching/removing the nozzle attachment

- Take a nozzle attachment **2** and insert it from above into the holder **4** of the dental water jet **12** so that it clicks into place.
- To remove it, hold down the release button **5** and carefully pull the nozzle attachment **2** upwards.

NOTE: it is advisable to replace the nozzle attachments **2** every 12 months (see "Ordering accessories" on page 25).

8. Operation

CAUTION:

- ⊙ In the first few days of use of a dental water jet, some users may experience slight gum bleeding. If the bleeding lasts longer than 14 days, a dentist should be consulted. Please also note the remarks on the water jet intensity (see "Setting the water jet intensity" on page 22).
- ⊙ If you suffer from severe periodontitis, have injuries or ulcers in the oral mucosa or have had dental surgery during the last two months, please consult a dentist before using the dental water jet.

8.1 Filling the water tank

You can fill the water tank with water without removing it from the device.

WARNING! Risk of material damage!

- ⊙ Use only clean tap water without additives to fill the water tank. Additives (such as mouthwashes) or other liquids can damage the device.
 - ⊙ The water should have a pleasant temperature. Do not use water that is hotter than 40 °C.
1. Hold the switched-off dental water jet **12** horizontally.
 2. **Figure B:** open the cover over the filling opening **18**.
 3. Fill the water tank **19** through the filling opening **18**.
 4. Close the cover over the filling opening **18**.

8.2 Setting the water jet intensity

You can set three different water jet intensities on the dental water jet:

LED	Jet intensity
soft	soft water jet: at the beginning, for sensitive gums, for gentle gum massage
normal	normal water jet
jet	strong water jet: for advanced users, strong gums, intensive cleaning of the spaces between the teeth

NOTE: if you have not used a dental water jet before, please start on the **soft** setting so that your gums can get used to the water jet. At a later time you can increase the water jet intensity to **normal** and **jet**. The pressure should always remain pleasant.

- Press the **mode 7** button repeatedly until the desired LED **10**, **9** or **8** is lit.

8.3 Before initial use

1. Rinse the nozzle attachments **2** under running water.
2. Clean the water tank **19** and the dental water jet **12** (see "Cleaning / Hygiene requirements" on page 23).
3. Fill the water tank **19** with clean tap water (see "Filling the water tank" on page 21).
4. Insert a nozzle attachment **2** into the holder **4** so that it clicks into place.
5. Point the nozzle **1** at a washbasin and switch the dental water jet **12** on with the on/off switch **6**.

6. Let the dental water jet **12** run until the water tank **19** is empty in order to rinse through the pipes inside. The dental water jet is now ready for use.

8.4 Using the dental water jet

1. Take the dental water jet **12** off the battery charger **15**.
2. Insert a nozzle attachment **2** into the holder **4** if necessary so that it clicks into place.
3. Fill the water tank **19** with fresh tap water.
4. To set the water jet intensity, press the **mode 7** button repeatedly until the desired LED is lit:
 - LED **soft 10**:
soft water jet
 - LED **normal 9**:
normal water jet
 - LED **jet 8**:
strong water jet
5. Insert the nozzle **1** into your slightly open mouth and lean over a washbasin so that water and rinsed-out food particles can flow out.
6. Switch on the dental water jet **12** with the on/off switch **6**.
7. To start the cleaning process, direct the water jet at the spaces between your teeth and the edges of the teeth.
Figure C: hold the nozzle **1** in such a way that the water jet can flush the food residues away from your gums.
8. Move the nozzle **1** slowly from one space between two teeth to the next. Use this method to clean all of the external and internal surfaces.
9. To massage your gums, direct the water jet at your gums.

10. When no more water is coming out of the nozzle **1**, the water tank **19** is empty. Switch off the dental water jet **12** and replenish the water tank. When the device is switched on again, the device will automatically start on the last selected water jet setting.
11. To switch off the dental water jet **12**, press the on/off switch **6**. Do not take the nozzle **1** out of your mouth until the dental water jet has been switched off.

8.5 After every use

1. Open the cover over the filling opening **18**, pour out the remaining water from the water tank **19** and close the cover again.
2. Switch on the dental water jet **12** for a few seconds with the on/off switch **6** to remove any remaining water from the pipes inside the dental water jet.
3. Hold down the release button **5** and gently pull the nozzle attachment **2** upwards.
4. Rinse the nozzle attachment **2** under running water and insert it into a holder **13**.
5. Place the dental water jet **12** on the battery charger **15**.

NOTE: it is advisable to replace the nozzle attachments **2** every 12 months (see "Ordering accessories" on page 25).

9. Cleaning / Hygiene requirements

DANGER! Risk of electric shock!

- ⊙ Before cleaning the battery charger **15**, pull the mains plug **17** out from the mains socket.
- ⊙ Do not immerse the dental water jet **12** and/or the battery charger **15** in water!
The battery charger **15** must be completely dry before you reconnect it to the wall socket.

Nozzle attachments

1. Hold down the release button **5** and gently pull the nozzle attachment **2** upwards.
2. Rinse the nozzle attachment **2** under running water and insert it into a holder **13**.

Water tank

Clean the water tank **19** once a week or before you are not going to use the device for more than one week.

1. **Figure D:** push the water tank **19** downwards off the dental water jet **12**.
2. Rinse the water tank **19** thoroughly both inside and out with warm water.
3. Dry the water tank **19** on the outside and pour out the remaining water.

NOTES:

- if you are not going to use the dental water jet **12** for more than one week, allow the water tank **19** to dry completely before you put it back into the dental water jet.
- It is also recommended that the device be disinfected at regular intervals. Special disinfectant powders and prosthesis cleaners are suitable. Dilute the solution with water and place it into the

water tank, and then leave it to work for approx. 2 hours. Then let the dental water jet **12** run until the water tank **19** is empty in order to rinse through the pipes inside. Then rinse with plenty of water.

4. Push the water tank **19** from below into the dental water jet **12**. Make sure that
 - the suction hose **20** is located in the water tank **19**,
 - the water tank **19** is positioned in the recesses on both sides and
 - the water tank **19** is firmly in place so that the seal provides a tight fit.

Dental water jet and battery charger

- Rinse off the dental water jet **12** and rub it dry.
- Wipe the battery charger **15** with a damp cloth if necessary, and rub it dry.

10. Storage

- Remove any excess water from the device and let all of the parts dry completely before you store them away.
- Disconnect the device from the mains and store it in a dry location.

11. Disposal

DANGER from the battery!

- ⊙ The device contains a lithium-ion battery.
 - The battery cannot be removed!
 - This device must not be opened!
 - The device with the installed battery must be disposed of properly.

The lithium-ion battery used for this device may not be treated as household waste.

The battery must be disposed of properly, together with the device. Contact local waste disposal points for this.

This product is subject to the provisions of European Directive 2012/19/EC. The symbol showing a wheelie bin crossed through indicates

that the product requires separate refuse collection in the European Union. This applies to the product and all accessories marked with this symbol. Products identified with this symbol may not be discarded with normal household waste, but must be taken to a collection point for recycling electric and electronic appliances.

This recycling symbol marks e.g. an object or parts of the material as valuable for recycling. Recycling helps to reduce the consumption of raw materials and protect the environment.

Packaging

When disposing of the packaging, make sure you comply with the environmental regulations applicable in your country.

12. Troubleshooting

If your device fails to function as required, please try this checklist first. Perhaps there is only a minor problem, and you can solve it yourself.

DANGER! Risk of electric shock!

- ⊙ Do not attempt to repair the device yourself under any circumstances.

Fault	Possible causes / Action
The device cannot be switched on.	<ul style="list-style-type: none">• Is the battery empty?
Charging indicator light charge 11 is not lit.	<ul style="list-style-type: none">• Is the battery charger 15 connected to the mains?

13. Ordering accessories

Accessories for the dental water jet NMD 3.7 B4 can be reordered.

Order online
shop.hoyerhandel.com

1. Scan the QR code with your smartphone/ tablet.
2. With the QR code you will access a website where you can submit your re-order.

14. Technical specifications

Model:	NMD 3.7 B4
Mains voltage:	100 - 240 V ~ 50/60 Hz
Protection class:	II 
Power rating:	1.4 W
Battery:	1x3.7 V === (DC voltage) Li-Ion 800 mAh 2.96 Wh
Charging time (first charging):	14-16 hours
Protection class of battery charger:	IPX7 (protected against the effects of temporary immersion)
Protection class of hand unit:	IPX4 (protected against water splashes)
LEDs:	RG 0
Ambient conditions:	10 °C to 40 °C Relative humidity 40 to 85 %

Symbols used

	Double insulation
	Certified Safety. Devices must comply with the generally acknowledged rules of technology and the German Product Safety Act (Produktsicherheitsgesetz - ProdSG).
	With the CE marking, HOYER Handel GmbH declares the conformity with EU guidelines.
	This symbol reminds you to dispose of the packaging in an environmentally friendly manner.
	Recyclable materials: cardboard (except corrugated cardboard)

Subject to technical modifications.

15. Warranty of the HOYER Handel GmbH

Dear Customer,
your device is provided with a 3 year warranty starting with the purchase date. In the event of product defects, you are entitled to statutory rights against the vendor. These statutory rights are not restricted by our warranty presented in the following.

Warranty conditions

The warranty period starts with the purchase date. Please keep the original purchase receipt in a safe place. This document is required to verify the purchase.

If within three years from the purchase date of this product a material or factory defect occurs, the product will be repaired or replaced by us – at our discretion – free of charge to you, or the purchase price will be refunded. This warranty implies that within the period of three years the defective device and the purchase receipt are presented, including a brief written description of the defect and the time it occurred.

If the defect is covered by our warranty, the repaired or a new product will be returned to you. No new warranty period starts with a repair or replacement of the product.

Warranty period and statutory claims for defects

The warranty period is not extended when the warranty has been claimed. This also applies to replaced and repaired parts. Any damages and defects already existing at the time of purchase must be reported immediately upon unpacking. Repairs arising after the expiration of the warranty period are subject to a charge.

Warranty coverage

The device was produced carefully according to strict quality guidelines and tested diligently prior to delivery.

The warranty applies to material or factory defects.

Excluded from the warranty are wear parts subject to normal wear and damages to fragile parts, e.g. switches, batteries, lamps or other parts manufactured from glass.

This warranty expires if the product is damaged, not used as intended or not serviced. For the proper operation of the product, all instructions listed in the operating instructions must be observed carefully. Any form of use and handling that is advised against in the operating instructions or warned against must always be avoided.

The product is only intended for private and not for commercial use. In the case of incorrect and improper treatment, use of force and interventions not performed by our authorised Service Centre, the warranty shall cease.

Handling in case of a warranty claim

In order to ensure prompt processing of your matter, please observe the following notes:

- Please keep the article number **IAN: 326029_1904** and the purchase receipt as a purchase verification for all inquiries.
- The article number can be found on the rating plate, an engraving, the title page of your instructions (in the bottom left), or as a label on the rear or underside of the device.
- If faulty operation or other defects occur, first contact the Service Centre listed in the following by **telephone** or **email**.
- Then, you are able to send a product reported as defective free of charge to the service address specified to you, including the purchase receipt and the information on the defect and when it occurred.

At www.lidl-service.com you can download this and many other manuals, product videos and installation software.

This QR code takes you directly to the Lidl Service website (www.lidl-service.com), where you can enter the article number (IAN) **326029_1904** to open your user instructions.

Supplier

Please note that the following address is **no service address**. First contact the aforementioned Service Centre.

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