## SHIATSU NECK MASSAGER SSNM 2 E2

(GB)

### SHIATSU NECK MASSAGER

Operation and safety notes

(SI)

### NAPRAVA ZA SHIATSU MASAŽO ZATILJA

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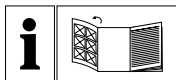
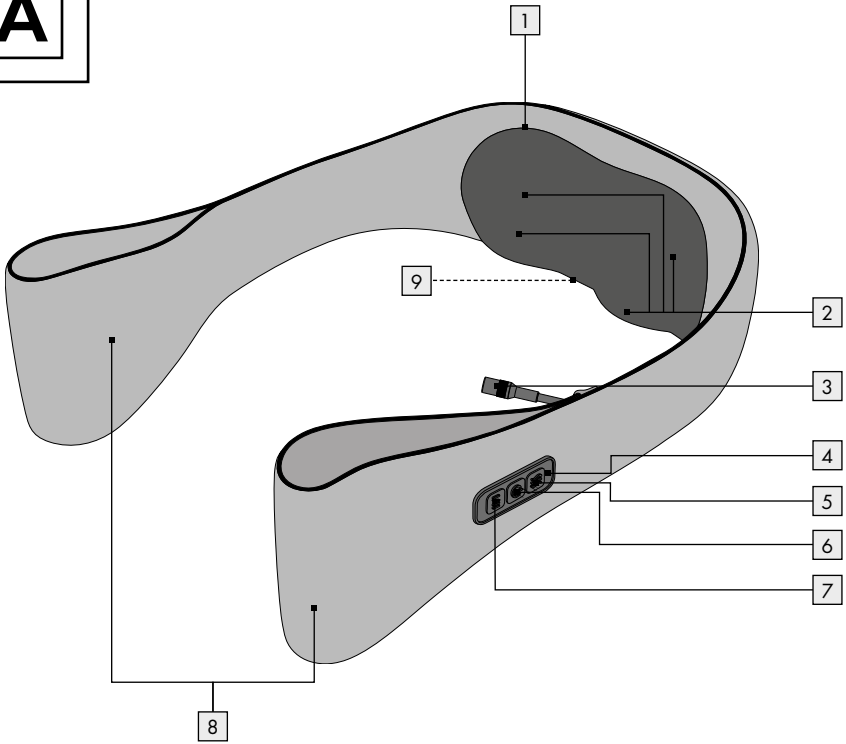
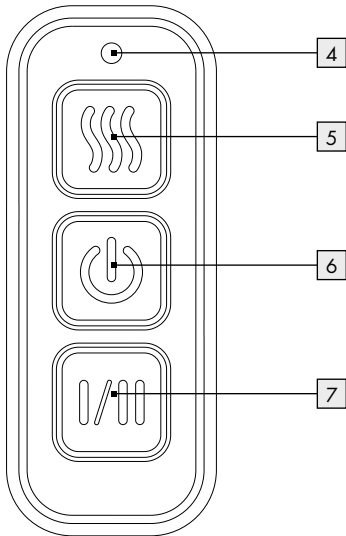
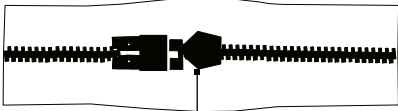
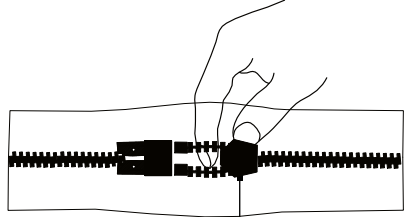
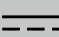
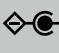
9

9

|                                      |         |
|--------------------------------------|---------|
| <b>List of pictograms used</b> ..... | Page 6  |
| <b>Introduction</b> .....            | Page 6  |
| Intended use .....                   | Page 6  |
| Description of parts .....           | Page 7  |
| Technical data .....                 | Page 7  |
| Scope of delivery .....              | Page 7  |
| <b>Safety advice</b> .....           | Page 7  |
| <b>Before use</b> .....              | Page 9  |
| <b>Operation</b> .....               | Page 9  |
| Warming function .....               | Page 9  |
| Switch off the product .....         | Page 9  |
| <b>Cleaning and care</b> .....       | Page 10 |
| <b>Storage</b> .....                 | Page 10 |
| <b>Disposal</b> .....                | Page 10 |
| <b>Warranty</b> .....                | Page 10 |
| Warranty claim procedure .....       | Page 11 |
| Service .....                        | Page 11 |

| List of pictograms used   |   |   |   |
|---|---|---|---|
|   | Observe caution and safety notes!   |  | Fire hazard!  |
|   | Only use in dry indoor rooms.   | <b>Hz</b>   | Hertz (mains frequency)   |
|   | Alternating current/voltage   |  | Do not use bleach   |
|   | Direct current/voltage  |  | Do not dry in a tumbler dryer   |
|   | The cover is machine-washable.<br>Set the washing machine to a delicate wash cycle at 40 °C |  | Do not iron   |
|   | Short-circuit proof safety transformer  |  | Do not dry clean  |
|   | Polarity  |  | Safety class II   |
|   | SMPS (Switch mode power supply unit)<br>(switching power supply)                            |  | Power   |
|   | Massage function  |  | Heat function   |
| <b>IP20</b>   | Protection against penetration of solid objects (Ø > 12.5 mm)                               |  | Do not use the charger if the pins of the plug are damaged                            |
|    | Danger - Risk of electric shock!  | <b>CE</b>   | CE mark indicates conformity with relevant EU directives applicable for this product. |
|   | Safety information<br>Instructions for use  |   |   |

## Shiatsu neck massager

### ● Introduction

We congratulate you on the purchase of your new product. You have chosen a high quality product. The instructions for use are part of the product. They contain important information concerning safety, use and disposal. Before using the product, please familiarise yourself with all of the safety information and instructions for use. Only use the product as

described and for the specified applications. If you pass the product on to anyone else, please ensure that you also pass on all the documentation with it.

### ● Intended use

The product is only designed for massaging the neck, shoulders, back, abdomen, thighs and lower legs of people with warming function and shiatsu style (shi = finger, astu = pressure in Japanese). The product is only for indoor household use. It is

not intended for medical or commercial purposes and is not a substitute for medical treatment. Any other form of use is not permitted and may result in injury or damage of the product. The manufacturer accepts no liability for any damage caused by improper use.

## ● Description of parts

- 1 Removable cover
- 2 Massage heads
- 3 Socket
- 4 Heat function indicator
- 5 Heat function button
- 6 Power button
- 7 Massage button
- 8 Hand loops
- 9 Zipper head on removable cover

## ● Technical data

### Mains power adapter

| Information   | Value          | Unit |
|---|----------------|------|
| Manufacturer's name or trade mark, commercial registration number and address:  |                |      |
| <b>Xiamen Keli Electronics Co., Ltd</b><br>91350211737865674J<br>No.1, Nantang Road, Guankou Town, Jimei District, Xiamen, 361023 Fujian, P.R. China. |                |      |
| Model identifier  | KL-WE120200-N1 |      |
| Input voltage   | 100-240        | V~   |
| Input AC frequency  | 50/60          | Hz   |
| Output voltage  | 12.0           | V=== |
| Output current  | 2.0            | A    |
| Output power  | 24.0           | W    |
| Average active efficiency   | 87.9           | %    |
| Efficiency at low load (10%)  | 84.9           | %    |
| No-load power consumption   | 0.1            | W    |

| Information       | Value | Unit |
|-------------------|-------|------|
| Protection class  | II    |      |
| Input current     | 0.7   | A    |
| Protection degree | IP20  |      |

### Shiatsu neck massager

|                             |                                  |
|-----------------------------|----------------------------------|
| Input voltage:              | 12 V===                          |
| Rated power:                | 24 W                             |
| Duration of use:            | max. 15 min (automatic shut-off) |
| Off mode power consumption: | 0.2 W                            |

## ● Scope of delivery

Immediately after unpacking, please check the package contents for completeness and if the product and all parts are in good condition. Remove all packaging materials before use.

- 1 Shiatsu neck massager
- 1 Mains power adaptor
- 1 Instruction manual

## Safety advice

Read all safety advices and instructions. Non-observance of the safety advices and instructions may cause electric shock, fire and/or bad injuries.

**KEEP ALL SAFETY NOTICES AND INSTRUCTIONS FOR FUTURE REFERENCE!**

## **⚠ WARNING!**

### **DANGER TO LIFE**

#### **AND RISK OF ACCIDENT FOR INFANTS AND CHILDREN!**

Never leave children unsupervised with the packaging materials. The packaging material poses a suffocation hazard. Children frequently underestimate the dangers. Please always keep the product out of the reach of children.

**⚠ CAUTION!** This product is not a toy for children! Children are not aware of the dangers associated with handling electrical products.

■ This product can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the product in a safe way and understand the hazards involved. Children should not clean or perform user maintenance unsupervised.

**⚠ CAUTION!** The product has a heated surface. Persons insensitive to heat must be careful when using it.

**⚠ WARNING!** Avoid the hazard of life by electric shock!

- Do not insert needles or pointed objects into the product.
- Never immerse the product into water or any liquid.
- Never touch the product, or the mains power adaptor with wet hands or use it in very humid environments.
- Only use the provided mains power adaptor as described in Technical Data section.
- Make sure that the details on the rating label are in conformity with your mains voltage.
- In the event of malfunctions, remove the mains power adaptor from the mains socket.
- Do not open the product or attempt to repair it yourself. Defective products should only be repaired by qualified personal.
- Always remove the mains power adaptor from the mains socket after each use and before clean the product.
- Keep the power cord away from hot surfaces, sharp edges and mechanical forces.
- Never let the power cord become a trip hazard, that no-one can get tangled in it or tread on it.

- Allow the product to cool down before moving it to a different location.

-  **WARNING! HAZARD OF FIRE!** Do not operate the product close to or below curtains, wall cupboards or other inflammable materials.
- Do not use the product on swollen, burnt, inflamed, diseased or injured regions of the skin and body areas. If in doubt, obtain medical advice before using the product.
- Never use the product if you have injuries or pain in the neck, shoulders, back, abdomen, thighs or lower legs.
- No action is needed from users to shift the product between 50 Hz and 60 Hz. The product adapts itself for both 50 Hz and 60 Hz.

## ● Before use

**Note:** Remove all packaging material from the product.

- Connect the connection plug at the wire of the mains power adaptor to the socket [3] at the inner side of the hand loops [8].
- Place the massager around your neck, shoulders, back, abdomen, thighs or lower legs. Align the massage heads [2] on the side of your neck, shoulders, back, abdomen, thighs or lower legs use the hand loops [8] to adjust it. The control buttons must be on the left.

- Plug the mains power adaptor to the mains power socket.
- Ensure the massager heads [2] aligned to your neck, shoulders, back, abdomen, thighs or lower legs.

## ● Operation

- To switch on the massager, press the power button [6] to start rotation of massage heads [2] in low speed.
- To change the rotation speed of the massage heads [2], use the massage button [7] to toggle between low speed rotation/high speed rotation/pause. To stop the massage function, press the power button [6].  
**Note:** Never use the product for more than 15 minutes. A longer massage may lead to tension caused by overstimulation of the muscles.  
**Note:** If you found yourself painful or uncomfortable, stop the massage immediately.

## ● Warming function


- To complement the massage with the warming function, press the heat function button [5] to turn on the heat function. The red heat function indicator [4] and the red light of the massage heads [2] switch on.

## ● Switch off the product


- To switch off the product, press the power button [6], the heat function indicator [4] and the light of the massage heads [2] will go off, the massage heads [2] stop rotating. The product switches off automatically after 15 minutes.
- To turn off the product completely, disconnect the mains power adaptor from the power supply.

## ● Cleaning and care
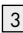
### ■ **⚠ WARNING! DANGER TO LIFE BY ELECTRIC SHOCK!**

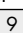Disconnect the mains power adaptor from the mains socket and unplug connection plug from the socket  before cleaning. Otherwise there is a risk of electric shock.

- **⚠ WARNING!** Never immerse the product, or mains power adaptor into water.
- Ensure the product is completely dry before connecting the mains power adaptor.
- Do not use abrasive, aggressive or chemical cleaning agents to clean the product. It could damage the product.
- Never allow any liquid to enter the product.
- Use a damp lint-free cloth to clean the product. Use a little mild detergent on the cloth if necessary.

-  To clean the removable cover , use finger to pull out the zipper head on the removable cover 

(see Fig. C) and open the zip completely. Clean the removable cover  in accordance with the cleaning symbols on the label. It is machine-washable. Set the washing machine to a delicate wash cycle at 40 °C.

**Note:** Failure to observe the following warnings may damage the product.

-  Do not bleach. Do not use any cleaning agents containing bleach (such as heavy duty detergents).
-  Do not tumble dry.
-  Do not iron.
-  Do not dry clean.

## ● Storage

- Let the product completely cool down before storage.
- Do not place other objects on top of the product.

- If the product is not used for a certain period of time, store the product in a clean and dry location.

## ● Disposal

The packaging is made entirely of recyclable materials, which you may dispose of at local recycling facilities.

Observe the marking of the packaging materials for waste separation, which are marked with abbreviations (a) and numbers (b) with following meaning: 1–7: plastics/20–22: paper and fibre-board/80–98: composite materials.

Contact your local refuse disposal authority for more details of how to dispose of your worn-out product.

To help protect the environment, please dispose of the product properly when it has reached the end of its useful life and not in the household waste. Information on collection points and their opening hours can be obtained from your local authority.

## ● Warranty

The product has been manufactured to strict quality guidelines and meticulously examined before delivery. In the event of material or manufacturing defects you have legal rights against the retailer of this product. Your legal rights are not limited in any way by our warranty detailed below.

The warranty for this product is 3 years from the date of purchase. The warranty period begins on the date of purchase. Keep the original sales receipt in a safe location as this document is required as proof of purchase.

Any damage or defects already present at the time of purchase must be reported without delay after unpacking the product.

Should the product show any fault in materials or manufacture within 3 years from the date of purchase, we will repair or replace it – at our choice – free of charge to you. The warranty period is not extended as a result of a claim being granted. This also applies to replaced and repaired parts.

This warranty becomes void if the product has been damaged, or used or maintained improperly.

The warranty covers material or manufacturing defects. This warranty does not cover product parts subject to normal wear and tear, thus considered consumables (e.g. batteries, rechargeable batteries, tubes, cartridges), nor damage to fragile parts, e.g. switches or glass parts.

## ● Warranty claim procedure

To ensure quick processing of your claim, observe the following instructions:

Make sure to have the original sales receipt and the item number (IAN 487603\_2501) available as proof of purchase.

You can find the item number on the rating plate, an engraving on the product, on the front page of the instruction manual (bottom left), or as a sticker on the rear or bottom of the product.

If functional or other defects occur, contact the service department listed below either by telephone or by e-mail.

Once the product has been recorded as defective you can return it free of charge to the service address that will be provided to you. Ensure to enclose the proof of purchase (sales receipt) and a short, written description outlining the details of the defect and when it occurred.

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**OWIM GmbH & Co. KG**