**ADVANCED ELECTRIC TOOTHBRUSH /
AKKU-ZAHNBÜRSTE „ADVANCED“ /
BROSSE À DENTS ÉLECTRIQUE SANS FIL « ADVANCED »
NZAOR 600 A2**

For EU market:
HOYER Handel GmbH
Kühnehöfe 12
22761 Hamburg
GERMANY

Imported for GB market by:
Lidl Great Britain Ltd
Lidl House · 14 Kingston Road
Surbiton · KT5 9NU

Stand der Informationen · Status of information · Version des informations · Stand van de informatie
Stan informacii · Stav informací · Stav informácií · Estado de las informaciones · Informationernes stand:

02/2022 ID: NZAOR 600 A2_22_V1.2

DE AT CH Kurzanleitung

AKKU-ZAHNBÜRSTE „ADVANCED“

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IAN: 377265_2110

Short manual

This document is an abridged print version of the complete user instructions. Scanning the QR code will take you directly to the Lidl service page (www.lidl-service.com) and by entering the article number (IAN) **377265_2110** you can view and download the complete user instructions.

WARNING! Observe the complete user instructions and the safety information to avoid personal injury and damage to property. The quick reference guide is a component of this product. Familiarise yourself with all user and safety information before using the product. Keep the quick reference guide in a safe place and hand over all documents when passing the product on to third parties.

Items supplied

- 1 toothbrush **10**
- 2 toothbrush heads **3**
- 1 battery charger **12**
- 1 travel case **14**
- 1 short manual

Intended purpose

The electric toothbrush is intended solely for tooth cleaning and oral hygiene in humans. The device is designed for private, domestic use. This device must not be used for commercial purposes.

Overview

- 1** Brush head
- 2** Indicator bristles
- 3** Toothbrush head (replaceable)
- 4** Marking ring
- 5** Peg of the toothbrush
- 6** LED:
 - is lit white after switching on
 - is lit red when pressure too high
 - flashes red when the transport lock is switched on
 - flashes 3x red when switching the transport lock on/off
- 7**  On/off switch + transport lock for on and off switch
- 8** Function displays (illuminated):
 -  Clean
 -  Sensitive
 -  Massage
- 9** Charging indicator light
- 10** Toothbrush
- 11** Peg of the battery charger
- 12** Battery charger
- 13** Mains power cable and mains plug
- 14** Travel case
- 15** Pressure cap

Instructions for safe operation

- ⊙ This device can be used by children and people with restricted physical, sensory or intellectual abilities or people without adequate experience and/or understanding, if they are supervised or instructed on how to use this device and if they are aware of the resulting hazards. Cleaning and user maintenance must not be performed by children unless they are supervised.
- ⊙ Children must not play with the device.
-  ⊙ The battery charger is connected to the power voltage, and must therefore not be used in the shower, in the bath, or over washbasins filled with water.
- ⊙ The mains power cable of this device cannot be replaced. If the cable should become damaged, the device itself must be disposed of.
- ⊙ The battery must be charged only with the original accessory (battery charger) supplied.
- ⊙ This device contains a battery that cannot be replaced.
- ⊙ Protect the battery against mechanical damages. **Fire hazard!**
- ⊙ Do not expose the device to direct sunshine or heat. The ambient temperature should not be lower than 10 °C or higher than 40 °C.
- ⊙ Should the electrolyte solution leak out of the battery, avoid contact with eyes, mucous membranes and skin. Rinse immediately with plenty of clean water and consult a doctor. The electrolyte solution can cause irritation.
- ⊙ The device is provided with a lithium-ion battery.
 - The battery cannot be removed!
 - This device must not be opened!
 - The device with the installed battery must be disposed of properly.
- ⊙ The toothbrush, the battery charger, the mains power cable and the mains plug must never be immersed in water or any other liquids.
- ⊙ Never touch the device with wet hands.
- ⊙ If the battery charger has fallen into water, pull out the mains plug immediately, and only then remove the battery charger from the water.

If this happens, do not use the battery charger again until it has been checked by a specialist workshop.
- ⊙ The installation of a residual current device with a nominal trip current not exceeding 30 mA in the building wiring offers additional protection. Contact a skilled electrician in this regard.
- ⊙ Do not start using the device if the device, the battery charger or the mains power cable shows visible damage, or if the device has previously been dropped.
- ⊙ In order to avoid any hazards, do not make modifications to the device. Repairs must be carried out by a specialist workshop or the Service Centre.

Charging the battery

- Never set up the battery charger **12** anywhere it could fall into water.

NOTES:

- Fully charge the battery before using it for the first time.
- It takes approx. 16 hours to fully charge the battery.

1. Plug the mains plug **13** of the battery charger **12** into a suitable wall socket. The wall socket must still be easily accessible after connection.
2. Place the switched-off toothbrush **10** onto the peg of the battery charger **11**. The charging indicator light **9** flashes white throughout the entire charging process. Lit segments of the charging indicator light show the progress of the charging process. When the battery is fully charged, first all segments of the charging indicator light **9** are lit, and then they go out.

Battery signals: charging indicator light **9**

Toothbrush **10** on the battery charger **12**:

- All segments of the charging indicator light **9** are lit white and then go out: the battery is fully charged.
- One or two segments of the charging indicator light **9** briefly light up white: the battery is being charged.

Toothbrush **10** switched on:

- The charging indicator light **9** flashes red: the battery is weak. Charge the battery.

Toothbrush heads

The items supplied include 2 toothbrush heads **3**.

The various marking rings **4** allow you to assign the toothbrush heads **3** to different family members.

- Take a toothbrush head **3** and slide it onto the peg of the toothbrush **5**.
- In order to remove, carefully pull the toothbrush head **3** off upwards.
- The indicator bristles **2** on the brush head **1** slowly lose their colour during use. The decolouration shows when the toothbrush heads **3** should be replaced.

NOTE: we recommend that you replace the toothbrush heads **3** every 3 months or after infections.

Operation

CAUTION:

- ⊙ Users could possibly experience mild bleeding of the gums during the first days of electric toothbrush usage. Should the bleeding continue over several days, a dentist should be consulted.
 - ⊙ If you have had dental surgery in the last two months, consult a dentist before using the toothbrush.
 - ⊙ Exert only light pressure when brushing to prevent damage to the gums and to increase the service life of the bristles.
-

Modes

The toothbrush **10** has 3 different modes. Depending on the desired mode, press the on/off switch **Ⓛ 7** repeatedly when switched on, until the corresponding symbol is lit:

Clean mode

- for cleaning teeth with high pulsation

Sensitive mode

- gentle brushing for cleaning teeth with sensitive gums

NOTE: if you are not used to using an electric toothbrush and you have comparatively sensitive gums, you can use this mode for the first 1 to 2 weeks, and then switch to Clean mode later.

Massage mode

- pulsating brushing to massage your gums

Brushing teeth

1. Remove the toothbrush **10** from the battery charger **12**.
 2. Switch on the toothbrush using the on/off switch **Ⓛ 7**.
 3. Press on/off switch **Ⓛ 7** repeatedly until the function display **8** for the desired mode is lit.
-

NOTE: switching over to another mode must be done quickly after the last button press of the on/off switch **Ⓛ 7** (**within approx. 1 second**). Pressing the on/off switch **Ⓛ 7** approx. 2 seconds after the last button press will switch the toothbrush off.

4. Switch off the toothbrush using the on/off switch **Ⓛ 7**.
5. Moisten the bristles, apply some toothpaste, and only switch on the toothbrush **10** again when the brush head **1** is already in your mouth.
6. Slowly guide the brush head **1** from tooth to tooth without exerting too much pressure. Carefully clean all tooth surfaces, the inner surfaces of the teeth and the tooth edges.

Figure A: the recommended brushing time is at least 2 minutes, divided up into 30-second periods for each area of the mouth (bottom right, bottom left, top right, top left). The toothbrush **10** indicates these intervals by two brief interruptions.

NOTE: if the brush head **1** is pressed against the teeth too hard, the LED **6** lights up red. You should then reduce the pressure.

7. The toothbrush signals the end of a cycle of 2 minutes with three short interruptions. Switch off the toothbrush using the on/off switch **Ⓛ 7**. Remove the brush head **1** from the mouth only when the toothbrush **10** has been switched off.

Cleaning

For additional information on cleaning and looking after the device and accessories, please see the online user instructions.

- To ensure the hygienic functionality of the toothbrush heads **3** they must be rinsed thoroughly after every use and allowed to dry.

Disposal

Packaging and device must be disposed of according to the environmental regulations of your country.

Trouble-shooting


If your device fails to function as required, please read the corresponding chapter in the online user instructions.

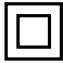Technical specifications

Model:	NZAOR 600 A2
Toothbrush	
Protection class:	IPX4 (protected against water splashes)
Battery:	1 x 3.7 V  Li ion 600 mAh 2.22 Wh
Charging time:	16 hours
Operating time:	approx. 90 min
Battery charger	
Mains voltage:	100 - 240 V ~ 50/60 Hz
Power rating:	1.4 W Standby: 0.5 W
Protection class:	II 
Protection class:	IPX7 (protected against the effects of temporary immersion)
General	
Dimensions:	approx. 38 x 190 x 29 mm (W x H x D)
Weight:	approx. 150 g
LEDs:	RG 0
Ambient conditions:	10 °C to 40 °C humidity 40 to 85% relative

NOTE: no action by the user is required to switch the product between 50 and 60 Hz. The product adjusts itself for both 50 and 60 Hz.

Symbols used

	Double insulation
	Geprüfte Sicherheit (certified safety): devices must comply with the generally acknowledged rules of technology and the German Product Safety Act (Produktsicherheitsgesetz - ProdSG). (Not for UK)
	With the CE marking, HOYER Handel GmbH declares the conformity with EU guidelines.
	This symbol reminds you to dispose of the packaging in an environmentally friendly manner.
	Recyclable materials are marked with the recycling symbol (3 arrows). The material can be specified by the recycling number in the centre (here: 21) and/or an abbreviation (here: PAP).
	Direct current
	AC voltage
	It is a recyclable product which is subject to an extended producer responsibility and for which the rules of waste separation apply.
	With the UKCA marking, HOYER Handel GmbH declares that the product is in conformity with the UKCA requirements.
	The symbol identifies the product manufacturer.

Subject to technical modifications.

Service Centre

If you require service, please contact our Service Centre:

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Med CE-mærkning erklærer HOYER Handel GmbH EU-konformiteten.
Med UKCA-mærkningen erklærer HOYER Handel GmbH UK-konformiteten.