SHIATSU NECK MASSAGE CUSHION / SHIATSU BACK MASSAGE CUSHION SSNR 12 B2

(GB)

SHIATSU NECK MASSAGE CUSHION / SHIATSU BACK MASSAGE CUSHION

Operation and safety notes

(SE)

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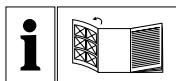
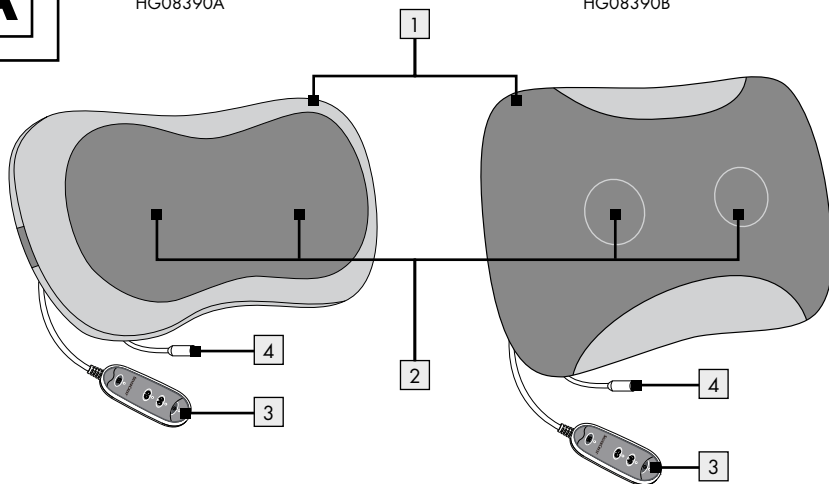
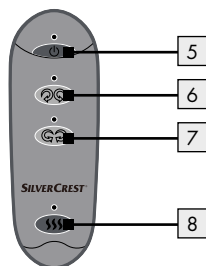
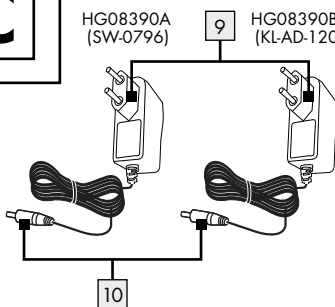
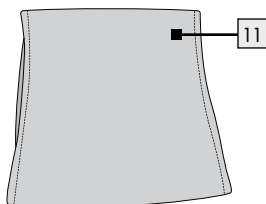
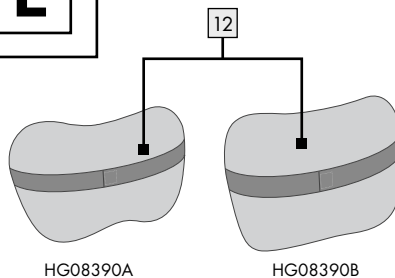
HG08390B

**B****C**HG08390A
(SW-0796)HG08390B
(KL-AD-120100)**D****E**

List of pictograms used	Page 6
Introduction	Page 6
Intended use	Page 6
Description of parts	Page 7
Technical data	Page 7
Scope of delivery	Page 8
Safety advice	Page 8
Before use	Page 9
Operation	Page 9
Warming function	Page 10
Switch off the product	Page 10
Cleaning and care	Page 10
Storage	Page 10
Disposal	Page 10
Warranty	Page 11
Warranty claim procedure	Page 11
Service	Page 11

List of pictograms used			
	Observe caution and safety notes!		Danger - Risk of electric shock!
	Only use in dry indoor rooms.		Fire hazard!
	Alternating current/voltage	Hz	Hertz (mains frequency)
	Direct current/voltage		Do not use bleach
	The cover is machine-washable. Set the washing machine to a delicate wash cycle at 30 °C		Do not dry in a tumbler dryer
	WARNING: Do not use this product near water.		Do not iron
	Short-circuit proof safety transformer		Do not dry clean
	Polarity		Safety class II
	SMPS (Switch mode power supply unit) (switching power supply)		Outward rotation
	Inward rotation		Light and heat function
IP20	Protection against contact with the fingers	t_a	Highest rated ambient temperature
	CE mark indicates conformity with relevant EU directives applicable for this product.		

Shiatsu Neck Massage Cushion/ Shiatsu Back Massage Cushion SSNR 12 B2

● Introduction

We congratulate you on the purchase of your new product. You have chosen a high quality product. The instructions for use are part of the product. They contain important information concerning safety, use and disposal. Before using the product, please familiarise yourself with all of the safety information

and instructions for use. Only use the product as described and for the specified applications. If you pass the product on to anyone else, please ensure that you also pass on all the documentation with it.

● Intended use

The product is only designed for massaging the neck and/or back and legs of people with warming function and shiatsu style (shi = finger, astu = pressure in Japanese). The product is only for indoor

household use. It is not intended for medical or commercial purposes and is not a substitute for medical treatment. Any other form of use is not permitted and may result in injury or damage of the product. The manufacturer accepts no liability for any damage caused by improper use.

Protection class: II/□
 Input current: 0.7 A
 Highest rated ambient temperature: 40 °C
 Ingress protection: IP20
 Plug-in mains adapter (Model identifier: SW-0796)
 TÜV Rheinland/GS certified

● Description of parts

- 1 Massage cushion
- 2 Massage heads
- 3 Control handset
- 4 Socket
- 5 Power button
- 6 Outward rotation button
- 7 Inward rotation button
- 8 Light and heat function button
- 9 Mains power adaptor
- 10 Connection plug
- 11 Cover
- 12 Velcro strap

● Technical data

Mains power adaptor 9 for HG08390A

Manufacturer's name or trade mark, commercial registration number and address:

Xiamen Keli
 Electronics Co., Ltd.,
 91350211737865674J,
 No. 19, XingLin North
 3rd Road, Jimei District,
 Xiamen, Fujian, China

Model identifier: SW-0796
 Input voltage: 100-240V~
 Input AC frequency: 50/60 Hz
 Output voltage: 12.0V===
 Output current: 1.2 A
 Output power: 14.4 W
 Average active efficiency: 84.0 %
 Efficiency at low load (10 %): 73.8 %
 No-load power consumption: 0.10 W

Mains power adaptor 9 for HG08390B

Manufacturer's name or trade mark, commercial registration number and address:

Xiamen Keli
 Electronics Co., Ltd.,
 91350211737865674J,
 No. 19, XingLin North
 3rd Road, Jimei District,
 Xiamen, Fujian, China

Model identifier: KL-AD-120100
 Input voltage: 100-240V~
 Input AC frequency: 50/60 Hz
 Output voltage: 12.0V===
 Output current: 1.0 A
 Output power: 12.0 W
 Average active efficiency: 83.0 %
 Efficiency at low load (10 %): 73.8 %
 No-load power consumption: 0.10 W
 Protection class: II/□
 Input current: 0.7 A
 Highest rated ambient temperature: 40 °C
 Ingress protection: IP20
 Plug-in mains adapter (Model identifier: KL-AD-120100)
 TÜV Rheinland/GS certified

Shiatsu Neck Massager Cushion/ Back Massager Cushion

Input voltage: 12.0V===
 Rated power: 14.4 W (HG08390A)
 12.0 W (HG08390B)
 Polarity: 
 Duration of use: max. 15 min (automatic shut-off)

● Scope of delivery

Immediately after unpacking, please check the package contents for completeness and if the product and all parts are in good condition. Remove all packaging materials before use.

- 1 Shiatsu Neck Massage Cushion/Back Massage Cushion
- 1 Mains power adaptor
- 1 Cover
- 1 Instruction Manual

Safety advice

Read all safety advices and instructions. Non-observance of the safety advices and instructions may cause electric shock, fire and/or bad injuries.

KEEP ALL SAFETY NOTICES AND INSTRUCTIONS FOR FUTURE REFERENCE!

-  **⚠ WARNING!**
DANGER TO LIFE AND RISK OF ACCIDENT FOR INFANTS AND CHILDREN! Never leave children unsupervised with the packaging materials. The packaging material poses a suffocation hazard. Children frequently underestimate the dangers. Please always keep the product out of the reach of children.

⚠ CAUTION! This product is not a toy for children! Children are not aware of the dangers associated with handling electrical products.

- This product can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the product in a safe way and understand the hazards involved. Children shall not play with the product. Do not allow children to clean or perform user maintenance without supervision.

⚠ CAUTION! The product has a heated surface. Persons insensitive to heat must be careful when using it.

⚠ WARNING! Avoid the hazard of life by electric shock:

- Do not insert needles or pointed objects into the massage cushion **[1]**.
- Never immerse the product into water or any liquid.
- Never touch the product, or the mains power adaptor **[9]** with wet hands or use it in very humid environments.

- Only use the provided mains power adaptor [9] as described in Technical Data section.
- Make sure that the details on the rating label are in conformity with your mains voltage.
- In the event of malfunctions, remove the mains power adaptor [9] from the mains socket.
- Do not open the product or attempt to repair it yourself. Defective products should only be repaired by qualified personnel.
- Always remove the mains power adaptor [9] from the mains socket after each use and before clean the product.
- Keep the power cord away from hot surfaces, sharp edges and mechanical forces.
- Never let the power cord becomes a trip hazard, that no-one can get tangled in it or tread on it.
- Allow the product to cool down before moving it to a different location.

WARNING! **HAZARD OF FIRE!**

Do not operate the product close to or below curtains, wall cupboards or other inflammable materials.

- Do not use the product on swollen, burnt, inflamed, diseased

or injured regions of the skin and body areas. If in doubt, obtain medical advice before using the product.

- Never use the product if you have injuries or pain in the neck/back or spine.
- No action is needed from users to shift the product between 50 Hz and 60 Hz. The product adapts itself for both 50 Hz and 60 Hz.

● **Before use**

Note: Remove all packaging material from the product.

- Always keep the cover [11] on the massage cushion [1] before use. Only the cover [11] is intended to be washed according to "Cleaning and care" section.
- Place the massage cushion [1] on a high chair and fasten it to the back of the high chair with the Velcro strap [12].
- Connect connection plug [10] to the socket [4] at the back of the massage cushion [1].
- Plug the mains power adaptor [9] to the mains power socket.
- Sit on the chair and ensure the massage cushion [1] is in the desired position.

● **Operation**

- Switch on the massage cushion [1] by pressing the power button [5] on the control handset [3], the blue light indicator switches on.
- Press the outward rotation button [6] to start rotation of massage heads [2]. The blue light indicator switches on above the button.
- To reverse the rotation direction, press the inward rotation button [7]. The other blue light indicator switches on instead.

Note: Never use the product for more than 15 minutes. A longer massage may lead to tension caused by overstimulation of the muscles.

● Warming function

- To complement the massage with the warming function, press the light and heat function button [8] to turn on the heat function. The red light indicator switches on.

Note: The warming function can only work when the massage heads [2] are in operation.

● Switch off the product

- To switch off the product, press the power button [5], all indicators will go off and the massage heads [2] stop rotating. The product switches off automatically after 15 minutes.
- To turn off completely, disconnect the mains power adaptor [9] from the mains socket.

● Cleaning and care

-  **⚠ WARNING! DANGER TO LIFE BY ELECTRIC SHOCK!**
Disconnect the mains power adaptor [9] from the mains socket and unplug connection plug [10] from the socket [4] before cleaning.
- **⚠ WARNING!** Never immerse the product, or mains power adaptor [9] into water.
- Ensure the massage cushion [1] is completely dry before connecting the mains power adaptor [9].
- Do not use abrasive, aggressive or chemical cleaning agents to clean the massage cushion [1]. It could damage the product.
- Never allow any liquid to enter the massage cushion [1].
- Use a damp lint-free cloth to clean the product. Use a little mild detergent on the cloth if necessary.

-  Clean the cover [11] in accordance with the cleaning symbols on the label. It is machine-washable. Set the washing machine to a delicate wash cycle at 30 °C.

Note: Failure to observe the following warnings may damage the cover [11].

-  Do not bleach. Do not use any cleaning agents containing bleach (such as heavy duty detergents).
-  Do not tumble dry.
-  Do not iron.
-  Do not dry clean.

● Storage

- Let the massage cushion [1] completely cool down before storage.
- Do not place other objects on top of the massage cushion [1].
- If the product is not used for a certain period of time, store the product in a clean and dry location.

● Disposal

The packaging is made entirely of recyclable materials, which you may dispose of at local recycling facilities.

Observe the marking of the packaging materials for waste separation, which are marked with abbreviations (a) and numbers (b) with following meaning:
1-7: plastics/20-22: paper and fibre-board/80-98: composite materials.

The product and packaging materials are recyclable, dispose of it separately for better waste treatment. The Triman logo is valid in France only.

Contact your local refuse disposal authority for more details of how to dispose of your worn-out product.

To help protect the environment, please dispose of the product properly when it has reached the end of its useful life and not in the household waste. Information on collection points and their opening hours can be obtained from your local authority.

● **Warranty**

The product has been manufactured to strict quality guidelines and meticulously examined before delivery. In the event of product defects you have legal rights against the retailer of this product. Your legal rights are not limited in any way by our warranty detailed below.

The warranty for this product is 3 years from the date of purchase. The warranty period begins on the date of purchase. Please keep the original sales receipt in a safe location. This document is required as your proof of purchase.

Should this product show any fault in materials or manufacture within 3 years from the date of purchase, we will repair or replace it – at our choice – free of charge to you. This warranty becomes void if the product has been damaged, or used or maintained improperly.

The warranty applies to defects in material or manufacture. This warranty does not cover product parts subject to normal wear, thus possibly considered consumables (e.g. batteries) or for damage to fragile parts, e.g. switches, rechargeable batteries or glass parts.

● **Warranty claim procedure**

To ensure quick processing of your case, please observe the following instructions:

Please have the till receipt and the item number (e.g. IAN 123456_7890) available as proof of purchase.

You will find the item number on the rating plate, an engraving on the front page of the instructions for use (bottom left), or as a sticker on the rear or bottom of the product.

If functional or other defects occur, please contact the service department listed either by telephone or by e-mail.

You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

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ja osoite: Xiamen Keli Electronics
Co., Ltd.,
91350211737865674J,
No. 19, XingLin North
3rd Road, District Jimei,
Xiamen, Fujian, China

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ja osoite: Xiamen Keli Electronics
Co., Ltd.,
91350211737865674J,
No. 19, XingLin North 3rd
Road, District Jimei,
Xiamen, Fujian, China

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OWIM GmbH & Co. KG