SILVER CREST®

MASSAGE GUN / PISTOLET DE MASSAGE SMGB D3

(GB) (IE)

MASSAGE GUN

Operating instructions

(FR) (BE)

PISTOLET DE MASSAGE

Mode d'emploi

(CZ)

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(SK)

MASÁŽNY PRÍSTROJ

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(DK)

MASSAGEPISTOL

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(HU)

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(DE) (AT) (CH)

MASSAGE GUN

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(NL) (BE)

MASSAGE GUN

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(PL)

MASSAGE GUN - PISTOLET DO MASAŻU

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(ES)

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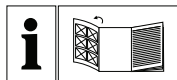
(IT)

PISTOLA MASSAGGIANTE

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IAN 444831_2307

(DE) (FR) (BE) (NL)
(CZ) (PL) (SK)

GB IE

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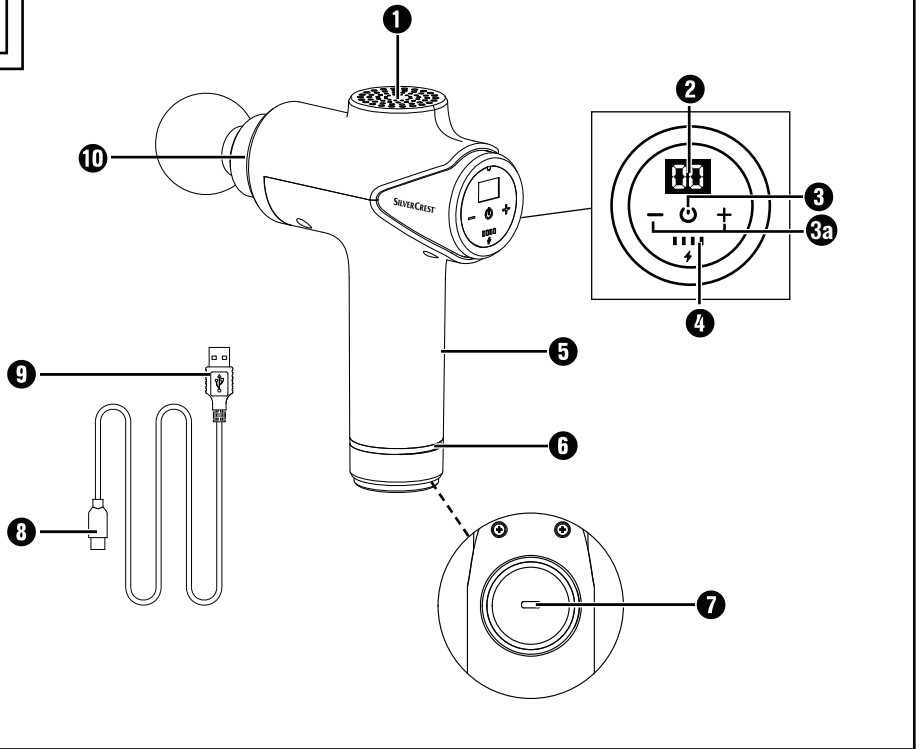
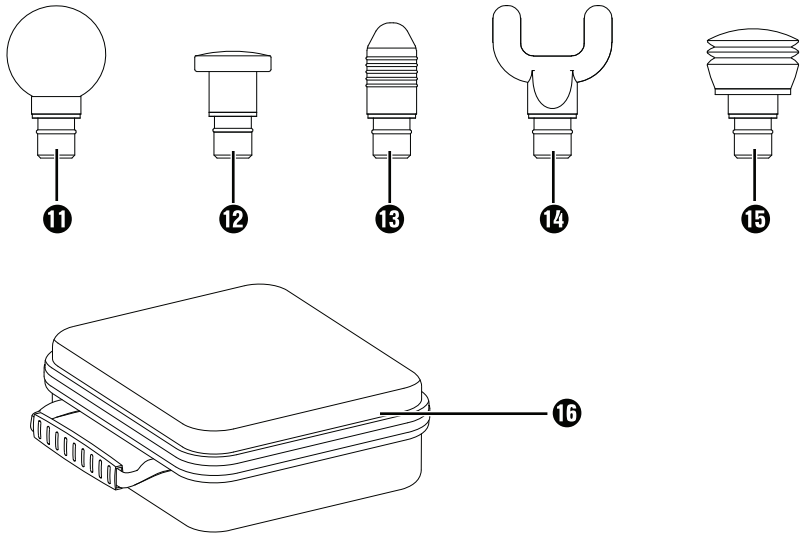
A**B**

Contents

Introduction	2
Intended use	2
Warning symbols used	2
Notes on trademarks	3
Safety information	3
Safety guidelines for rechargeable batteries	6
Package contents	7
Parts description	7
Before first use	8
Charging the integrated battery	8
Using the appliance	9
Selecting/fitting attachments	10
Examples of use	11
Switching the appliance on/off and selecting intensity	12
Intensity levels	13
Cleaning and care	13
Storage	14
Troubleshooting	14
Disposal	14
Disposal of the appliance	15
Disposal of the packaging	15
Technical specifications	15
Kompernass Handels GmbH warranty	16
Service	17
Importer	17

Introduction

Congratulations on the purchase of your new appliance.

You have selected a high-quality product. The operating instructions are part of this product. They contain important information about safety, usage and disposal. Before using the product, please familiarise yourself with all operating and safety instructions. Use the product only as described and for the specified range of applications. Please also pass these operating instructions on to any future owner.

Intended use

This appliance and its attachments are intended exclusively for massaging the muscles of the neck, shoulders, arms, back (excluding the spine), thighs, calves and buttocks. Only use the attachments intended for this purpose. Never use the appliance directly on bones.

This appliance is intended solely for personal, private use. It is not to be used for medical or commercial applications and is not a substitute for medical treatment. It is not intended for use in medical/therapeutic or commercial fields. Any other use or use beyond that specified is considered improper use.

Warning symbols used

The following warnings are used in these operating instructions (where applicable):

DANGER

A warning notice at this hazard level indicates an imminently hazardous situation.

If the hazardous situation is not avoided, it could result in death or serious physical injury.

- ▶ Follow the instructions in this warning to prevent a fatality or serious personal injury.

WARNING

A warning at this hazard level indicates a potentially hazardous situation.

Failure to avoid this hazardous situation could result in injuries.

- ▶ Follow the instructions in this warning notice to prevent injuries.

CAUTION

A warning at this hazard level indicates a risk of property damage.

Failure to avoid this situation could result in property damage.

- ▶ Follow the instructions on this warning label to prevent property damage.

NOTE

- ▶ A note provides additional information which will assist you in using the appliance.

Notes on trademarks

- USB® is a registered trademark of USB Implementers Forum, Inc. All other names and products may be trademarks or registered trademarks of their respective owner

Safety information

To ensure safe operation of the appliance, follow the safety guidelines set out below.

DANGER! RISK OF ELECTROCUTION!

- ▶ Make sure that the charging cable does not become wet or damp. Route the cable in such a way that it cannot be caught or damaged.
- ▶ Keep the appliance and the charging cable away from hot surfaces.
- ▶ Do not perform any repairs to the appliance. All repairs must be performed by our Customer Service Unit or by a qualified specialist technician.
- ▶ Never open the appliance housing. There are no parts inside the appliance which require user maintenance. In addition, this will invalidate the warranty.
- ▶ Switch off the appliance after each use and before cleaning it.
- ▶ Use the appliance in dry interior rooms only; never use it outside.
- ▶ Never immerse the appliance in water or other liquids!
- ▶ Do not use the appliance while taking a bath or shower. Do not store or keep the appliance in places from where it can fall into the bath tub, sink or other containers filled with water or other liquids.

⚠ DANGER! RISK OF ELECTROCUTION!

- ▶ Never reach for the appliance if it has fallen into water.
- ▶ Do not use the appliance if the power cable is damaged, if the appliance is not working properly, if the attachments show signs of damage, if it has sustained a fall, if it has fallen into water or if it is otherwise damaged.

⚠ WARNING! RISK OF INJURY!

- ▶ This appliance may be used by children aged 8 years and above and by persons with limited physical, sensory or mental capabilities or lack of experience and knowledge, provided that they are under supervision or have been told how to use the appliance safely and are aware of the potential risks.
- ▶ Cleaning and user maintenance tasks must not be carried out by children unless they are supervised.
- ▶ Children must not be allowed to play with the appliance.
- ▶ Never leave children unattended with the packaging material. Keep children away from the packaging. There is a risk of suffocation!
- ▶ People with pacemakers, artificial joints or electronic implants should consult their doctor before using the appliance.
- ▶ Do not use the appliance around your heart, face or genital area.
- ▶ Do not use the appliance if you suffer from one or more of the following diseases or medical conditions: circulatory disorders, varicose veins, open wounds, bruises, cracked skin, phlebitis, thrombosis, diabetes or slipped disc.
- ▶ Never use the appliance on children, animals, during pregnancy, on your head, throat or other parts of the body that are exhibiting swelling, burns, inflammation, wounds or other sensitive areas.

⚠ WARNING! RISK OF INJURY!

- ▶ Do not use the appliance under blankets or pillows. Never use the appliance in the vicinity of petrol or other inflammable materials.
- ▶ Do not use the appliance while sleeping or after taking substances that may reduce your ability to react quickly (e.g. painkillers or alcohol).
- ▶ Never use the appliance directly on bones, e.g. on your spine, shins, the tops of your feet or similar places on your body.
- ▶ Use the appliance only on dry and clean parts of your body.
- ▶ Stop using the appliance immediately if you find the massage unpleasant or if you experience pain.
- ▶ Uninterrupted application of the appliance should not last any longer than 15 minutes. After use, switch off the appliance for at least 15 minutes and allow it to cool down.
- ▶ Do not insert your fingers or any other objects into the opening for the attachments. Keep long hair away from moving parts.
- ▶ Use the appliance only for its intended purpose.

CAUTION! PROPERTY DAMAGE!

Use the appliance only indoors.

- ▶ Use the appliance only with the supplied accessories. Make sure that the attachments can always move freely. Use only the supplied charging cable to charge the appliance.
- ▶ If you are using a USB power supply, the mains socket must always be easily accessible so that the USB power supply can be quickly disconnected from the socket in a hazardous situation. You should also observe the operating instructions for the USB power supply.

CAUTION! PROPERTY DAMAGE!

- ▶ Protect the appliance and attachments from shocks, dust, chemicals or electromagnetic fields. Keep the appliance away from heat sources (radiators, ovens, etc.) and do not expose it to extreme fluctuations of temperature.
- ▶ Do not use any aggressive cleaning agents or solvents. They could damage the surface.
- ▶ Keep the ventilation slits free of dust and do not cover them.

Safety guidelines for rechargeable batteries

The product is equipped or supplied with a lithium-ion battery. The battery in this appliance cannot be replaced.

⚠ DANGER! RISK OF EXPLOSION!

- ▶ Do not open or short-circuit the battery. This could cause the battery to overheat and explode. There is a risk of fire!
- ▶ Never throw the battery into a fire or water.
- ▶ Do not subject the battery to mechanical stress.
- ▶ Never leave the appliance unattended during charging.

⚠ WARNING! RISK OF INJURY DUE TO BATTERY LEAKAGE!

- ▶ High temperatures can affect the battery. Avoid direct sunlight and proximity to heat sources such as radiators.
- ▶ Under extreme conditions, the battery cells may leak. Avoid contact with skin, eyes and mucous membranes. In the event of accidental contact, rinse the affected areas with plenty of clean water and contact a doctor.
- ▶ Leaking or damaged batteries can cause chemical burns to the skin. Always wear protective gloves if you have to touch the batteries.

Package contents

The appliance is supplied with the following components as standard:

- Massage gun
- 5 attachments (1 ball, 1 conical, 1 U-shaped, 1 flat and 1 air cushion head)
- Storage bag
- USB charging cable
- Operating instructions

- 1) Remove all parts of the appliance and the operating instructions from the box.
- 2) Remove all packaging materials and any stickers from the appliance.
- 3) Remove the protective film from the display **2**.

NOTE

- ▶ Check contents for completeness and damage immediately after unpacking. Contact Customer Service if necessary.

Parts description

Figure A (see the fold-out page):

- 1** Ventilation slits
- 2** Display / intensity display
- 3** ON/OFF button 
- 3a** **+** and **-** buttons for setting the intensity level
- 4** Battery level indicator
- 5** Handle
- 6** Ring display
- 7** Connection socket
- 8** USB plug (type C)
- 9** USB charging cable (USB type A to USB type C)
- 10** Socket for attachments

Figure B:

- 11** Ball head
- 12** Flat head
- 13** Conical head
- 14** U-shaped head
- 15** Air cushion head
- 16** Storage bag

Before first use

- Charge the battery as described in the section **Charging the integrated battery**.

Charging the integrated battery

⚠ DANGER! ELECTRIC SHOCK!

- ▶ Only charge the appliance indoors in dry locations away from water. For example, do not charge the appliance next to or above a sink filled with water.

NOTE

- ▶ The USB charging cable ⑨ and connection socket ⑦ are used only to charge the integrated battery. They are not suitable for data transmission.
- ▶ Only charge the appliance with a USB power adapter with protection class II that is approved for use with household appliances and has an output voltage of 5 V $\overline{=}$ and an output current of max. 2 A.
- ▶ Use only the supplied USB charging cable ⑨ to charge the integrated battery. Disconnect the USB charging cable ⑨ from the appliance when charging is complete.
- ▶ The operating time with a fully charged battery is about 2.5 to 4 hours, depending on the selected intensity level.
- ▶ The appliance cannot be operated during charging. If you switch on the appliance during charging, CH will light up in the display ②.

- 1) Make sure that the appliance is switched off.
- 2) Insert the USB plug ⑥ into the USB connection socket ⑦ on the underside of the handle ⑤ (fig. 1):

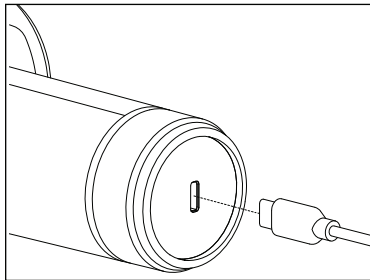
Fig. 1

- 3) Connect the USB charging cable ⑨ to a suitable mains adapter and plug the mains adapter into a mains socket.
- ✓ When the ring display ⑥ lights up green, the appliance is fully charged.

The ring display **6** has the following meanings:

Ring display 6	Meaning
Lights up red	Low battery
Flashes green	Appliance is charging
Lights up green	During charging: battery fully charged During operation: battery level sufficient

- 4) Disconnect the appliance from the power supply and unplug the USB charging cable **9** from the appliance.

The battery can be recharged at any time and at any battery level.

Each of the 4 LEDs of the battery level indicator **4** indicates 25% of a full battery charge.

- If, for example, 3 of 4 LEDs light up, this means that the battery is still 75% charged.
- If only 1 LED is still on, the battery has only 25% capacity.
If the battery charge drops below 5%, the last LED starts flashing.
Charge the appliance immediately to avoid damage to the battery.
- When the appliance is charging, the respective LED flashes until the corresponding battery charge is reached. After that, it lights up continuously.
When the battery is fully charged, all 4 LEDs light up continuously.

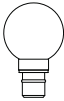
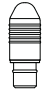
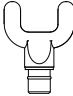
Using the appliance

You can use this appliance to give yourself or others a soothing and effective trigger point massage. Trigger point massages can relax or stimulate the muscles to relieve tension, pain and fatigue. The massage gun provides a powerful, intensive massage for the thighs, calves, shoulders, arms and the back region (excluding the spine). Depending on the extent of the tension in the part of the body to be treated, a maximum treatment time of approx. 15 minutes is recommended.

A massage should always feel pleasant. A slight reddening of the treated areas and even slight bruising are a normal side effect, as the massage specifically promotes blood circulation. However, if you experience excessive skin reactions, discontinue use and consult a doctor.

Selecting/fitting attachments

- 1) Select the appropriate attachment for the desired treatment and body part as follows:

Attachment	Suitable for:
 <p>11</p>	Uniform and gentle massage to relax smaller muscle groups, e.g. upper arms, shoulders, waist and legs.
 <p>12</p>	Large areas all over the body, such as buttocks, thighs and back (excluding spine).
 <p>13</p>	Targeted and selective massage of deep muscle tissue, e.g. soles of the feet, palms of the hands, muscles of the shoulder blades.
 <p>14</p>	Uniform, 2-point stimulation of deep muscles in the lower back (except spine), neck, Achilles tendon.
 <p>15</p>	Gentle, soothing tapping massage of sensitive parts of the body, e.g. neck, shoulders, wrists and feet.

- 2) Insert the desired attachment into the socket 10 and press it as far as it will go (Fig. 2):

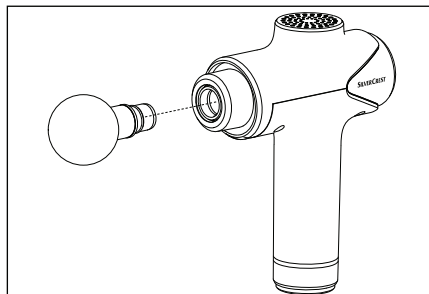
Fig. 2

- 3) To change the attachment, carefully pull it out of the socket 10 horizontally, turning it slightly if necessary.

Examples of use

NOTE

- Do not apply the attachments directly to a bone, only to muscle tissue.

Examples of application on arms (fig. 3), thighs (fig. 4), sole of the foot (fig. 5), Achilles tendon (fig. 6) and neck (fig. 7).

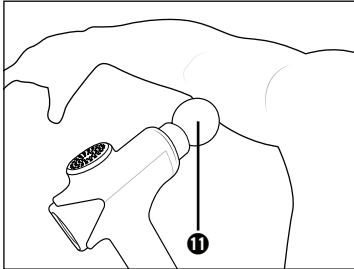
Fig. 3

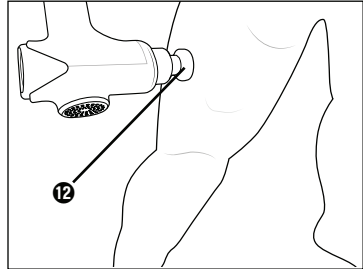
Fig. 4

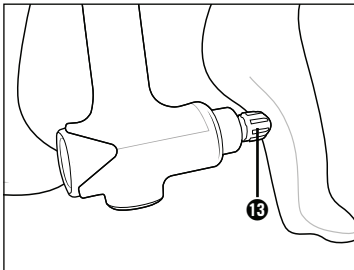
Fig. 5

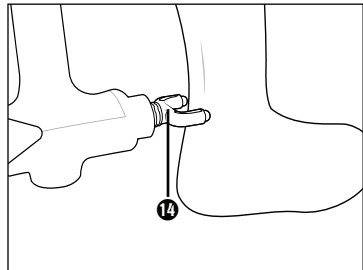
Fig. 6

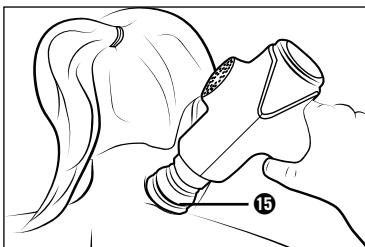
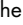
Fig. 7

NOTE

- The appliance may only be used with massage oil or gel when the silicone ball head 11 is attached. All other attachments are unsuitable for the use of oil or gel.

Switching the appliance on/off and selecting intensity

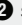
- 1) Switch on the appliance by pressing and holding the ON/OFF button  3 for about 2 seconds.
 - ✓ The ring display  6 lights up. The display  2 shows 00.
- 2) Press the +  3a button to start the appliance.
 - ✓ The display  2 shows the intensity level 0 1. The appliance starts on intensity level 0 1.
- 3) To increase the intensity level, briefly press the +  3a button. To decrease the intensity level, briefly press the -  3a button.
 - ✓ The intensity changes by one level each time you press it. This is shown on the display  2.
- 4) Massage the desired body part by gently moving the attachment back and forth over the muscle area.

NOTE

- ▶ Do not apply too much pressure, otherwise the appliance will automatically switch off after a few seconds to protect the motor from a possible overload.
If this safety device has been triggered, E will appear on the display  2. You need to restart the appliance.

- 5) If you no longer wish to use the appliance, switch it off by pressing and holding the ON/OFF button  3 for about 2 seconds.

NOTE

- ▶ The appliance is equipped with an automatic switch-off to prevent overheating. The appliance switches off automatically after 15 minutes of continuous operation. Allow the appliance to cool down for at least 15 minutes before using it again.
- ▶ If you do not select an intensity level after switching on (display  2 shows 00) or if you do not use the appliance for more than 5 minutes, it switches off automatically after 5 minutes.

Intensity levels

The intensity of the massage varies depending on the level selected, as follows:

Level	Intensity
00	Standby
01	1800 rpm
02	1935 rpm
03	2065 rpm
04	2200 rpm
05	2335 rpm
06	2470 rpm
07	2600 rpm
08	2735 rpm
09	2870 rpm
10	3000 rpm

Cleaning and care

DANGER! ELECTRIC SHOCK!

- ▶ Do not clean the appliance while it is charging. Switch off the appliance before cleaning.
- ▶ Ensure that no moisture can penetrate the device during cleaning.

ATTENTION! PROPERTY DAMAGE!

- ▶ Do not use abrasive or aggressive cleaning materials. These can damage the surface of the appliance!
- 1) Make sure that the appliance is switched off.
 - 2) Make sure that the appliance is not connected to a mains adapter. Unplug the USB charging cable  from the appliance if it is plugged in.
 - 3) Wipe off the appliance using a damp cloth. Do not use any solvents or cleansers which can damage plastic. If necessary, add a little mild detergent to the cloth and wipe off afterwards with fresh water. Make sure that no water or moisture can get into the appliance!

- 4) To hygienically clean the attachments, wash them with water and a mild detergent and then rinse them with clear water.
- 5) Allow everything to dry thoroughly before storing the appliance or putting it back into operation.

Storage

Store the cleaned appliance in the storage bag **16** in a dust-free and dry location out of direct sunlight. Charge the battery at regular intervals to ensure a long battery life.

Troubleshooting

Fault	Cause	Solution
The appliance is not working.	The battery is empty.	Charge the appliance.
	The appliance is defective.	Contact the Customer Service department.
The display 2 indicates E .	The motor overload protection has been triggered.	Restart the appliance. Do not press too hard while massaging.
The display 2 indicates CH .	The appliance is charging and cannot be used.	Wait until the battery is fully charged.

If the malfunction cannot be corrected using one of the suggestions above, or if you have found some other kind of malfunction, contact Customer Service.

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The product, its packaging and the operating instructions are recyclable. They are subject to an extended manufacturer responsibility and will be collected separately.

Disposal of the appliance

The adjacent symbol of a crossed-out dustbin means that this appliance is subject to Directive 2012/19/EU. This directive states that this appliance may not be disposed of in the normal household waste at the end of its useful life, but must be taken to specially set-up collection locations, recycling depots or disposal companies.

The disposal is free of charge for the user. Protect the environment and dispose of this appliance properly.

If your old appliance has stored any personal data, you are responsible for deleting it yourself before returning it.

If it is possible to do so without destroying the old appliance, remove the old batteries or rechargeable batteries before returning the appliance for disposal and take them to a separate collection point. In the case of permanently installed rechargeable batteries, you must indicate during disposal that the appliance contains a battery.

The integrated rechargeable battery in this appliance cannot be removed for disposal.

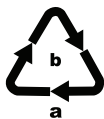
Your local community or municipal authorities can provide information on how to dispose of the worn-out product.

Disposal of the packaging

The packaging materials have been selected for their environmental friendliness and ease of disposal and are therefore recyclable. Dispose of packaging materials that are no longer needed in accordance with applicable local regulations.

Dispose of the packaging in an environmentally friendly manner. Note the labelling on the packaging and separate the packaging material components for disposal, if necessary. The packaging material is labelled with abbreviations (a) and numbers (b) with the following meanings: 1–7: plastics, 20–22: paper and cardboard, 80–98: composites.

Technical specifications

Appliance	
Input voltage	5.0 V \equiv (DC)
Input current	2.0 A
Integrated battery (Li-Ion)	3000 mAh / 7.2 V \equiv (DC) / 21.6 Wh
Protection class	III / \blacklozenge (protection through low voltage)
	detachable supply unit

Kompernass Handels GmbH warranty

Dear Customer,

This appliance has a 3-year warranty valid from the date of purchase. If this product has any faults, you, the buyer, have certain statutory rights. Your statutory rights are not restricted in any way by the warranty described below.

Warranty conditions

The warranty period starts on the date of purchase. Please keep your receipt in a safe place. This will be required as proof of purchase.

If any material or manufacturing fault occurs within three years of the date of purchase of the product, we will either repair or replace the product for you or refund the purchase price (at our discretion). This warranty service requires that you present the defective appliance and the proof of purchase (receipt) within the three-year warranty period, along with a brief written description of the fault and of when it occurred.

If the defect is covered by the warranty, your product will either be repaired or replaced by us. The repair or replacement of a product does not signify the beginning of a new warranty period.

Warranty period and statutory claims for defects

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking. Repairs carried out after expiry of the warranty period shall be subject to a fee.

Scope of the warranty

This appliance has been manufactured in accordance with strict quality guidelines and inspected meticulously prior to delivery.

The warranty covers material faults or production faults. The warranty does not extend to product parts subject to normal wear and tear or to fragile parts which could be considered as consumable parts such as switches, batteries or parts made of glass.

The warranty does not apply if the product has been damaged, improperly used or improperly maintained. The directions in the operating instructions for the product regarding proper use of the product are to be strictly followed. Uses and actions that are discouraged in the operating instructions or which are warned against must be avoided.

This product is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications / repairs which have not been carried out by one of our authorised Service centres.

Warranty claim procedure

To ensure quick processing of your case, please observe the following instructions:

- Please have the till receipt and the item number (IAN) 444831_2307 available as proof of purchase.
- You will find the item number on the type plate on the product, an engraving on the product, on the front page of the operating instructions (below left) or on the sticker on the rear or bottom of the product.
- If functional or other defects occur, please contact the service department listed either by **telephone** or by **e-mail**.
- You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

You can download these instructions along with many other manuals, product videos and installation software at www.lidl-service.com.

This QR code will take you directly to the Lidl service page (www.lidl-service.com) where you can open your operating instructions by entering the item number (IAN) 444831_2307.

Service

GB Service Great Britain

Tel.: 0800 404 7657

E-Mail: kompernass@lidl.co.uk

IE Service Ireland

Tel.: 1800 101010

E-Mail: kompernass@lidl.ie

IAN 444831_2307

Importer

Please note that the following address is not the service address. Please use the service address provided in the operating instructions.

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