WATER JET FLOSSER / MUNDDUSCHE HYDROPULSEUR NMD 3.7 D8

(GB) (IE)

WATER JET FLOSSER

Operating instructions

(FR) (CH)

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Mode d'emploi

(CZ)

ÚSTNÍ SPRCHA

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(SK)

ÚSTNA SPRCHA

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(DK)

MUNDSKYLLER

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(HU)

SZÁJZUHANY

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(DE) (AT) (CH)

MUNDDUSCHE

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(NL) (BE)

MONDDOUCHE

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(PL)

IRYGATOR DO JAMY USTNEJ

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(ES)

IRRIGADOR BUCAL

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
(IT)

IDROPULSORE

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
IAN 435653_2304

(DE) (FR) (BE) (NL)
(CZ) (PL) (SK)

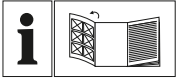GB/IE	Operating instructions	Page	1
DE/AT/CH	Bedienungsanleitung	Seite	19
FR/BE	Mode d'emploi	Page	39
NL/BE	Gebruiksaanwijzing	Pagina	61
CZ	Návod k obsluze	Strana	79
PL	Instrukcja obsługi	Strona	97
SK	Návod na obsluhu	Strana	117
ES	Instrucciones de uso	Página	135
DK	Betjeningsvejledning	Side	153
IT	Istruzioni per l'uso	Pagina	171
HU	Használati utasítás	Oldal	189

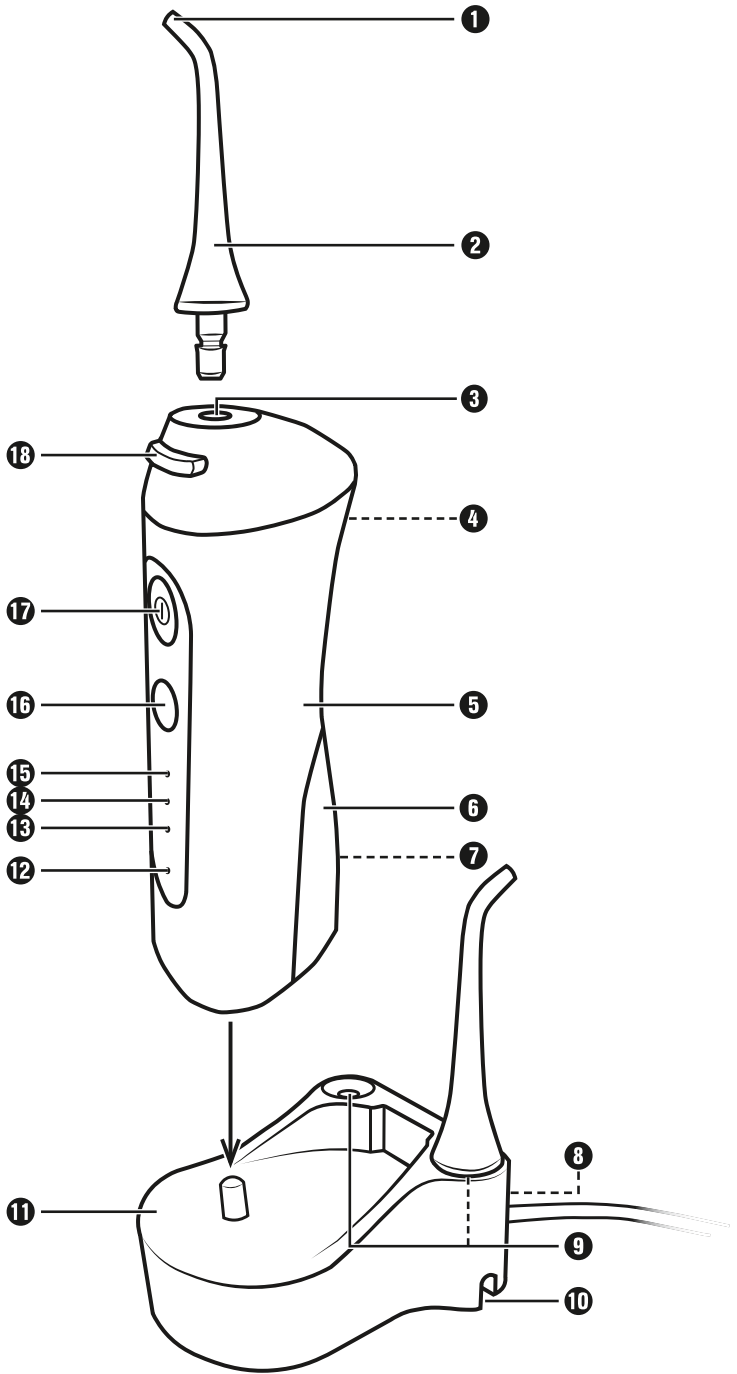Contents

Information about these operating instructions.	2
Intended use	2
Warnings and symbols used	2
Safety information	3
Technical data	6
Package contents	6
Appliance description	7
Setting up the appliance / wall mounting.	8
Charging the appliance	9
Filling the water tank	9
Fitting/removing attachments	10
Before first use	10
Power levels	11
Operation	11
After use	13
Cleaning and maintenance	13
Storage	14
Disposal	14
Disposal of the appliance	14
Disposal of the packaging.	15
Kompass Handels GmbH warranty.	15
Service.	17
Importer.	17
Ordering spare parts	17

Information about these operating instructions

Congratulations on the purchase of your new appliance. You have selected a high-quality product.

The operating instructions are part of this product. They contain important information about safety, usage and disposal. Before using the product, please familiarise yourself with all operating and safety instructions. Use the product only as described and for the specified range of applications. Please also pass on all documents to any future owner.


Intended use

This appliance is intended exclusively for personal dental care at home. It is intended solely for use in private households. This appliance is not intended for use in industrial, commercial or medical environments. It is not suitable for animal care.

It is not intended for any other purpose nor for use beyond the scope described.

Warnings and symbols used

The following warnings and symbols are used in these operating instructions, on the packaging and on the device (where applicable):

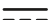	DANGER! A warning with this symbol and the signal word "DANGER", indicates an imminently hazardous situation which will result in death or serious injury if not avoided.
	WARNING! A warning with this symbol and the signal word "WARNING", indicates a potentially hazardous situation which could result in death or serious injury if not avoided.
	CAUTION! A warning with this symbol and the signal word "CAUTION", indicates a potentially hazardous situation which could result in a minor or moderate injury if not avoided.
	ATTENTION! A warning with this symbol and the signal word "ATTENTION", indicates a potential situation which could result in property damage if not avoided.
	Note: A note provides additional information which makes handling the device easier for you.
	Alternating current (AC)/alternating voltage
	Direct current (DC)/direct voltage
	Removable power supply

Safety information

⚠ DANGER! ELECTRIC SHOCK!

- ▶ Connect the appliance only to a properly installed and easily accessible mains power socket supplying a mains power voltage of 100–240 V ~, 50/60 Hz.
- ▶ Ensure that the power cable is not damaged. Keep it out of hot areas and route it in such a way that it cannot be damaged.
- ▶ Before use, check the appliance for visible external damage. Never use the appliance if the cable, plug or housing are damaged.
- ▶ In the event of malfunctions or obvious damage, please contact our Customer Service department.
- ▶ The power cable of this device cannot be replaced. If the cable is damaged, the device must be scrapped.
- ▶ Have all repairs carried out by a specialist workshop. Under no circumstances should you open the appliance yourself. Repairs that are not carried out by a specialist workshop could lead to physical injuries.
- ▶ Do not use the appliance in the vicinity of a bathtub/shower or a sink filled with water. The proximity of moisture presents a danger even when the appliance is switched off.
- ▶ To avoid damage, do not allow liquids to penetrate the appliance and do not submerge it in water.
- ▶ Ensure that no foreign objects enter the appliance.

- ▶ Ensure that the power cable can never become wet or damp during operation. Route the power cable in such a way that it cannot become trapped or damaged in any other way. Risk of electric shock!
- ▶ Never touch the charging station if it has fallen into water. Remove the plug from the wall socket first! Risk of electric shock!
- ▶  Do not submerge the handpiece or charging station in water or other liquids.
- ▶ Never touch the power cable or the plug with wet or damp hands.
- ▶ Do not use the power cable with an extension cable; connect the power cable directly to a power socket.

 WARNING! RISK OF INJURY!

- ▶ This appliance may be used by children from the age of 8 and by persons with reduced physical, sensory or mental abilities or lack of experience and/or knowledge if they are supervised.
- ▶ Do not allow children to play with the appliance.
- ▶ Always store the appliance indoors to avoid accidents. After use, store the appliance in a dry place away from direct sunlight.
- ▶ Under extreme conditions, leakage can occur from the battery cells and electrolyte solution can leak out. If the liquid comes into contact with skin or eyes, immediately rinse the affected area with plenty of clean water and seek medical attention.

- ▶ Switch off the appliance before you put it down or change accessories.
- ▶ Lay the power cable so that nobody can step on it or trip over it.
- ▶ Always pull out the plug before cleaning the appliance. Do not pull on the power cable. Always grip the power plug to disconnect the appliance from the power supply.


⚠ ATTENTION! PROPERTY DAMAGE!

- ▶ Charge the appliance only with the supplied charging station.
- ▶ Use only accessories recommended by the manufacturer.
- ▶ Do not make any unauthorised modifications or alterations to the appliance.
- ▶ Never put down the appliance next to or on radiators, ovens or other heated appliances or surfaces.
- ▶ The appliance may only be used with cold or lukewarm tap water up to a maximum of 40°C. If necessary, add a few drops of mouthwash. The water tank should never be completely filled with additives such as mouthwash or mouth rinse. These can damage the appliance.
- ▶ This appliance contains a battery that cannot be replaced.
- ▶ The product is equipped or supplied with a lithium-ion battery. Do not throw the rechargeable battery into a fire and do not subject it to high temperatures. There is a risk of an explosion!

- ▶ Do not charge or use the appliance outdoors.

Technical data

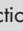Charging station	
Power supply	100–240 V ~, 50/60 Hz
Power consumption	3.5 W
Protection class	II /  (double insulation)
Protection type	IPX7 (protection against temporary submersion)
Battery	
Capacity	1200 mAh
Battery (lithium ion)	3.7 V — /4.44 Wh
Handpiece	
Water tank volume	130 ml
Protection type	IPX7 (protection against temporary submersion)

Note

- ▶ No user action is required to switch the product between 50 and 60 Hz. The product switches automatically to either 50 or 60 Hz.

Package contents

The appliance is supplied with the following components as standard:

- Handpiece with water tank
- Charging station
- 2 attachments
- 2 screws
- 2 wall plugs (∅ 5 mm)
- Operating instructions

Note

- ▶ Check the package for completeness and signs of visible damage.
- ▶ If the delivery is incomplete or damage has occurred as a result of defective packaging or during transport, contact the Customer Service hotline (see **Service** section).

Appliance description

(See fold-out page for illustrations)

- ❶ Nozzle
- ❷ Attachment
- ❸ Attachment socket
- ❹ Air supply openings
- ❺ Handpiece
- ❻ Water tank
- ❼ Filling opening, water tank/water tank lid
- ❽ Hangers
- ❾ Holder for attachments
- ❿ Recess for power cable
- ⓫ Charging station
- ⓬ Charge indicator 
- ⓭ Indicator light **soft**
- ⓮ Indicator light **normal**
- ⓯ Indicator light **jet**
- ⓰ Power level selection button
- ⓱ On/Off button 
- ⓲ Release button

Setting up the appliance / wall mounting

i Note

- ▶ The charging station **11** may only be operated free-standing if used in dry rooms. Place the charging station **11** on a horizontal, level and stable surface.
- ▶ In damp rooms, especially in the vicinity of water (bath), the charging station **11** must be firmly fixed to the wall to protect it from falling into the water.
- ▶ Before wall mounting, make sure that there are no cables or pipes in the area where you want to hang the charging station **11**.

- 1) Mark two drill holes on the wall 3.3 cm apart and drill the marked holes with a 5 mm bit. Ensure that there is a mains power socket in the immediate vicinity.
- 2) Push the supplied wall plugs into the drill holes and screw the screws into the wall plugs. Take care that the screws still protrude about 3 mm out of the holes.
- 3) Feed the power cable through one of the recesses **10** on the side of the charging station **11**, depending on which side the power socket is located.
- 4) Push the charging station **11** onto the screws so that the screws protrude into the hangers **8** at the back of the charging station **11**. Then press the charging station **11** down so that the screws slide into the rails of the hangers **8** and the charging station **11** sits firmly on the screws.

Charging the appliance

i Note

- ▶ Before starting to use the appliance, check that it is in fault-free condition. You must not use the appliance if it has any faults.
- ▶ Charge the appliance only when there is no water in the water tank **6**.
- ▶ In the following cases, the appliance must be charged for about 15 hours:
 - Before first use.
 - If the battery has been drained during use.
 - If you have not used the appliance for more than 3 months.
- ▶ A fully charged battery is enough for about 20 minutes of use.

- 1) Plug the mains plug into the mains power socket.
- 2) Place the handpiece **5** on the charging station **11** so that it slides onto the pin on the charging station **11**.
- 3) Allow the handpiece **5** to charge for around 15 hours before using it for the first time. The charge indicator  **12** on the front of the handpiece **5** lights up red during the charging process. When the battery is fully charged, the charge indicator  **12** lights up green.
- 4) Place the handpiece **5** on the charging station **11** after each use so that full power is always available.

Filling the water tank

! ATTENTION!

- ▶ Fill the water tank **6** with pure tap water. If necessary, add a few drops of mouthwash. The water tank **6** should never be completely filled with additives such as mouthwash or mouth rinse.
- ▶ Use only cold or lukewarm water with a maximum temperature of 40°C. Do not fill the water tank **6** with hot water.

You can fill/top up the water tank **6** at any time without removing it from the appliance:

- 1) Take the handpiece **5** with the attached empty water tank **6** from the charging station **11**. Hold the switched-off handpiece **5** horizontally.
- 2) Open the water tank lid **7** and fill the tank with water through the filling opening **7**.
- 3) Close the water tank lid **7**.

Fitting/removing attachments

i Note

- ▶ Make certain that the handpiece **5** is switched off before fitting or removing the attachments **2**.
- ▶ The appliance comes with two attachments **2** with different coloured marking rings for different users.
- ▶ If used daily, replace the attachments **2** after no more than six months.
- ▶ Use only the attachments **2** recommended by the manufacturer. Information on ordering replacement attachments **2** can be found in the section **Ordering replacement parts**. The attachment nozzles of Panasonic ©* water jet flossers (nozzle model EW0950) can also be used for this appliance (*Panasonic is a registered trademark of Panasonic Corporation).

- ◆ To fit an attachment **2** to the handpiece **5**, insert it into the attachment socket **3**. Press the attachment **2** down until it completely engages in the attachment socket **3**.
- ◆ To remove an attachment **2** from the handpiece **5**, press the release button **18** and pull the attachment **2** upwards off the handpiece **5**.
- ◆ Attachments **2** that are not in use can be placed in one of the holders **9** on the charging station **11**.

Before first use

- 1) Rinse the attachments **2** in clean water. One attachment **2** can be placed into one of the holders **9** on the charging station **11**. Fit the other attachment **2** onto the handpiece **5**.
- 2) Pull the water tank **6** downwards off the handpiece **5** and rinse the water tank **6** in clean water.
- 3) Push the water tank **6** back onto the handpiece **5** until it clicks into place completely. Ensure that the intake tube is in the water tank **6**.
- 4) Fill the water tank **6** as described in the section **Filling the water tank**.
- 5) Point the nozzle **1** of the attachment **2** downwards into a washbasin. When doing so, make sure that you do not tilt the handpiece **5** too far forwards. Otherwise, the intake tube will no longer reach water or water could leak out of the air supply openings **4**.
- 6) Press the On/Off button **1** **17**. The handpiece **5** starts pumping water and discharging it via the nozzle **1**. Let the handpiece **5** run until the water tank **6** is empty. This flushes all internal lines.
- 7) Press the On/Off button **1** **17** to switch off the handpiece **5**.

Power levels

This water jet flosser has three different power levels:

- **jet**: powerful water jet, indicator light **jet 15** lights up.
For cleaning left over food from your interdental spaces.
- **normal**: water jet mixed with air, indicator light **normal 14** lights up.
For cleaning left over food from your interdental spaces, for massaging your gums, for rinsing between teeth and gums.
- **soft**: water jet mixed with a lot of air, indicator light **soft 13** lights up.
For gently massaging your gums.

Note

- ▶ If the appliance was completely discharged during the last use, it always starts on the power level **jet**.
- ◆ Press the power level selection button **16** repeatedly until the desired power level is set.
- ◆ Press the On/Off button **17**. The appliance starts with the set power level.
- ◆ To switch off the water jet flosser, press the On/Off button **17** again.

Operation

CAUTION!

- ▶ The unfamiliar stimulation may cause some gum bleeding at first, which is not unusual. In this case, use a less intensive power level, such as **soft** or **normal**. If the bleeding gums last longer than two weeks, you should contact your dentist.
- ▶ Consult your dentist before using the appliance if you suffer from severe periodontitis, have injuries or abscesses in the mucous membrane of the oral cavity or underwent dental surgery within the last two months.

Note

- ▶ Use of the appliance is not a replacement for daily tooth cleaning with a toothbrush and toothpaste. Use the appliance only as a supplement to daily tooth cleaning with toothbrush and toothpaste.
- ▶ Make sure that you do not cover up the air supply openings **4** on the back of the handpiece **5** with your fingers during use.
- ▶ If you hold the handpiece **5** at an angle during use, water may leak out of the air supply openings **4**.
- ▶ The handpiece **5** switches off automatically after approx. 50 seconds, even if there is still water in the water tank **6**.

- 1) Fill the water tank **6**.
- 2) Place an attachment **2** on the handpiece **5**.
- 3) Press the power level selection button **16** to set the desired power level.
- 4) Place the attachment **2** in your mouth to prevent spray. Bend over the washbasin and hold the handpiece **5** as upright as possible (see Fig. 1). Direct the nozzle **1** of the attachment **2** towards your teeth/interdental spaces.

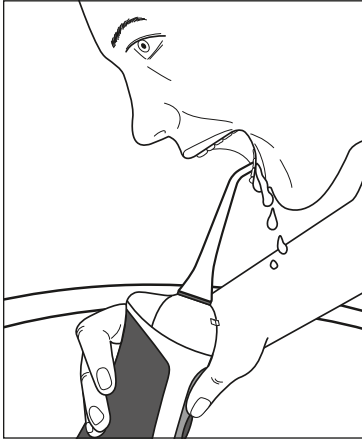
Fig. 1

- 5) Press the On/Off button **1** **17** to start the handpiece **5**. Close your mouth a little to prevent water from spraying around.
- 6) To clean the interdental spaces, slowly guide the nozzle **1** of the attachment **2** from interdental space to interdental space on the inner and outer sides of your teeth. To massage your gums, point the nozzle **1** at the gums. Move the water jet slowly over your gums.
The length of the process will depend on how much cleaning the interdental spaces require. Continue cleaning until your teeth feel clean.
When massaging your gums, use the water jet flosser for as long as you find the massage comfortable.

i Note

- ▶ The water in the water tank **6** is enough for about 30 seconds of use.
- ▶ If the water in the water tank **6** is used up in the meantime, stop the handpiece **5**. Remove it from your mouth and refill the water tank **6** through the filling opening **7**.

- 7) When the treatment is finished, press the On/Off button **①** **⑰** to switch off the handpiece **⑤**.
- 8) Only then should you take the attachment **②** out of your mouth and rinse it under clean, running water.

After use

i Note

- ▶ Always empty the water tank **⑥** and all tubes after every use to prevent the build-up of bacteria.

When the handpiece **⑤** is switched off and you have finished the treatment:

- 1) Open the water tank lid **⑦** with the attachment **②** still fitted and pour out the rest of the water.
- 2) Close the water tank lid **⑦**.
- 3) Point the nozzle **①** of the attachment **②** downwards into a washbasin. Switch the handpiece **⑤** on for a couple of seconds to flush any remaining water from the water tank **⑥**.
- 4) Remove the attachment **②** from the handpiece **⑤** and place it in the holder **⑨** on the charging station **⑩**.
- 5) Place the handpiece **⑤** on the charging station **⑩**.

Cleaning and maintenance

⚠ DANGER!

- ▶ Always unplug the charging station **⑩** from the power socket before cleaning. There is a risk of electric shock!

- ▶ Never immerse the handpiece **⑤** or charging station **⑩** in water or other liquids!

! ATTENTION!

- ▶ Never open the housing of the charging station **⑩** or the handpiece **⑤**! They cannot be resealed in a watertight manner.
- ▶ Do not use abrasive, aggressive or chemical cleaning agents. They could damage the outer surfaces of the appliance.

- ◆ Clean the charging station **⑩** only with a soft, lightly moistened cloth. Do not use detergents or solvents. They can damage the plastic surfaces.
- ◆ Remove the attachment **②** from the handpiece **⑤**. After use, clean the attachments **②** thoroughly under running water.

- ◆ To clean the handpiece **5**, remove the water tank **6** and wipe off the handpiece **5** with a damp cloth. For more stubborn soiling, add a little detergent to the cloth and rinse it off afterwards with a cloth dampened with clean water. Dry the handpiece **5** afterwards.
- ◆ Rinse the water tank **6** thoroughly with warm water, including from the outside. Empty the water tank **6** completely into a sink after cleaning. If you are not using the appliance for over a week, leave the water tank **6** open to ensure it can dry fully.

Storage

- ◆ Disconnect the appliance from the power supply by pulling the plug of the switched-off appliance from the socket.
- ◆ Remove any water from the appliance and allow the appliance to dry completely (see section **Cleaning and maintenance**).
- ◆ Store the appliance in a dry and dust-free location out of direct sunlight.

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The product, its packaging and the operating instructions are recyclable. They are subject to an extended manufacturer responsibility and will be collected separately.

Disposal of the appliance

The adjacent symbol of a crossed-out dustbin means that this appliance is subject to Directive 2012/19/EU. This directive states that this appliance may not be disposed of in the normal household waste at the end of its useful life, but must be taken to specially set-up collection locations, recycling depots or disposal companies.

The disposal is free of charge for the user. Protect the environment and dispose of this appliance properly.

If your old appliance has stored any personal data, you are responsible for deleting it yourself before returning it.

If it is possible to do so without destroying the old appliance, remove the old batteries or rechargeable batteries before returning the appliance for disposal and take them to a separate collection point. In the case of permanently installed rechargeable batteries, you must indicate during disposal that the appliance contains a battery.

The integrated rechargeable battery in this appliance cannot be removed for disposal.


Your local community or municipal authorities can provide information on how to dispose of the worn-out product.

Disposal of the packaging

The packaging materials have been selected for their environmental friendliness and ease of disposal and are therefore recyclable. Dispose of packaging materials that are no longer needed in accordance with applicable local regulations.

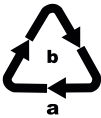Dispose of the packaging in an environmentally friendly manner. Note the labelling on the packaging and separate the packaging material components for disposal, if necessary. The packaging material is labelled with abbreviations (a) and numbers (b) with the following meanings: 1–7: plastics, 20–22: paper and cardboard, 80–98: composites.

Kompernass Handels GmbH warranty

Dear Customer,

This appliance has a 3-year warranty valid from the date of purchase. If this product has any faults, you, the buyer, have certain statutory rights. Your statutory rights are not restricted in any way by the warranty described below.

Warranty conditions

The warranty period starts on the date of purchase. Please keep your receipt in a safe place. This will be required as proof of purchase.

If any material or manufacturing fault occurs within three years of the date of purchase of the product, we will either repair or replace the product for you or refund the purchase price (at our discretion). This warranty service requires that you present the defective appliance and the proof of purchase (receipt) within the three-year warranty period, along with a brief written description of the fault and of when it occurred.

If the defect is covered by the warranty, your product will either be repaired or replaced by us. The repair or replacement of a product does not signify the beginning of a new warranty period.

Warranty period and statutory claims for defects

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking. Repairs carried out after expiry of the warranty period shall be subject to a fee.

Scope of the warranty

This appliance has been manufactured in accordance with strict quality guidelines and inspected meticulously prior to delivery.

The warranty covers material faults or production faults. The warranty does not extend to product parts subject to normal wear and tear or to fragile parts which could be considered as consumable parts such as switches, batteries or parts made of glass.

The warranty does not apply if the product has been damaged, improperly used or improperly maintained. The directions in the operating instructions for the product regarding proper use of the product are to be strictly followed. Uses and actions that are discouraged in the operating instructions or which are warned against must be avoided.

This product is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications/repairs which have not been carried out by one of our authorised Service centres.

Warranty claim procedure

To ensure quick processing of your case, please observe the following instructions:

- Please have the till receipt and the item number (IAN) 435653_2304 available as proof of purchase.
- You will find the item number on the type plate on the product, an engraving on the product, on the front page of the operating instructions (below left) or on the sticker on the rear or bottom of the product.
- If functional or other defects occur, please contact the service department listed either by **telephone** or by **e-mail**.
- You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

You can download these instructions along with many other manuals, product videos and installation software at www.lidl-service.com.

This QR code will take you directly to the Lidl service page (www.lidl-service.com) where you can open your operating instructions by entering the item number (IAN) 435653_2304.

Service

GB Service Great Britain

Tel.: 0800 404 7657

E-Mail: kompernass@lidl.co.uk

IE Service Ireland

Tel.: 1800 101010

E-Mail: kompernass@lidl.ie

IAN 435653_2304

Importer

Please note that the following address is not the service address. Please use the service address provided in the operating instructions.

KOMPERNASS HANDELS GMBH

BURGSTRASSE 21

44867 BOCHUM

GERMANY

www.kompernass.com

Ordering spare parts

You can order replacement parts for this product on the Internet at www.kompernass.com.

Scan this QR code with your smartphone or tablet.

You can use this QR code to go directly to our website

www.kompernass.com to view and order the available replacement parts.

i Note

- ▶ If you have problems with your online order, you can contact our service centre by phone or e-mail (see section **Service**).
- ▶ Please always quote the article number (e.g. **435653_2304**) when ordering. You can find it on the title page of these operating instructions.
- ▶ Please note that online ordering of replacement parts is not possible for all countries.

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