SILVER CREST®

TAPPING MASSAGER / KLOPFMASSAGEGERÄT / APPAREIL DE MASSAGE PAR TAPOTEMENT SKM 20 B3

(GB) (IE)

TAPPING MASSAGER

Operating instructions

(FR) (BE)

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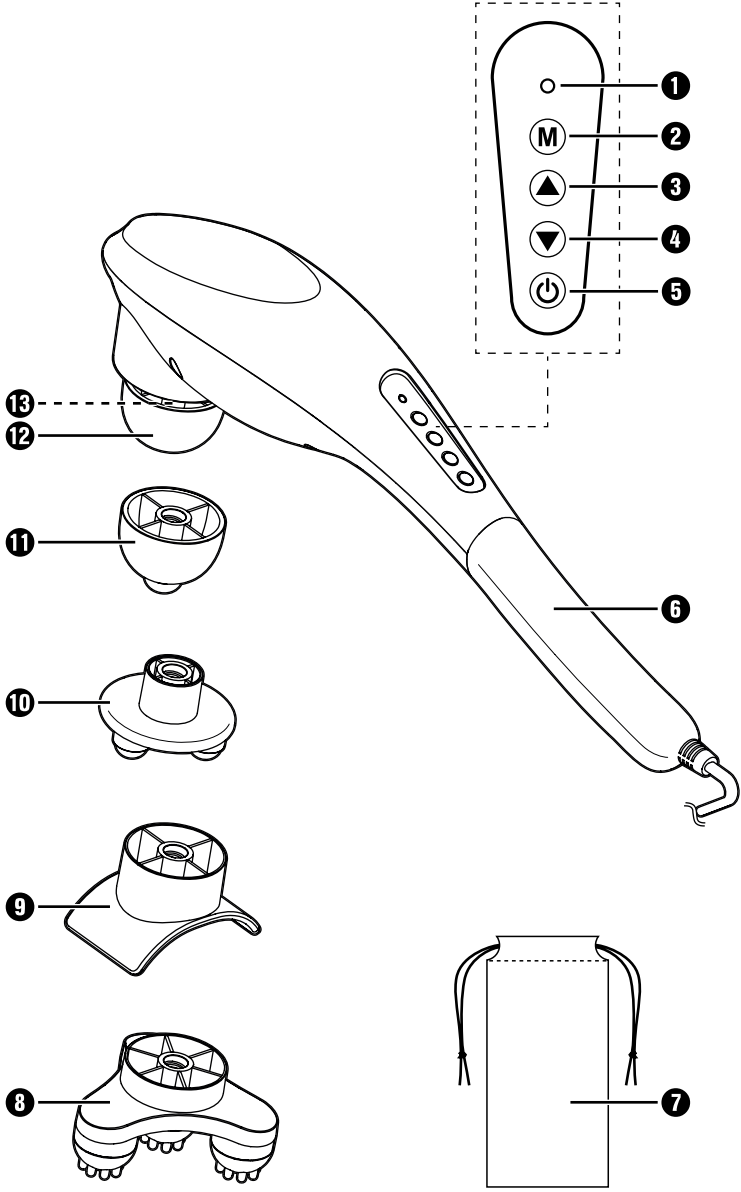
Contents

Introduction	2
Information about these operating instructions	2
Proper use	2
Warnings and symbols used	2
Safety information	3
Description of components	6
Before first use	6
Check package contents	6
Using the appliance	7
Choosing/attaching the massage attachment	7
Switching the appliance on/off	8
Selecting a mode and intensity level	8
Cleaning	9
Storage	9
Troubleshooting	9
Disposal	10
Disposal of the appliance	10
Disposal of the packaging	10
Appendix	10
Technical data	10
Kompernass Handels GmbH warranty	11
Service	12
Importer	12

Introduction

Information about these operating instructions

Congratulations on the purchase of your new appliance.

You have selected a high-quality product. The operating instructions are part of this product. They contain important information about safety, usage and disposal. Before using the product, please familiarise yourself with all operating and safety instructions. Use the product only as described and for the specified range of applications. Please also pass these operating instructions on to any future owner.

Proper use

This appliance is intended for massaging the skin, connective tissue and body muscles. The appliance is only intended for use with the massage attachments supplied. Do not use the appliance directly on the bones, including the spine.

It is not intended for any other purpose nor for use beyond the scope described. The appliance is not intended for use in medical/therapeutic or commercial settings and cannot replace medical treatment. Claims of any type due to damage from improper use are excluded. The risk is borne solely by the user.

Warnings and symbols used

The following warnings and symbols are used in these operating instructions, on the packaging and on the appliance:

	DANGER! A warning with this symbol and the signal word "DANGER" indicates an imminent hazardous situation which will result in death or serious injury if not avoided.
	WARNING! A warning with this symbol and the signal word "WARNING" indicates a potentially hazardous situation which could result in death or serious injury if not avoided.
	CAUTION! A warning with this symbol and the signal word "CAUTION" indicates a potentially hazardous situation which could result in a minor or moderate injury if not avoided.
	ATTENTION! A warning with this symbol and the signal word "ATTENTION" indicates a potential situation which could result in property damage if not avoided.
	Note: A note provides additional information that will assist you in using the appliance.

	AC current/voltage
	Protection class II (double insulation)

Safety information

⚠ DANGER!

- The appliance should only be connected to correctly installed and earthed mains power sockets. Ensure that the specifications for your local power supply correspond to the details on the rating plate of the appliance.
- Make sure that the power cable does not become wet or moist. Route the cable in such a way that it cannot be caught or damaged.
- Keep the appliance and the power cable away from hot surfaces.
- Should the power cable of the appliance be damaged, to avoid potential risks it must be replaced by the manufacturer, his customer service or by a qualified technician.
- Do not perform any repairs to the appliance. All repairs must be performed by our customer service department or by a qualified specialist technician.
- Never open the appliance housing. There are no parts inside the appliance which require user maintenance. In addition, this will invalidate the warranty.
- Switch off the appliance after each use and before cleaning it. For additional protection, we recommend that you install a residual current device with a rated tripping current of no more than 30 mA in the electrical circuit. Ask your electrician for advice.
- Use the appliance in dry interior rooms only; never use it outside.

- Never immerse the appliance in water or other liquids!
- Do not use the appliance while taking a bath or shower. Do not store or keep the appliance in places from where it can fall into the bath tub, sink or other containers filled with water or other liquids.
- Never reach for the appliance if it has fallen into water.
- Do not use the appliance if the power cable is damaged, if the appliance is not working properly, if the attachments show signs of damage, if it has sustained a fall, if it has fallen into water or if it is otherwise damaged.

⚠ WARNING!

- This appliance may be used by children aged 8 years or older and by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, as long as they are being supervised or have been instructed on how to use the appliance safely and are aware of the potential risks.
- Do not allow children to use the appliance as a toy.
- Cleaning and user maintenance may not be carried out by children unless they are supervised.
- People with pacemakers, artificial joints or electronic implants should consult their doctor before using the appliance.
- Do not use the appliance around your heart, face or genital area.
- Do not use the appliance if you suffer from one or more of the following diseases or medical conditions: circulatory disorders, varicose veins, open wounds, bruises, cracked skin, phlebitis, thrombosis, diabetes or slipped disc.

- Never use the appliance on children, animals, during pregnancy, on your head, throat or other parts of the body that are exhibiting swelling, burns, inflammation, wounds or other sensitive areas.
- Do not use the appliance under blankets or pillows. Never use the appliance in the vicinity of petrol or other inflammable materials.

⚠ CAUTION!

- Do not use the appliance while sleeping or after taking substances that may reduce your ability to react quickly (e.g. painkillers or alcohol).
- Never use the appliance directly on bones, e.g. on your spine, shins, the tops of your feet or similar places on your body.
- Use the appliance only on dry and clean parts of your body.
- Stop using the appliance immediately if you find the massage unpleasant or if you experience pain.
- Uninterrupted application of the appliance should not last any longer than 20 minutes.
- Do not insert your fingers or any other objects into the opening for the attachments. Keep long hair away from moving parts.
- Use the appliance only for its intended purpose.

⚠ ATTENTION!

- Use the appliance only indoors.
- Use the appliance only with the supplied accessories.
- Do not use any aggressive cleaning liquids or solvents to clean the appliance. They could damage the surface.

Description of components

(See fold-out page for illustrations)

- ❶ Operating indicator
- ❷ **M** Mode selector button
- ❸ ▲ button (increase intensity level)
- ❹ ▼ button (decrease intensity level)
- ❺ ⏻ On/Off button
- ❻ Handle
- ❼ Storage bag
- ❽ Massage attachment (3 points intensive)
- ❾ Massage attachment (curved)
- ❿ Massage attachment (3 points)
- ⓫ Massage attachment (1 point)
- ⓬ Massage attachment (hemisphere)
- ⓭ Thread

Before first use

Check package contents

 **DANGER!** Risk of suffocation! Do not allow children to play with packaging materials.

- ◆ Remove all parts and the operating instructions from the box.
- ◆ Remove all packaging material.

The following components are included in the delivery:

- 1 × tapping massager
- 5 × massage attachments
- 1 × storage bag (for accessories)
- These operating instructions

 **Note:** Check the package for completeness and signs of visible damage. If the delivery is incomplete or damage has occurred as a result of defective packaging or during transport, contact the service hotline (see the section **Service**).

Using the appliance

The appliance offers a relaxing or even intensive massage for arms, shoulders, neck, back (except spine), buttocks, legs and feet. Depending on the extent of the tension in the part of the body to be treated, a maximum treatment time of about 20 minutes is recommended. The ergonomically shaped handle **6** ensures safe and comfortable handling.

A massage should always feel pleasant. A slight reddening of the treated areas and even slight bruising are a normal side effect, as the massage specifically promotes blood circulation. However, if you experience excessive skin reactions, discontinue use and consult a doctor.

Choosing/attaching the massage attachment

! **ATTENTION!** Damage to the appliance! Switch off the appliance before attaching or removing a massage attachment.

◆ Select the appropriate massage attachment for the desired treatment and body part as follows:

Massage attachment	Treatment/body part
 8	For deep tissue and to relax tense muscles. Large areas all over the body, such as buttocks, legs and back (excluding spine).
 9	For relaxing tense muscles. Especially suitable for arms and legs.
 10	For deep tissue and to relax tense muscles. Particularly suitable for shoulders, neck, back (except spine), legs and feet.
 11	For deep tissue and acupressure points all over the body.
 12	For relaxation purposes all over the body.

◆ Screw the selected massage attachment clockwise onto the thread **13** of the appliance.

◆ Unscrew the attached massage attachment anticlockwise from the thread **13** of the appliance to attach another massage attachment.

Switching the appliance on/off

- ◆ Plug the power plug into a mains socket.
- ◆ Press and hold the  On/Off button **5** for about two seconds to switch on the appliance. A beep sounds, and the operating indicator **1** lights up continuously in white.
- ◆ Press the  On/Off button **5** to switch off the appliance. A beep sounds, and the operating indicator **1** goes out.
- ⓘ **Note:** The appliance switches off automatically after approx. 20 minutes to avoid overheating. Let the appliance cool down for about 10 minutes before switching it on again.
- ◆ Unplug the mains plug from the mains socket when you are no longer using the appliance.

Selecting a mode and intensity level

- ⓘ **Note:** The appliance is in mode 1 and intensity level 1 when switched on.
- ◆ Press the **M** mode selector button **2** to change the mode.

Mode	Function
Mode 1	Continuous pulse
Mode 2	1 pulse per second
Mode 3	1–2 pulses per second
Mode 4	2–3 pulses per second
Mode 5	Combination of modes 1–4

- ◆ Press the  button **3** to increase the intensity level.
- ◆ Press the  button **4** to decrease the intensity level.
- ⓘ **Note:** (1) The appliance offers 5 intensity levels in each of the 5 modes. (2) Each time the mode or intensity level is changed, a beep sounds. (3) Start in mode 1 and with intensity level 1, as preset, to avoid unpleasant feelings or pain. When you find the treatment satisfactory and you are familiar with how to use the appliance, you can switch to a higher intensity level and to another mode.

Cleaning

- ⚠ **DANGER!** Risk of fatal electrical shock! Switch off the appliance before cleaning. Always unplug the appliance from the mains socket before cleaning. Ensure that no moisture can penetrate the appliance during cleaning.
- ⚠ **ATTENTION!** Damage to the appliance! Do not use caustic, abrasive or solvent-based cleaning agents. They can damage the surfaces of the appliance.
- ◆ Wipe off the appliance using a damp cloth. Do not use any solvents or cleansers which can damage plastic. If necessary, add a little mild detergent to the cloth and wipe off afterwards with fresh water. Make sure that no water or moisture can get into the appliance!
- ◆ For hygienic cleaning, wash the massage attachments **8/9/10/11/12** with water and a mild detergent and then rinse them with clean water.
- ◆ Allow everything to dry thoroughly before storing the appliance or putting it back into operation.

Storage

- ◆ Unplug the appliance from the mains socket when not in use.
- ◆ Store the appliance in a clean, dry place out of direct sunlight when not in use.
- ◆ Store the massage attachments **8/9/10/11/12** in the storage bag **7**.

Troubleshooting

Fault	Cause	Remedy
The appliance does not switch on after the  On/Off button 5 is pressed.	The mains plug is not connected to a mains socket.	Insert the plug into a mains socket.
The appliance switches itself off during operation.	The appliance switches off automatically after approx. 20 minutes to avoid overheating.	Let the appliance cool down for about 10 minutes before switching it on again.

Disposal

Applies only to France:

The product, its packaging and the operating instructions are recyclable. They are subject to an extended manufacturer responsibility and will be collected separately.

Disposal of the appliance

The adjacent symbol of a crossed-out dustbin means that this appliance is subject to Directive 2012/19/EU. This directive states that this appliance may not be disposed of in the normal household waste at the end of its useful life, but must be taken to specially set-up collection locations, recycling depots or disposal companies.

The disposal is free of charge for the user. Protect the environment and dispose of this appliance properly.

If your old appliance has stored any personal data, you are responsible for deleting it yourself before returning it.

Your local community or municipal authorities can provide information on how to dispose of the worn-out product.

Disposal of the packaging

The packaging materials have been selected for their environmental friendliness and ease of disposal and are therefore recyclable. Dispose of packaging materials that are no longer needed in accordance with applicable local regulations.

Dispose of the packaging in an environmentally friendly manner. Note the labelling on the packaging and separate the packaging material components for disposal, if necessary. The packaging material is labelled with abbreviations (a) and numbers (b) with the following meanings: 1-7: plastics, 20-22: paper and cardboard, 80-98: composites.

Appendix

Technical data

Operating voltage	220-240 V~, 50-60 Hz
Power consumption	20 W
Protection class	II/□

Kompernass Handels GmbH warranty

Dear Customer,

This appliance has a 3-year warranty valid from the date of purchase. If this product has any faults, you, the buyer, have certain statutory rights. Your statutory rights are not restricted in any way by the warranty described below.

Warranty conditions

The warranty period starts on the date of purchase. Please keep your receipt in a safe place. This will be required as proof of purchase.

If any material or manufacturing fault occurs within three years of the date of purchase of the product, we will either repair or replace the product for you or refund the purchase price (at our discretion). This warranty service requires that you present the defective appliance and the proof of purchase (receipt) within the three-year warranty period, along with a brief written description of the fault and of when it occurred.

If the defect is covered by the warranty, your product will either be repaired or replaced by us. The repair or replacement of a product does not signify the beginning of a new warranty period.

Warranty period and statutory claims for defects

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking. Repairs carried out after expiry of the warranty period shall be subject to a fee.

Scope of the warranty

This appliance has been manufactured in accordance with strict quality guidelines and inspected meticulously prior to delivery.

The warranty covers material faults or production faults. The warranty does not extend to product parts subject to normal wear and tear or to fragile parts which could be considered as consumable parts such as switches, batteries or parts made of glass.

The warranty does not apply if the product has been damaged, improperly used or improperly maintained. The directions in the operating instructions for the product regarding proper use of the product are to be strictly followed. Uses and actions that are discouraged in the operating instructions or which are warned against must be avoided.

This product is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications/repairs which have not been carried out by one of our authorised Service centres.

Warranty claim procedure

To ensure quick processing of your case, please observe the following instructions:

- Please have the till receipt and the item number (IAN) 430498_2301 available as proof of purchase.
- You will find the item number on the type plate on the product, an engraving on the product, on the front page of the operating instructions (below left) or on the sticker on the rear or bottom of the product.
- If functional or other defects occur, please contact the service department listed either by **telephone** or by **e-mail**.
- You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

You can download these instructions along with many other manuals, product videos and installation software at www.lidl-service.com.

This QR code will take you directly to the Lidl service page (www.lidl-service.com) where you can open your operating instructions by entering the item number (IAN) 430498_2301.

Service

GB Service Great Britain

Tel.: 0800 404 7657

E-Mail: kompennass@lidl.co.uk

IE Service Ireland

Tel.: 1800 101010

E-Mail: kompennass@lidl.ie

IAN 430498_2301

Importer

Please note that the following address is not the service address. Please use the service address provided in the operating instructions.

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