AKKU-ZAHNBÜRSTE CBZB 800 A1 **ELECTRIC TOOTHBRUSH CBZB 800 A1** **BROSSE À DENTS ÉLECTRIQUE CBZB 800 A1**

For EU market:
HOYER Handel GmbH
Kühnehöfe 12
22761 Hamburg
GERMANY

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Contents

- 1. Overview 19
- 2. Intended purpose 20
- 3. Safety information 20
- 4. Items supplied 22
- 5. Setting up the battery charger 23
- 6. Charging the battery 23
- 7. Toothbrush heads 23
- 8. Operation 24
 - 8.1 Modes 24
 - 8.2 Brushing teeth 24
- 9. Transport when traveling 25
 - 9.1 Transport lock 25
 - 9.2 Travel case 25
- 10. Cleaning / Hygiene requirements 25
- 11. Storage 25
- 12. Disposal 26
- 13. Trouble-shooting 26
- 14. Ordering toothbrush heads 26
- 15. Technical specifications 27
- 16. Warranty of the HOYER Handel GmbH 28

1. Overview

- 1** Brush head
- 2** Indicator bristles
- 3** Toothbrush head (replaceable)
- 4** Marking ring
- 5** Peg of the toothbrush
- 6** LED:
 - is lit white after switching on
 - is lit red when pressure too high
 - flashes red when the transport lock is switched on
 - flashes 3x red when switching the transport lock on/off
- 7**  On/off switch + transport lock for on and off switch
- 8** Function displays (illuminated):
 -  Clean
 -  Sensitive
 -  Massage
- 9** Charging indicator light
- 10** Toothbrush
- 11** Peg of the battery charger
- 12** Battery charger
- 13** Mains power cable and mains plug
- 14** Travel case
- 15** Pressure cap
- 16** Toothbrush head holder

Thank you for your trust!

Congratulations on the purchase of your new electric toothbrush.

For safe handling of the device and in order to get to know the entire range of features:

- **Thoroughly read these user instructions prior to initial use.**
- **Above all, observe the safety information!**
- **The device should only be used as described in this copy of the user instructions.**
- **Keep this copy of the user instructions for reference.**
- **If you pass the device on to someone else, please include this copy of the user instructions. The user instructions are a part of the product.**

We hope you enjoy using your new electric toothbrush!

All brand and manufacturer names are the property of their respective owners. They are specified solely in order to simplify the explanation of the points illustrated.

2. Intended purpose

The electric toothbrush is intended solely for tooth cleaning and oral hygiene in humans. The device is designed for private, domestic use. This device must not be used for commercial purposes.

3. Safety information

Warnings

If necessary, the following warnings are used in this copy of the user instructions:

DANGER! High risk: failure to observe this warning may result in injury to life and limb.

WARNING! Moderate risk: failure to observe this warning may result in injury or serious material damage.

CAUTION: low risk: failure to observe this warning may result in minor injury or material damage.

NOTE: circumstances and specifics that must be observed when handling the device.

Instructions for safe operation

- ⦿ This device can be used by children and people with restricted physical, sensory or intellectual abilities or people without adequate experience and/or understanding, if they are supervised or instructed on how to use this device and if they are aware of the resulting hazards. Cleaning and user maintenance must not be performed by children unless they are supervised.
- ⦿ Children must not play with the device.
-  ⦿ The battery charger is connected to the power voltage, and must therefore not be used in the shower, in the bath, or over washbasins filled with water.
- ⦿ The mains power cable of this device cannot be replaced. If the cable should become damaged, the device itself must be disposed of.
- ⦿ The battery must be charged only with the original accessory (battery charger) supplied.
- ⦿ This device contains a battery that cannot be replaced.
- ⦿ Misuse may result in injuries.

DANGER for children!

- ⦿ Packing materials are not children's toys. Children should not be allowed to play with the plastic bags: there is a risk of suffocation.

DANGER to and from pets and livestock!

- ⦿ Electrical devices can represent a danger to pets and livestock. In addition, animals can also cause damage to the device. For this reason you should keep animals away from electrical devices at all times.

DANGER! Risk of electric shock due to moisture!

- ⦿ The toothbrush, the battery charger, the mains power cable and the mains plug must never be immersed in water or any other liquids.

- ⦿ Protect the device from moisture, water drops and water splashes: failure to observe this instruction may result in an electrical shock hazard.
- ⦿ Never touch the device with wet hands.
- ⦿ Make sure that the surface under the device is dry.
- ⦿ If the battery charger has fallen into water, pull out the mains plug immediately, and only then remove the battery charger from the water.
If this happens, do not use the battery charger again until it has been checked by a specialist workshop.

DANGER! Risk of electric shock!

- ⦿ The installation of a residual current device with a nominal trip current not exceeding 30 mA in the building wiring offers additional protection.
Contact a skilled electrician in this regard.

- ⊙ Do not start using the device if the device, the battery charger or the mains power cable shows visible damage, or if the device has previously been dropped.
- ⊙ Only connect the mains plug into a properly installed and easily accessible wall socket whose voltage corresponds to the specifications on the rating plate. The wall socket must still be easily accessible after the device is plugged in.
- ⊙ The hand unit of the toothbrush must only be charged using the battery charger supplied.
- ⊙ Ensure that the mains power cable cannot be caught or crushed, and that they cannot be damaged by sharp edges or hot points.
- ⊙ To fully disconnect the battery charger from the power, pull out the mains plug.
- ⊙ When pulling the mains plug out of the wall socket, always pull by the plug and never by the cable.
- ⊙ Disconnect the mains plug from the wall socket:
 - if there is a fault
 - when you are not using the device
 - before you clean the device
 - during thunderstorms
- ⊙ In order to avoid any hazards, do not make modifications to the device. Repairs must be carried out by a specialist workshop or the Service Centre.
- ⊙ When using an extension cable, this must correspond to the applicable safety regulations.

DANGER from batteries!

- ⊙ Protect the battery against mechanical damages. **Fire hazard!**
- ⊙ Do not expose the device to direct sunshine or heat. The ambient temperature should not be lower than 10 °C or higher than 40 °C.
- ⊙ The battery must be charged only with the original accessory (battery charger) supplied.
- ⊙ Should the electrolyte solution leak out of the battery, avoid contact with eyes, mucous membranes and skin. Rinse immediately with plenty of clean water and consult a doctor. The electrolyte solution can cause irritation.
- ⊙ The device is provided with a lithium-ion battery.
 - The battery cannot be removed!
 - This device must not be opened!
 - The device with the installed battery must be disposed of properly.

CAUTION

- ⊙ Users could possibly experience mild bleeding of the gums during the first days of electric toothbrush usage. Should the bleeding continue over several days, a dentist should be consulted.
- ⊙ Exert only light pressure when brushing to prevent damage to the gums and to increase the service life of the bristles.

4. Items supplied

- 1 toothbrush **10**
- 2 toothbrush heads **3**
- 1 battery charger **12**
- 1 travel case **14**
- 1 toothbrush head holder **16**
- 1 operating instructions

- Remove all packing materials.
- Check that the device is undamaged.

5. Setting up the battery charger

DANGER! Risk of electric shock due to moisture!

- ⊙ The battery charger **12** is connected to the power voltage, and must therefore not be used in the shower, in the bath, or over washbasins filled with water.
- Never set up the battery charger **12** anywhere it could fall into water.

6. Charging the battery

NOTES:

- Fully charge the battery before using it for the first time.
 - It takes approx. 10 hours to fully charge the battery.
 - Place the toothbrush **10** on the peg of the battery charger **11** after use. This ensures that the toothbrush **10** always works at full power. The battery is not damaged by regular charging.
1. Plug the mains plug **13** of the battery charger **12** into a suitable wall socket. The wall socket must still be easily accessible after connection.
 2. Place the switched-off toothbrush **10** onto the peg of the battery charger **11**. All segments of the charging indicator light **9** briefly light up one after the other during the charging process. Permanently lit segments of the charging indicator light show the progress of the charging process. When the battery is fully charged, first all segments of the charging indicator light **9** are lit, and then they go out.

Battery signals: charging indicator light 9

Toothbrush **10** on the battery charger **12**:

- All segments of the charging indicator light **9** are lit and then go out: the battery is fully charged.
- All segments of the charging indicator light **9** briefly light up one after the other: the battery is being charged.

Toothbrush **10** switched on:

- The last segment of the charging indicator light **9** flashes: the battery is weak. Charge the battery.

7. Toothbrush heads

The items supplied include 2 toothbrush heads **3**.

The various marking rings **4** allow you to assign the toothbrush heads **3** to different family members.

- Take a toothbrush head **3** and slide it onto the peg of the toothbrush **5**.
- In order to remove, carefully pull the toothbrush head **3** off upwards.
- The indicator bristles **2** on the brush head **1** slowly lose their colour during use. The decolouration shows when the toothbrush heads **3** should be replaced.

NOTE: we recommend that you replace the toothbrush heads **3** every 3 months or after infections (see "Ordering toothbrush heads" on page 26).

8. Operation

CAUTION:

- ⊙ Users could possibly experience mild bleeding of the gums during the first days of electric toothbrush usage. Should the bleeding continue over several days, a dentist should be consulted.
 - ⊙ If you have had dental surgery in the last two months, consult a dentist before using the toothbrush.
 - ⊙ Exert only light pressure when brushing to prevent damage to the gums and to increase the service life of the bristles.
-

8.1 Modes

The toothbrush **10** has 3 different modes. Depending on the desired mode, press the on/off switch **7** repeatedly when switched on, until the corresponding symbol is lit:

Clean mode

- for cleaning teeth with high pulsation

Sensitive mode

- gentle brushing for cleaning teeth with sensitive gums
-

NOTE: if you are not used to using an electric toothbrush and you have comparatively sensitive gums, you can use this mode for the first 1 to 2 weeks, and then switch to Clean mode later.

Massage mode

- pulsating brushing to massage your gums

8.2 Brushing teeth

1. Remove the toothbrush **10** from the battery charger **12**.
 2. Switch on the toothbrush using the on/off switch **7**.
 3. Press on/off switch **7** repeatedly until the function display **8** for the desired mode is lit.
-

NOTE: switching over to another mode must be done quickly after the last button press of the on/off switch **7** (**within approx. 1 second**). Pressing the on/off switch **7** approx. 2 seconds after the last button press will switch the toothbrush off.

4. Switch off the toothbrush using the on/off switch **7**.
5. Moisten the bristles, apply some toothpaste, and only switch on the toothbrush **10** again when the brush head **1** is already in your mouth.
6. Slowly guide the brush head **1** from tooth to tooth without exerting too much pressure. Carefully clean all tooth surfaces, the inner surfaces of the teeth and the tooth edges.

Figure A: the recommended brushing time is at least 2 minutes, divided up into 30-second periods for each area of the mouth (bottom right, bottom left, top right, top left). The toothbrush **10** indicates these intervals by two brief interruptions.

NOTE: if the brush head **1** is pressed against the teeth too hard, the LED **6** lights up red. You should then reduce the pressure.

7. The toothbrush signals the end of a cycle of 2 minutes with four short interruptions. Remove the brush head **1** from the mouth only when the toothbrush **10** has been switched off.

NOTES:

- We recommend that you replace the toothbrush heads **3** every 3 months or after infections (see "Ordering toothbrush heads" on page 26).
 - The indicator bristles **2** on the brush head **1** slowly lose their colour during use. The decolouration shows when the toothbrush heads **3** should be replaced.
 - After use, the toothbrush heads **3** can be placed on the toothbrush head holder **16**.
-

9. Transport when traveling

9.1 Transport lock

To secure the toothbrush **10** against accidentally being switched on, activate the transport lock.

- Switching on the transport lock:
Hold down the on/off switch  **7** until the LED **6** flashes red three times.
- If you try to switch on the toothbrush **10** with the transport lock activated, the LED **6** will flash red.
- Switching off the transport lock:
Hold down the on/off switch  **7** until the LED **6** flashes red three times.

9.2 Travel case

The travel case **14** has space for the toothbrush **10** and 2 toothbrush heads **3**.

- Open the travel case **14** by pressing on the pressure cap **15**.
- Close the travel case **14** by carefully pressing the lid down until it audibly clicks into place.

10. Cleaning / Hygiene requirements

DANGER! Risk of electric shock!

- ⊙ Before cleaning the battery charger **12**, pull the mains plug **13** out from the mains socket.
 - ⊙ Do not immerse the toothbrush **10** and battery charger **12** in water! The battery charger **12** must be completely dry before you reconnect it to the wall socket.
-

- To ensure the hygienic functionality of the toothbrush heads **3** they must be rinsed thoroughly after every use and allowed to dry.
- Remove the toothbrush head **3** once per week and clean it in hot water.
- Rinse the toothbrush **10** as well and wipe it dry.
- Wipe the battery charger **12**, the travel case **14** and the toothbrush head holder **16** with a damp cloth if necessary and rub it dry.

11. Storage


- Disconnect the device from the mains and store it in a dry location.

12. Disposal

The battery integrated into this device must not be disposed of in household waste. The device must be disposed of properly with the battery installed.

This product is subject to the provisions of European Directive 2012/19/EC. The symbol showing a wheeled bin crossed through indicates that the product requires separate refuse collection in the European Union. This applies to the product and all accessories marked with this symbol. Products marked as such may not be discarded with normal domestic waste, but must be taken to a collection point for recycling electrical and electronic devices.

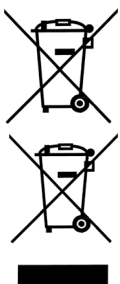This recycling symbol marks e.g. an object or parts of the material as valuable for recycling. Recycling helps to reduce the consumption of raw materials and protect the environment.

Packaging

When disposing of the packaging, make sure you comply with the environmental regulations applicable in your country.

13. Trouble-shooting

If your device fails to function as required, please try this checklist first. Perhaps there is only a minor problem, and you can solve it yourself.

DANGER! Risk of electric shock!

- ⊙ Do not attempt to repair the device yourself under any circumstances.

Fault	Possible causes / Action
The device cannot be switched on.	- Is the battery empty? - Is the transport lock activated? Press the on/off switch 7 repeatedly until the toothbrush 10 is switched on.
Charging indicator light 9 is not lit.	The battery charger 12 is not connected to the mains power.

14. Ordering toothbrush heads

You can order additional toothbrush heads **3** on the internet.

Ordering online

<https://shop.hoyerhandel.com/gb/>

15. Technical specifications

Model:	CBZB 800 A1
Toothbrush	
Protection class:	IPX7 (protected against the effects of temporary immersion)
Battery:	1x 3.7 V === Li ion 800 mAh 2.96 Wh
Charging time:	10 hours
Operating time:	approx. 90 min
Battery charger	
Mains voltage:	100 - 240 V ~ 50/60 Hz
Power rating:	1.4 W Standby: 0.5 W
Protection class:	II 
Protection class:	IPX7 (protected against the effects of temporary immersion)
General	
LEDs:	RG 0
Ambient conditions:	10 °C to 40 °C humidity 40 to 85% relative

NOTE: no action by the user is required to switch the product between 50 and 60 Hz. The product adjusts itself for both 50 and 60 Hz.

Information on batteries/rechargeable batteries according to EU 2023/1542

Name / trade-mark of the manufacturer	Shenzhen Bofuneng Battery Co.,Ltd. A Building, Chun yang Industrial Park, Zhugushi Road, Wulian Community, Long Cheng Street, Longgang District, Shenzhen, P.R.China
Battery type	Li-ion 14500 Batch number: 2605
Brand name	BFN
The declaration of conformity is issued under the sole responsibility of:	Shenzhen Bofuneng Battery Co.,Ltd.
The mentioned battery was produced in accordance with EU Regulation EU2023/1542.	

Symbols used

	Double insulation
	G eprüfte S icherheit (certified safety): devices must comply with the generally acknowledged rules of technology and the German Product Safety Act (Produktsicherheitsgesetz - ProdSG). (Not for UK)
	With the CE marking, HOYER Handel GmbH declares the conformity with EU guidelines.
	This symbol reminds you to dispose of the packaging in an environmentally friendly manner.
	Recyclable materials are marked with the recycling symbol (3 arrows). The material can be specified by the recycling number in the centre (here: 21) and/or an abbreviation (here: PAP).
	Direct current
	AC voltage
	It is a recyclable product which is subject to an extended producer responsibility and for which the rules of waste separation apply.
	With the UKCA marking, HOYER Handel GmbH declares that the product is in conformity with the UKCA requirements.
	The symbol identifies the product manufacturer.

Subject to technical modifications.

16. Warranty of the HOYER Handel GmbH

Dear Customer,
your device is provided with a 3 year warranty starting with the purchase date. In the event of product defects, you are entitled to statutory rights against the vendor. These statutory rights are not restricted by our warranty presented in the following.

Warranty conditions

The warranty period starts with the purchase date. Please keep the original purchase receipt in a safe place. This document is required to verify the purchase.

If within three years from the purchase date of this product a material or factory defect occurs, the product will be repaired or replaced by us – at our discretion – free of charge to you, or the purchase price will be refunded. This warranty implies that within the period of three years the defective device and the purchase receipt are presented, including a brief written description of the defect and the time it occurred.

If the defect is covered by our warranty, the repaired or a new product will be returned to you. No new warranty period starts with a repair or replacement of the product.

Warranty period and statutory claims for defects

The warranty period is not extended when the warranty has been claimed. This also applies to replaced and repaired parts. Any damages and defects already existing at the time of purchase must be reported immediately upon unpacking. Repairs arising after the expiration of the warranty period are subject to a charge.

Warranty coverage

The device was produced carefully according to strict quality guidelines and tested diligently prior to delivery.

The warranty applies to material or factory defects.

Excluded from the warranty are wear parts subject to normal wear and damages to fragile parts, e.g. switches, lamps or other parts manufactured from glass.

This warranty expires if the product is damaged, not used as intended or not serviced. For the proper operation of the product, all instructions listed in the operating instructions must be observed carefully. Any form of use and handling that is advised against in the operating instructions or warned against must always be avoided.

The product is only intended for private and not for commercial use. In the case of incorrect and improper treatment, use of force and interventions not performed by our authorised Service Centre, the warranty shall cease.

Handling in case of a warranty claim

In order to ensure prompt processing of your matter, please observe the following notes:

- Please keep the article number **IAN: 516700_2510** and the purchase receipt as a purchase verification for all inquiries.
- The article number can be found on the rating plate, an engraving, the title page of your instructions (in the bottom left), or as a label on the rear or underside of the device.
- If faulty operation or other defects occur, first contact the Service Centre listed in the following by **telephone** or **email**.
- Then, you are able to send a product reported as defective free of charge to the service address specified to you, including the purchase receipt and the information on the defect and when it occurred.

At www.lidl-service.com you can download this and many other manuals, product videos and installation software.

This QR code takes you directly to the Lidl Service website (www.lidl-service.com), where you can enter the article number (IAN) **516700_2510** to open your user instructions.

Service Centre

GB Service Great Britain
Tel.: 0800 404 7657
E-Mail: hoyer@lidl.co.uk

IE Service Ireland
Tel.: 1800 101010
E-Mail: hoyer@lidl.ie

MT Service Malta
Tel.: 80062230
E-Mail: hoyer@lidl.com.mt

IAN: 516700_2510

Supplier

Please note that the following address is **no service address**. First contact the aforementioned Service Centre.

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