WATER JET FLOSSER / MUNDDUSCHE HYDROPULSEUR CBMDA 3.7 A1

GB IE NI MT

WATER JET FLOSSER

Operating instructions

FR BE CH

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Mode d'emploi

CZ

ÚSTNÍ SPRCHA

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SK

ÚSTNA SPRCHA

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DK

MUNDSKYLLER

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HU

SZÁJZUHANY

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DE AT BE CH

MUNDDUSCHE

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NL BE

MONDDOUCHE

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PL

IRYGATOR

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ES

IRRIGADOR BUCAL

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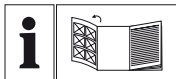
IT CH MT

IDROPULSORE

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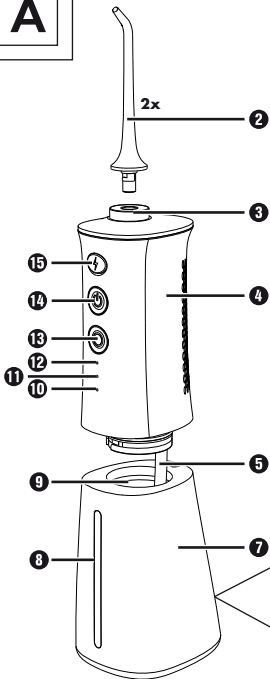
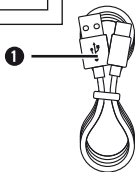
IAN 508148_2507

DE FR BE NL
CZ PL SK ES

English	Operating instructions	Page	1
Deutsch	Bedienungsanleitung	Seite	47
Français	Mode d'emploi	Page	101
Nederlands	Gebruiksaanwijzing	Pagina	159
Česky	Návod k obsluze	Strana	207
Polski	Instrukcja obsługi	Strona	253
Slovenčina	Návod na obsluhu	Strana	303
Español	Instrucciones de uso	Página	349
Dansk	Betjeningsvejledning	Side	399
Italiano	Istruzioni per l'uso	Pagina	445
Magyar	Használati utasítás	Oldal	493

A**B**

Contents

1. Introduction	3
2. Intended use	3
3. Package contents	4
4. Warnings and symbols used	5
5. Safety information	10
5.1. Safety guidelines for rechargeable batteries	16
6. Appliance description	19
7. Preparations	20
7.1. Charging the appliance	21
7.2. Filling the water tank	23
7.3. Fitting/removing attachments	25
8. Power levels	27
9. Using the appliance	28
9.1. After use	32
10. Cleaning and care	33
11. Storage	35

12. Ordering replacement parts	36
13. Technical specifications	37
14. Disposal	38
14.1. Disposal of the appliance	38
14.2. Disposal of batteries	39
14.3. Disposal of the packaging	40
15. Kompernass Handels GmbH warranty	40
15.1. Service	45
15.2. Importer	46

1. Introduction

Congratulations on the purchase of your new appliance. You have selected a high-quality product.

The instructions for use are part of this product. They contain important information about safety, usage and disposal. Before using the product, please familiarise yourself with all instructions regarding safety and use. Use the product only as described and for the specified range of applications. Please also pass on all documents to any future owner.

2. Intended use

This appliance is intended exclusively for personal dental care at home. It can also be used to clean bridges, crowns or fixed braces. It is intended exclusively for use in domestic households. This appliance is not intended for use in industrial, commercial or medical environments.

It is not suitable for animal care. Any other use or use beyond that specified is considered improper use.

3. Package contents

The appliance is supplied with the following components as standard:

- Handpiece with water tank
- USB-C charging cable
- Instructions for use

Note

- ▶ Check the package for completeness and signs of visible damage.
- ▶ If the delivery is incomplete or damage has occurred as a result of defective packaging or during transport, contact the Customer Service hotline (see section **15.1. Service**).

4. Warnings and symbols used

The following warnings and symbols are used in these instructions for use, on the packaging and on the appliance (where applicable):

	<p>DANGER! A warning with this symbol and the signal word "DANGER" indicates an imminently hazardous situation that will result in death or serious injury if not avoided.</p>
	<p>WARNING! A warning with this symbol and the signal word "WARNING" indicates a potentially hazardous situation that could result in death or serious injury if not avoided.</p>

	<p>CAUTION! A warning with this symbol and the signal word "CAUTION" indicates a potentially hazardous situation that could result in a minor or moderate injury if not avoided.</p>
	<p>ATTENTION! A warning with this symbol and the signal word "ATTENTION" indicates a potential situation that could result in property damage if not avoided.</p>
	<p>Note: A note identifies additional information that facilitates the use of the appliance.</p>
	<p>DC current/voltage</p>
	<p>Protection class III: Protection through extra low voltage</p>

	Removable power supply
	Read the instructions.
	USB [®] is a registered trademark of USB Implementers Forum, Inc.
	Do not dispose of electrical appliances in the household waste!
	Dispose of the packaging in an environmentally-responsible manner.

The packaging is made of recyclable materials. Observe the labelling on the packaging material when sorting for disposal: The material is labelled with abbreviations (a) and numbers (b) with the following meanings: 1-7: plastics, 20-22: paper and cardboard, 80-98: composites.

The packaging contains paper and/or cardboard components.

The packaging contains plastic and/or metal components.

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FR

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FR

FR: The product, its
packaging and the
operating instructions

are recyclable. They are subject to an
extended manufacturer responsibility and
will be collected separately.

5. Safety information

DANGER!

- Make sure that the charging cable does not become wet or damp. Route the cable in such a way that it cannot be caught or damaged.
- Keep the appliance and the charging cable away from hot surfaces.
- Do not perform any repairs to the appliance. All repairs must be performed by our customer service department or by a qualified specialist technician.

- Never open the handpiece housing. There are no parts inside the appliance that require user maintenance. In addition, this will invalidate the warranty.
- Switch off the appliance after each use and before cleaning it.
- Have all repairs carried out by a specialist workshop. Under no circumstances should you open the appliance yourself. Repairs that are not carried out by a specialist workshop could lead to physical injuries.

WARNING!

- This appliance may be used by children over the age of 8 and by persons with reduced physical, sensory or mental capabilities or lack of experience and/or knowledge if they are supervised.
- Do not allow children to play with the appliance.
- Always store the appliance indoors to avoid accidents. After use, store the appliance in a dry place away from direct sunlight.

- Switch off the appliance before you put it down or change accessories.

ATTENTION!

- Use only accessories recommended by the manufacturer.
-  Do not charge or use the appliance outdoors.
- Use the appliance only with the supplied accessories. Use only the supplied charging cable to charge the appliance.

- If you are using a USB power supply, the mains socket must always be easily accessible so that the USB power supply can be quickly disconnected from the socket in a hazardous situation.
You should also observe the operating instructions for the USB power supply.
- Do not use any aggressive cleaning agents or solvents. They could damage the surface.
- The appliance may only be used with cold or lukewarm

tap water up to a maximum of 40°C. The water tank should never be filled with additives such as mouth-wash or mouth rinse. These can damage the appliance.

- This appliance contains a battery that cannot be replaced by the user.
- The product is equipped or supplied with a lithium-ion battery. Do not throw the rechargeable battery into a fire and do not subject it to high temperatures. There is a risk of an explosion!

5.1. Safety guidelines for rechargeable batteries

The product is equipped or supplied with a lithium-ion battery. The battery in this appliance cannot be replaced.

RISK OF EXPLOSION!

- Do not open or short-circuit the battery.
This could cause the battery to overheat and explode.
There is a risk of fire!
- Never throw the battery into a fire or water.
- Do not subject the battery to mechanical stress.

- Never leave the appliance unattended during charging.

WARNING!

RISK OF INJURY DUE TO BATTERY LEAKAGE!

- High temperatures can affect the battery. Avoid direct sunlight and proximity to heat sources such as radiators.
- Under extreme conditions, the battery cells may leak. Avoid contact with skin, eyes and mucous membranes. In the event of accidental contact, rinse the affected areas with plenty of clean water and contact a doctor.

- Leaking or damaged batteries can cause chemical burns to the skin. Always wear protective gloves if you have to touch the batteries.

ⓘ **ATTENTION!**

- Never store the appliance when it is discharged. Prolonged storage in a discharged state can cause permanent damage to the battery. If the appliance is stored for a longer period of time, the charge level of the battery must be checked regularly.

The optimum charge level is between 50% and 80%.

6. Appliance description

- 1 USB-C charging cable
- 2 Attachment
- 3 Socket
- 4 Handpiece
- 5 Suction tube
- 6 Small water tank opening
- 7 Water tank
- 8 Water level indicator
- 9 Large water tank opening
- 10 "**Jet**" indicator light
- 11 "**Normal**" indicator light
- 12 "**Soft**" indicator light
- 13 "Power level" button
- 14 On/Off button 
- 15 USB-C port (with protective cap)

7. Preparations

Always clean the appliance thoroughly before first use:

- 1) Wipe the handpiece ④ and the exterior surfaces of the water tank ⑦ with a damp cloth to remove any packaging or production residues.
- 2) Rinse the two attachments ② under flowing water and wipe them both dry.
- 3) Rinse the water tank ⑦ with clean water and then let it dry while the small and large water tank openings ⑥/⑨ are open.
- 4) Thoroughly charge the appliance as described in the section **7.1. Charging the appliance** before using it for the first time.

7.1. Charging the appliance

i Note

- ▶ The USB-C charging cable **1** and USB-C charging port **15** are used only to charge the integrated battery. They are not suitable for data transmission.
- ▶ Charge the appliance with a USB power adapter with protection class II that is approved for use with household appliances and has an output voltage of 5 V $\overline{=}$ and a maximum output current of 1 A.
- ▶ Use only the supplied USB-C charging cable **1** to charge the integrated battery. Disconnect the USB-C charging cable **1** from the appliance when charging is complete.
- ▶ The operating time with a fully charged battery is about 60 minutes, depending on the selected power level.

- ▶ The appliance cannot be used during charging.
- 1) Make sure that the appliance is switched off.
- 2) Open the protective cap over the USB-C charging port **15**.
- 3) Insert the USB-C plug of the USB-C charging cable **1** into the USB-C charging port **15**.
- 4) Connect the USB-C charging cable **1** to a suitable mains adapter and plug the mains adapter into a mains socket.
- 5) While the battery is charging, the indicator lights **10/11/12** light up in sequence. When the battery is fully charged, all indicator lights **10/11/12** light up at the same time.
- 6) Disconnect the appliance from the power supply and unplug the USB-C charging cable **1** from the appliance.
- 7) Close the protective cap.

The battery can be recharged at any time and at any battery level.

7.2. Filling the water tank

ATTENTION!

- ▶ Fill the water tank **7** with pure tap water. The water tank **7** should never be filled with additives such as mouth-wash or mouth rinse.
 - ▶ Use only cold or lukewarm water with a maximum temperature of 40°C. Do not fill the water tank **7** with hot water.
- 1) To easily fill the water tank **7**, remove it from the handpiece **4** by pulling them apart with a light turning motion.
 - 2) Make sure that the small water tank opening **6** is closed.
 - 3) Fill the water tank **7** under flowing water or fill it with a measuring cup with a maximum of 230 ml of water.

- 4) Attach the water tank **7** to the handpiece **4** again (see Fig. 1):
Insert the suction tube **5** into the water tank **7**, place the handpiece **4** on the water tank **7** with a slight offset and lock the two parts together with a light turning motion.

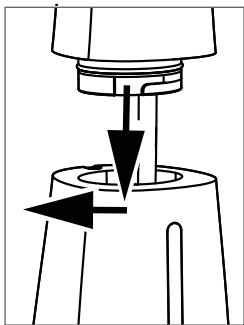
Fig. 1

7.3. Fitting/removing attachments

i Note

- ▶ Make certain that the handpiece ④ is switched off before fitting or removing the attachments ②.
- ▶ If used daily, replace the attachments ② after no more than six months.
- ▶ Use only the attachments ② recommended by the manufacturer. Information on ordering replacement attachments ② can be found in the section **12. Ordering replacement parts.**

- ◆ To fit an attachment **2** to the handpiece **4**, insert it into the socket **3** with a light turning motion (see Fig. 2):

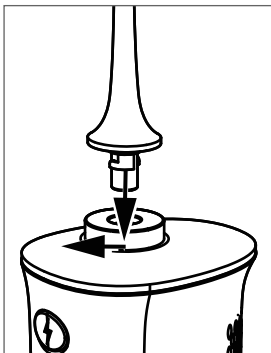
Fig. 2

- ◆ To remove an attachment **2** from the handpiece **4**, turn it slightly and pull it off upwards.

8. Power levels

This oral irrigator has three different power levels:

- **Soft:** Soft water jet, **Soft** indicator light **12** lights up.
For gently massaging your gums.
- **Normal:** Water jet mixed with air, **Normal** indicator light **11** lights up.
For cleaning left over food from your interdental spaces, for massaging your gums, for rinsing between teeth and gums.
- **Jet:** Pulsing cleaning, **Jet** indicator light **10** lights up.
For removal of stubborn food residue.

9. Using the appliance

CAUTION!

- ▶ The unfamiliar stimulation may cause some gum bleeding at first, which is not unusual. If this happens, use the lower power level **soft**. If the bleeding gums last longer than two weeks, you should contact your dentist.
- ▶ Consult your dentist before using the appliance if you suffer from severe periodontitis, have injuries or abscesses in the mucous membrane of the oral cavity or underwent dental surgery within the last two months.

Note

- ▶ Use of the appliance is not a replacement for daily tooth cleaning with a toothbrush and toothpaste.
Use the appliance only as a supplement to daily tooth cleaning with toothbrush and toothpaste.

The appliance is correctly assembled and the water tank **7** is filled:

- 1) Press the "Power level" button **13** repeatedly until the indicator light of the desired power level is lit.

i **Note**

- ▶ The appliance switches off automatically after about 2 seconds if no other button is pressed.
- ▶ The power level can also be changed during operation by pressing the "Power level" button **13** again.

- 2) Place the attachment **2** in your mouth to prevent spray. Bend over the wash-basin and hold the handpiece **4** as upright as possible (see Fig. 3). Direct the nozzle of the attachment **2** towards your teeth/interdental spaces.

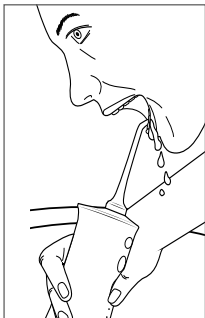
Fig. 3

- 3) Press the On/Off button **14** to switch on the appliance. Close your mouth a little to prevent water from spraying around.
- 4) To clean the interdental spaces, slowly guide the nozzle of the attachment **2** from interdental space to interdental space on the inner and outer sides of your teeth. To massage your gums, point the nozzle at the gums. Move the water jet slowly over your gums.

The length of the process will depend on how much cleaning the interdental spaces require. Continue cleaning until your teeth feel clean.

When massaging your gums, use the oral irrigator for as long as you find the massage comfortable.

i Note

- ▶ If the water in the water tank **7** is used up in the meantime, stop the appliance. Remove it from your mouth and refill the water tank **7** through the small water tank opening **6**.
 - ▶ The device switches off automatically if no button is pressed for 2 minutes.
- 5) If you want to exit the application before the automatic shutdown after 2 minutes, press the On/Off button **14** to switch off the appliance.
- 6) Only then should you take the attachment **2** out of your mouth and rinse it under clean, running water.

9.1. After use

i Note

- ▶ Always empty the water tank **7** and all tubes after every use to prevent the build-up of bacteria.

When the handpiece **4** is switched off and you have finished the treatment:

- 1) Open the small water tank opening **6** with the attachment **2** still fitted and pour out the rest of the water.
- 2) Close the small water tank opening **6**.
- 3) Point the nozzle of the attachment **2** downwards into a washbasin. Switch the appliance on for a couple of seconds to flush any remaining water from the water tank **7**.

10. Cleaning and care

DANGER!

- ▶  NEVER submerge the hand-piece **4** in water or other liquids!
- ▶ Do not clean the appliance while it is charging. Switch off the appliance before cleaning it.

ATTENTION!

- ▶ Do not use abrasive, aggressive or chemical cleaning agents. They could damage the outer surfaces of the appliance.
- ◆ Make sure that the appliance is switched off.
- ◆ Make sure that the appliance is not connected to a mains adapter. Unplug the USB-C charging cable **1** from the appliance if it is plugged in.

- ◆ Remove the attachment ② from the handpiece ④. After use, clean the attachments ② thoroughly under running water and empty the appliance as described in section **9.1. After use**.
- ◆ To clean the handpiece ④, remove the water tank ⑦ and wipe off the handpiece ④ with a damp cloth. For more stubborn soiling, add a little detergent to the cloth and rinse it off afterwards with a cloth dampened with clean water. Dry the handpiece ④ afterwards.
- ◆ Rinse the water tank ⑦ thoroughly with warm water, including from the outside. Empty the water tank ⑦ completely into a sink after cleaning. If you are not using the appliance for over a week, leave the water tank ⑦ open to ensure it can dry fully.

11. Storage

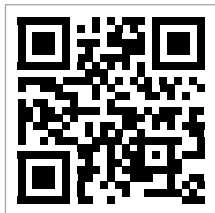
- ◆ Remove any water from the appliance and allow the appliance to dry completely (see section **10. Cleaning and care**).
- ◆ Store the appliance in a location that is clean, dry, dust-free and out of direct sunlight. Charge the battery at regular intervals to ensure a long battery life.

ⓘ **ATTENTION!**

- ▶ Never store the appliance when it is discharged. Prolonged storage in a discharged state can cause permanent damage to the battery. If the appliance is stored for a longer period of time, the charge level of the battery must be checked regularly. The optimum charge level is between 50% and 80%.

12. Ordering replacement parts

You can order replacement parts for this product on the Internet at www.kompernass.com.

Scan the QR code with your smart-phone/tablet. You can use this QR code to go directly to our website to view and order the available spare parts.

i Note

- ▶ If you have problems with your online order, you can contact our service centre by phone or e-mail.
- ▶ Always quote the article number (IAN) 508148_2507 in your order.

- ▶ Please note that online ordering of replacement parts is not possible for all countries.

13. Technical specifications

Input voltage	5.0 V 
Input current	1.0 A
Integrated battery (Li-ion)	1500 mAh / 3.7 V 
Protection class	IPX7 (protection against temporary submersion)

14. Disposal

14.1. Disposal of the appliance

The symbol of the crossed-out wheelie bin means that this appliance may not be disposed of in regular household waste at the end of its service life. The appliance must be deposited at an established collection point, recycling centre or disposal company. Please erase all personal data before returning an appliance.

Before returning an appliance, please remove any batteries that are not permanently enclosed within the old appliance as well as any light sources that can be removed without destroying the appliance. Then dispose of these items separately.

If the rechargeable batteries are permanently installed, you must indicate during disposal that the appliance contains a battery.

The built-in rechargeable battery in this appliance cannot be removed for disposal.

14.2. Disposal of batteries

Batteries/rechargeable batteries must be treated as hazardous waste and must therefore be disposed of in an environmentally sound manner by appropriate bodies (dealers, specialist dealers, public municipal bodies, commercial disposal companies). Batteries/rechargeable batteries may contain toxic heavy metals.

The heavy metals contained are identified by letters below the symbol: Cd = cadmium, Hg = mercury, Pb = lead.

For this reason, do not dispose of batteries/ rechargeable batteries in domestic waste. Take them to a specialist collection point. Only return batteries that are fully discharged.

14.3. Disposal of the packaging

The packaging is made from environmentally friendly material which can be disposed of at your local recycling plant.

Dispose of the packaging in an environmentally friendly manner.

15. Kompernass Handels GmbH warranty

Dear Customer,

This appliance has a 3-year warranty valid from the date of purchase. If this product has any faults, you, the buyer, have certain statutory rights. Your statutory rights are not restricted in any way by the warranty described below.

Warranty conditions

The warranty period starts on the date of purchase. Please keep your receipt in a safe place. This will be required as proof of purchase.

If any material or manufacturing fault occurs within three years of the date of purchase of the product, we will either repair or replace the product for you or refund the purchase price (at our discretion). This warranty service requires that you present the defective appliance and the proof of purchase (receipt) within the three-year warranty period, along with a brief written description of the fault and of when it occurred.

If the defect is covered by the warranty, your product will either be repaired or replaced by us. The repair or replacement of a product does not signify the beginning of a new warranty period.

Warranty period and statutory claims for defects

The warranty period is not prolonged by repairs effected under the warranty. This also applies to replaced and repaired components. Any damage and defects present at the time of purchase must be reported immediately after unpacking. Repairs carried out after expiry of the warranty period shall be subject to a fee.

Scope of the warranty

This appliance has been manufactured in accordance with strict quality guidelines and inspected meticulously prior to delivery.

The warranty covers material faults or production faults. The warranty does not extend to product parts subject to normal wear and tear or to fragile parts which could be considered as consumable parts such as switches or parts made of glass.

The warranty does not apply if the product has been damaged, improperly used or improperly maintained. The directions in the operating instructions for the product regarding proper use of the product are to be strictly followed. Uses and actions that are discouraged in the operating instructions or which are warned against must be avoided. This product is intended solely for private use and not for commercial purposes. The warranty shall be deemed void in cases of misuse or improper handling, use of force and modifications / repairs which have not been carried out by one of our authorised Service centres.

Warranty claim procedure

To ensure quick processing of your case, please observe the following instructions:

- Please have the till receipt and the item number (IAN) 508148_2507 available as proof of purchase.

- You will find the item number on the type plate on the product, an engraving on the product, on the front page of the operating instructions (below left) or on the sticker on the rear or bottom of the product.
- If functional or other defects occur, please contact the service department listed either by **telephone** or by **e-mail**.
- You can return a defective product to us free of charge to the service address that will be provided to you. Ensure that you enclose the proof of purchase (till receipt) and information about what the defect is and when it occurred.

You can download these instructions along with many other manuals, product videos and installation software at www.lidl-service.com.

This QR code will take you directly to the Lidl service page (www.lidl-service.com) where you can open your operating instructions by entering the item number (IAN) 508148_2507.

15.1. Service

- GB** **Service Great Britain**
Tel.: 0800 404 7657
E-Mail: kompennass@lidl.co.uk
- IE** **Service Ireland**
Tel.: 1800 101010
E-Mail: kompennass@lidl.ie
- MT** **Service Malta**
Tel.: 80062230
E-Mail: kompennass@lidl.com.mt
- IAN 508148_2507

15.2. Importer

Please note that the following address is not the service address. Please use the service address provided in the operating instructions.

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